97050180
rev. 005
12/2010

ITALIANO (Lingua originale).....	3
ENGLISH (translation from the original in Italian).....	77
FRANÇAIS (traduction de l'original en italien).....	151
DEUTSCH (Übersetzung aus dem Original in Italienisch).....	225

Indice

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TABLE OF CONTENTS

1. Safety guidelines.....	78	5.6.	Scaler.....	112	
1.1.	Symbol definition.....	78	5.6.1.	SURGISON 2 surgical detartarer device.....	114
1.2.	Intended use.....	78	5.6.1.1.	Normal operational mode.....	114
1.2.1.	Classification.....	78	5.6.1.2.	BOOST functioning mode.....	115
1.2.2.	Environmental conditions.....	78	5.7.	LEDA curing light.....	117
1.2.3.	Warranty.....	79	5.8.	C-U2 PRO intraoral camera system.....	118
1.2.4.	Disposing the equipment When no longer used.....	79	5.8.1.	OSD control panel (On Screen Display).....	122
1.3.	Safety rules.....	79	5.9.	Integrated laser L810.....	123
1.4.	Cleaning and disinfecting.....	80	5.10.	Peristaltic pump.....	124
			5.11.	Electronic apex locator (LAEC).	125
2. Description of the equipment.....	81	6. Assistant's board operation.....	127		
2.1.	Nameplate.....	81	6.1.	Assistant's touchpad.....	127
2.2.	Dental units.....	81	6.2.	Instruments on assistant's board.....	128
2.3.	Dental chair.....	82	6.3.	Suction tubes.....	128
			6.4.	Instrument tray on assistant's board.....	129
			6.5.	Hydraulic saliva ejector.....	129
3. Turning on the dental operatory.....	82	7. Hydrogroup operation.....	129		
4. NEW SKEMA Dental chair operation.....	82	7.1.	Fill cup and bowl.....	129	
4.1.	Safety devices.....	83	7.2.	Distilled water tank system.....	131
4.2.	Emergency devices.....	83	7.2.1.	Manual disinfecting cycle with hydrogen peroxide.....	132
4.3.	Adjustable headrest.....	83	7.3.	M.W.B. system.....	133
4.4.	Adjustable armrest (optional).....	83	7.4.	AUTOSTERIL automatic disinfection system.....	134
			7.5.	S.H.D. system for flushing the suction tubes.....	136
5. Instrument board operation.....	84		7.6.	TIME FLUSHING automatic cycle.....	137
5.1.	Doctor's control console.....	86	7.7.	Opening/closing the cover at the side of the hydrogroup.....	138
5.1.1.	User interface.....	87	8. Accessories.....	138	
5.1.1.1.	Operator selection.....	88	8.1.	Operating lamp.....	138
5.1.1.2.	General settings.....	88	8.2.	Monitor / Workstation with light pole.....	138
5.1.1.2.1.	Adjusting the display brightness.....	89	8.3.	Air/water/230V quick-connect couplers.....	138
5.1.1.2.2.	Language selection.....	89	9. Maintenance.....	140	
5.1.1.2.3.	Adjusting/engaging the pantograph arm's brake.....	90	9.1.	Instrument maintenance.....	140
5.1.1.2.4.	Entering the operator's personal identification.....	91	9.2.	Draining condensate.....	140
5.1.1.2.5.	Chronometer.....	91	9.3.	Cleaning the surgical suction filters.....	140
5.1.1.2.6.	Setting the date and time.....	92	9.4.	Surgical suction.....	141
5.1.1.2.7.	Selecting your favourite button.....	92	9.5.	CATTANI surgical separator.....	142
5.1.1.2.8.	LAEC setting.....	92	9.6.	Cleaning the turbine return air filter.....	143
5.1.1.3.	Setting the disinfecting systems.....	93	9.7.	METASYS amalgam separator.....	143
5.1.1.3.1.	Setting the AUTOSTERIL cycle (AS).....	93	9.8.	DURR amalgam separator.....	143
5.1.1.3.2.	Setting the Time Flushing cycle (TF).....	94	9.9.	CASTELLINI dental chair model NEW SKEMA.....	143
5.1.1.3.3.	Emptying the M.W.B. system tank.....	95	10. Packing.....	144	
5.1.1.4.	Hydrogroup settings.....	95	10.1.	Transport and storage.....	144
5.1.1.4.1.	Setting water to bowl.....	96	11. Troubleshooting.....	145	
5.1.1.4.2.	Setting water to cup.....	96	12. Specifications.....	146	
5.1.1.4.3.	Automatic bowl motion settings.....	97	12.1.	Overall dimensions: SKEMA 6 RS.....	147
5.1.2.	Setting the dental chair's "rinse" and "home" positions.....	97	12.2.	Overall dimensions: SKEMA 6 CP.....	148
5.1.3.	Setting dental chair positions 1, 2, 3 and 4.....	98	13. SKEMA 6 dental operatory disinfecting and maintenance plan	149	
5.1.4.	Emergency stop button.....	98			
5.2.	Foot control.....	98			
5.2.1.	"Multifunction" foot control.....	98			
5.2.2.	"Push-pedal" foot control.....	101			
5.2.3.	"Power Pedal" foot control.....	102			
5.3.	Syringe.....	104			
5.4.	Turbine.....	105			
5.5.	Micromotor.....	106			
5.5.1	CONSERVATIVE operating mode.....	108			
5.5.2	ENDO operating mode.....	109			
5.5.3	SURGERY operating mode.....	110			
5.5.4	Reduction ratio setup menu.....	111			

1. Safety guidelines

- These instructions describe how to use the SKEMA 6 series dental units correctly.
Carefully read and become familiar with the content of this manual before using the equipment.
- All versions of the SKEMA 6 fully loaded dental units are described herein therefore not all the paragraphs may pertain to the equipment purchased.
- No part of this manual is to be reproduced, stored in a retrieval system or transmitted in any form or by any means, i.e. electronic, mechanical, photocopying, translation or otherwise, without the prior written permission of CASTELLINI S.p.A.
- The information, technical specifications and illustrations contained in this publication are not binding.
CASTELLINI S.p.A reserves the right to make technical improvements and changes without modifying the instructions contained herein.
- The manufacturer has a company policy of continual development. Although every effort is made to keep technical documentation up-to-date at all times some of the instructions, specifications and figures given in this manual may slightly differ from the purchased product. The manufacturer reserves the right to make changes without prior notice.
- The original version of this manual is written in Italian.
- This equipment is equipped with a device that prevents liquid back up.

1.1. Symbol definition

List of symbols used in this document to denote certain conditions:

- | | | | | | | |
|--|---|--|----|--|----|--|
| 1) Type of protection against direct and indirect contact: Class I.
Type of protection against direct and indirect contact: Type B. | 1 | | 8 | | 15 | |
| 2) WARNING!
Failure to observe may result in equipment damage or injury to the user and/or patient. | 2 | | 9 | | 16 | |
| 3) OPERATING INSTRUCTIONS:
This symbol indicates that the user should read and become familiar with the content of the User's manual before attempting to use the relative part of the equipment. | 3 | | 10 | | 17 | |
| 4) NOTE:
Identifies information that is especially important for the user and/or assistant. | 4 | | 11 | | 18 | |
| 5) Earth ground. Functional earth connection. | 5 | | 12 | | 19 | |
| 6) Alternating current | 6 | | 13 | | | |
| 7) Part sterilized in autoclave. | 7 | | 14 | | | |
| 8) Equipment on | | | | | | |
| 9) Equipment off | | | | | | |
| 10) On (a part of the apparatus) | | | | | | |
| 11) Off (a part of the apparatus) | | | | | | |
| 12) Equipment in accordance with essential requirements of directive EEC 93/42 and subsequent changes (Class II equipment). | | | | | | |
| 13) Equipment in accordance with essential requirements of directive EEC 93/42 and subsequent changes . | | | | | | |
| 14) Symbol indicating separate collection for electrical and electronic equipment in compliance with Directives 2002/95/ EC, 2002/96/ EC e 2003/108/ EC. | | | | | | |
| 15) "Warning biological hazard".
It provides information about possible risks of contamination deriving from contact with fluids, storage of infected biological waste. | | | | | | |
| 16) Manufacturer | | | | | | |
| 17) Month and year of constructions | | | | | | |
| 18) Apparatus serial number | | | | | | |
| 19) DVGW mark (quality assurance kitemark regarding supply of drinking water). | | | | | | |

1.2. Intended use

- The SKEMA line 6, SKEMA 6 and SKEMA 6 CP model operatory units are medical devices intended for dental treatment.
- The instrument board may hold up to 6 instruments.
- The assistant's board may be equipped with 2 suction tubes and 3 instruments.
- This equipment must be used only by adequately trained personnel (doctors and paramedics).
- Equipment designed for continuous operation with intermittent loads (see the times for the individual parts in the dedicated sections).

1.2.1. Classification

MEDICAL DEVICES classification

Classification of the dental unit in accordance with the indications given in annex IX of directive 93/42 and subsequent changes/EEC: Class IIa

ELECTRICAL MEDICAL EQUIPMENT classification

Classification of the dental unit in accordance with standard EN 60601-1 for safety of medical equipment: Class I - Type B.

Reference standards

SKEMA 6 series operatory units are designed and constructed in compliance with CEI EN 60601-1, CEI EN 60601-1-2, ISO 7494, ISO 6875 and EN 1717 (type AA and AB) standards as far as the water mains safety devices are concerned.

1.2.2. Environmental conditions

The equipment is to be installed in rooms that satisfy the following requirements:

- temperature between 10 and 40 °C.
- relative humidity between 30 and 75%.
- atmospheric pressure ranging from 700 to 1060 hPa.
- air pressure entering equipment ranging from 6 to 8 bar.
- water hardness entering equipment not over 60 mg/l.
- water pressure entering equipment ranging from 3 to 5 bar.
- water temperature entering equipment not higher than 25 °C.

1.2.3. Warranty

CASTELLINI S.p.A. stands behind its products warranting safety, reliability and performance.

The warranty is valid only under the following terms:

- The conditions given on the warranty certificate are observed.
- Performing yearly scheduled maintenance.
- The equipment is used only as instructed in this manual.
- The electrical wiring in the room in which the equipment is installed must conform to I.E.C. 60364-7-710 (standards for electrical wiring in medical and dental offices).
- A 3x1.5 mm² line protected by a bi-polar cut-out that conforms to applicable standards (10 A, 250 V, cut off power at least 6000 A, distance between contacts at least 3 mm) must be used to feed the equipment.

ATTENTION!

The color of the three wires (POWER, NEUTRAL and EARTH) should satisfy the requirements of current standards.

- Installation, repairs and, in general, any other operations requiring the casing to be opened are to be performed exclusively by personnel authorized by CASTELLINI.

1.2.4. Disposing the equipment When no longer used

As set out in Directives 2002/95/ EC, 2002/96/ EC and 2003/108/ EC, on the restrictions of the use of certain hazardous substances in electrical and electronic equipment along with collection, treatment, recycling and disposal of waste electrical and electronic equipment the latter must be treated as municipal waste, therefore sorted and collected separately. When new equipment of equivalent type is purchased the waste equipment should be returned to the distributor on a one-to-one basis for disposal. As far as reuse, recycling and other forms of waste recovery mentioned above are concerned, the manufacturer is responsible for the actions specified by individual local laws.

Efficient collection of sorted waste separately to recycle and treat waste electrical and electronic equipment aids in preventing negative environmental impacts while protecting human health.

In addition it facilitates recycling of the materials used to construct the equipment.

Illegal waste disposal carries heavy fines defined by local laws.

WARNING!

The crossed out M.W.B.ed bin placed on the equipment indicates that the waste equipment must be collected separately from other waste.

1.3. Safety rules

WARNING!

- **All equipment is permanently installed.**

Depending on the type of chair the unit comes with, refer to the installation DATA given in paragraph "Specifications".

- **CASTELLINI S.p.A shall not be held liable for any personal injury or equipment damage resulting from failure to heed the precaution given above.**

- **Floor conditions.**

The floor (continuous) should meet the load-bearing capacity set forth by DIN 1055, sheet 3. The weight of the dental unit including a 190 kg patient is about 350 kg/sq.m.

- The technician not authorized by CASTELLINI who alters the product replacing parts or components with ones different from those used by the manufacturer shall assume liabilities in an identical manner as the manufacturer.

- **CASTELLINI S.p.A shall not be held liable for any personal injury or equipment damage resulting from failure to heed the precaution given above.**

- **Dental chair**

The maximum chair capacity is 190Kg. This weight must never be exceeded.

- **Tray holders**

The maximum weights that can be held must never be exceeded:

- Instrument tray attached to the instrument board maximum allowable load 2 Kg, evenly distributed.
- Instrument tray attached to the instrument board maximum allowable load 1 Kg, evenly distributed.

- **Connections to external instruments**

The equipment can be hooked up only to other instruments that bear the CE mark.

- **Electromagnetic interferences.**

Use of electrical equipment that does not comply to standard I.E.C. 60601-1-2 in the office or nearby may cause electromagnetic or other types of interferences resulting in dental unit malfunctions.

In these cases, shut off power to the dental unit before using this equipment.

- **Replacing the chucks**

Operate the turbine release and contra angle only once the chuck has come to a complete stop. On the contrary, the locking system will wear down and the chucks can slip off causing injury. Use only high quality chucks with gauged diameter attachment. To check the state of the locking system, make certain the chuck is firmly secured to the instrument every day before starting work. Locking system defects caused by misuse are easily identified and not covered by the warranty.

- **Patients with pace makers and/or hearing aids.**

When treating patients with pace makers and/or hearing aids, take into consideration the effects the instruments may have on pace makers and/or hearing aids. Carefully read technical-scientific information available on this subject.

- **Implants.**

If the dental unit is used for implant operations using separate equipment designed for this purpose, shut off power to the dental chair to avoid unwanted movements resulting from faults and/or accidental start up of the controls.

- Do not forget to turn off the office's water supply and master switch on the equipment before leaving the surgery.

- The equipment is not protected against liquid penetration (IPX 0).

- This equipment is not suitable for use in the presence of a mix of inflammable anaesthetic gas with oxygen or nitrous oxide.

- This equipment must be stored properly so that it is kept in top working order at all times. The manufacturer shall not be held responsible for misuse, carelessness or improper use of the equipment.

- This equipment is to be used exclusively by qualified personnel (doctors and paramedics) with the proper training.

- The user must be present at all times When the equipment is turned on or ready for start-up. In particular, never leave the equipment unattended in the presence of children/the mentally disabled or other unauthorised personnel in general.

Any companions must keep out of the area in which treatment is performed and in any case under the responsibility of the operator. The area in which treatment is performed refers to the space around the dental unit plus 1.5 meters.

- **Water quality**

CASTELLINI S.p.A recommends providing the dental unit with a water disinfecting system to satisfy national and medical legal requirements regarding water coming from dental units. The dental unit user is responsible for the quality of the water and must take alternative actions to maintain the required quality if the dental unit is not equipped with a disinfecting system.

Once the dental unit has been installed in the office and the ducts are filled with water, the equipment is at risk of contamination caused by aquatic microorganisms harbored in the municipal water supply capable of forming biofilm in the ducts. To effectively overcome this problem, it is advisable to install the dental unit only at the moment in which the professional will use it daily and perform the disinfecting procedures at the specified frequency starting from the day in which it is installed.

 **NOTE:** Contact your local dealer or Dental association for more detailed information about national laws and requirements.

1.4. Cleaning and disinfecting

Cleaning is the first step of any disinfecting process. Physically scrubbing with detergents and surface-active substances and rinsing with water removes a considerable amount of micro organisms. If a surface is not clean first, the disinfecting process cannot be successful. If a surface cannot be adequately cleaned, it should be protected with barriers. The external parts of the equipment must be cleaned and disinfected with hospital-grade tuberculocide products effective for HIV, HBV prevention (medium grade) and specific for small surfaces.

Pay strict attention to all the cleaning product manufacturer's instructions, warnings and cautions.

The wide range of drugs and chemical products used in the dentist's office can damage the coated surfaces and plastic parts. Research and tests run show that the surfaces cannot be fully protected against the harsh action of all products available on the market. We therefore recommend protecting with barriers Whenever possible. The harsh actions of chemical products also depend on the amount of time they are left on the surfaces. It is therefore important not to leave the product on the surfaces longer than the time specified by the manufacturer.

Due to the harsh effects of the active substances present in disinfectants, products that contain the following maximum percentage of chemical disinfectants should be used:

- **Ethanol 96%.** Concentration: maximum 30 g per 100 g of disinfectant.
- **Propanole.** Concentration: maximum 20 g per 100 g of disinfectant.
- **Combination of ethanol and propanole.** Concentration: the combination of the two should be maximum 40 g per 100 g of disinfectant.

CASTELLINI S.p.A. has conducted compatibility tests amongst the main disinfectants available on the market and its own products, the most suitable disinfectants are listed below:

- **Plastic surfaces, coated and non-coated surfaces**
 - STER 1 PLUS (Castellini).
- **Upholstery.**
 - RelyOn Virkosept (Du Pont).
 - Plastisept (Alpro).
 - Zeta 3 Foam (Zhermack).

ATTENTION!

The acid used for the mordant process stains VISCO-elastic upholstery. Immediately rinse with plenty of water if acid spatters on the upholstery.

ATTENTION!

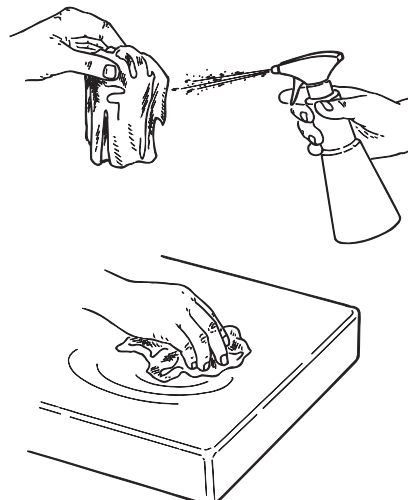
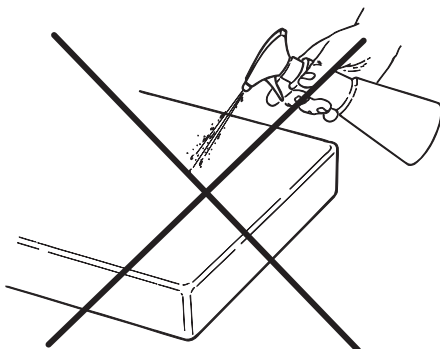
- Do not use products containing isopropyl alcohol (2-propanol, iso-propanol).
- Do not use products that contain sodium hypochlorite (bleach).
- Do not use cleaners that contain phenol.
- Do not spray the selected products directly on the surfaces.
- Never combine products with each other or with liquids other than the products listed above.
- All products must be used as directed by the manufacturer.

Cleaning and disinfecting instructions.

Clean and disinfect with single-use non-abrasive paper (avoid using recycled paper) or sterile gauze.
Do not use sponges or in any case, any material that can be reused.

ATTENTION!

- Shut off the dental unit prior to clean and disinfecting the external parts.
- All material used to clean and disinfect must be thrown away.

2. Description of the equipment

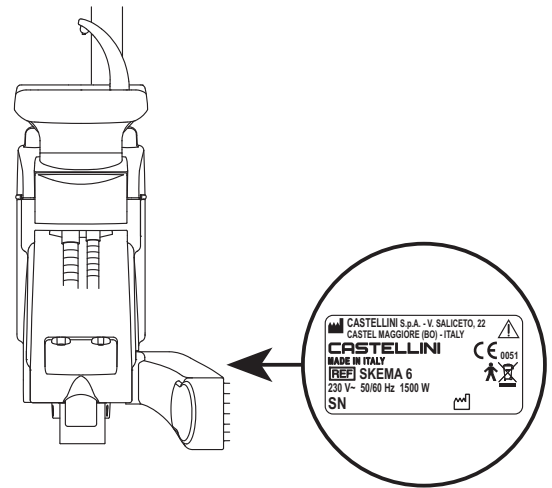
2.1. Nameplate

Dental operatory

The nameplate is found on the arm between the dental chair and hydro-group.

Data given on plate:

- Manufacturer's name
- Name of equipment
- Voltage
- Type of current
- Frequency
- Maximum power absorbed
- Serial number
- Month and year of construction

2.2. Dental units

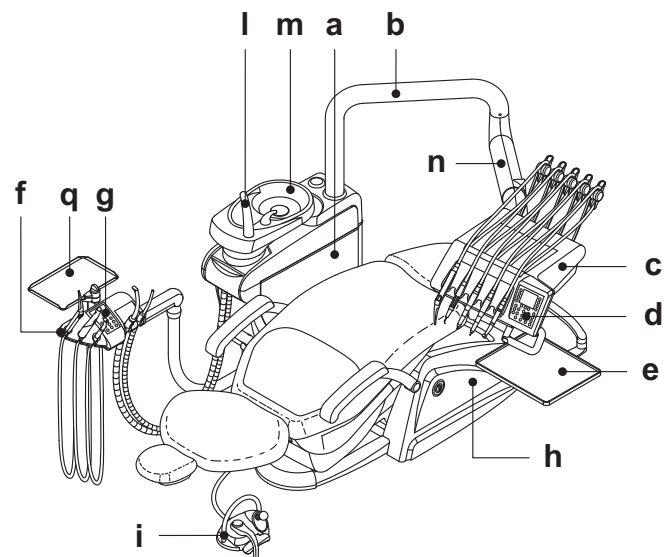
Dental units SKEMA 6 are available in the following versions:

RS version

"RS" version instrument board (instruments will return to their original position through the pulling action of the spring-operated arms) attached to a double supporting arm, one of which is articulated while the other is autobalancing.

Description of equipment:

- [a] Hydrogroup
- [b] Adjustable arm
- [c] Instrument board
- [d] Doctor's console
- [e] Tray holder
- [f] Assistant's board
- [g] Assistant's control console
- [h] Utility service center
- [i] Multifunction foot control
- [l] Water to cup
- [m] Bowl
- [n] Autobalancing arm
- [q] Instrument tray on assistant's board (optional)


CP version

Instrument board version CP (the instruments are placed vertically in position) mounted on double arm of which one is articulated and self-balanced.

Description of equipment:

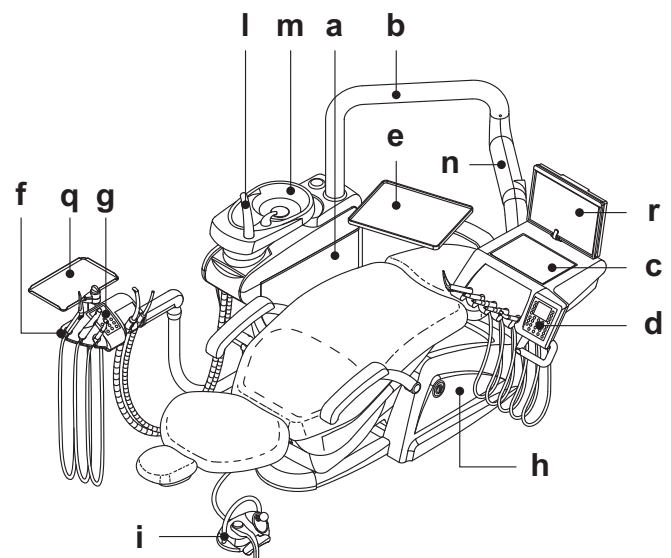- [a] Hydrogroup
- [b] Adjustable arm
- [c] Instrument board
- [d] Doctor's console
- [e] Tray holder (optional).
- [f] Assistant's board
- [g] Assistant's control console
- [h] Utility service center
- [i] Multifunction foot control
- [l] Water to cup
- [m] Bowl.
- [n] Autobalancing arm
- [q] Instrument tray on assistant's board (optional)
- [r] Viewer for panoramic x-rays (optional).

2.3. Dental chair

Description of the chair

- [a] Headrest
- [b] Back
- [c] Left mobile armrest (optional).
- [d] Right mobile armrest (optional).
- [e] Safety foot board
- [f] Sliding footrest

Operating times

The operating and rest times are as follows: **work 25 sec. - rest 10 min.**

Maximum weight capacity.

The maximum chair capacity is 190 Kg.

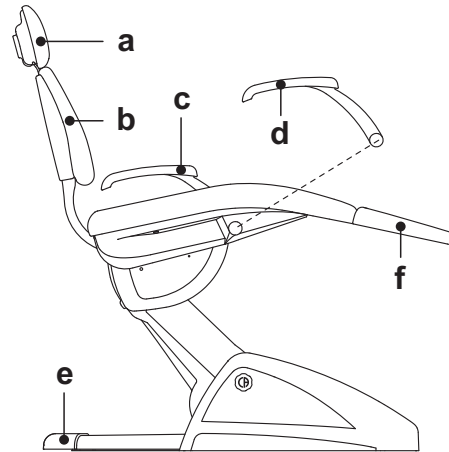
ATTENTION!

Do not exceed this value.

Note about biocompatibility.

For skai upholstery, only CASTELLINI original spare parts must be used.


3. Turning on the dental operatory

Use the main switch [f]:

Equipment on:

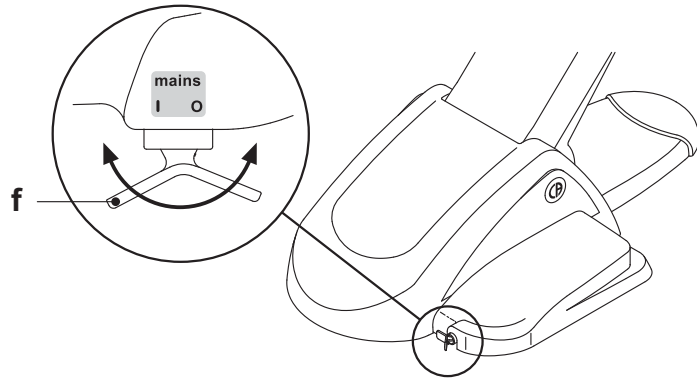- electrical power supplied
- pneumatic system connected
- hydrogroup connected

0

Equipment off:

- electrical power not supplied
- pneumatic system not connected
- hydrogroup not connected

4. NEW SKEMA Dental chair operation

The NEW SKEMA dental chair can be moved as follows:

- Chair seat up/down
- Back up/down with inclination of the chair seat (Trendelenburg compensated)

The dental chair can be operated from the following places:


- Instrument board [a] (see par. 5)
- Multifunction foot control [b] (see par. 5.2)
- Assistant's board [c] (see par. 6)

The footrest [f] can be extended approximately 10 cm.

Dental chair motion blocked.

With the instruments in place, hold down button MED for at least 5 seconds, to disable/enable dental chair movements.

The relative icon (A) appears at the top of the console display to signal movement has been disabled.

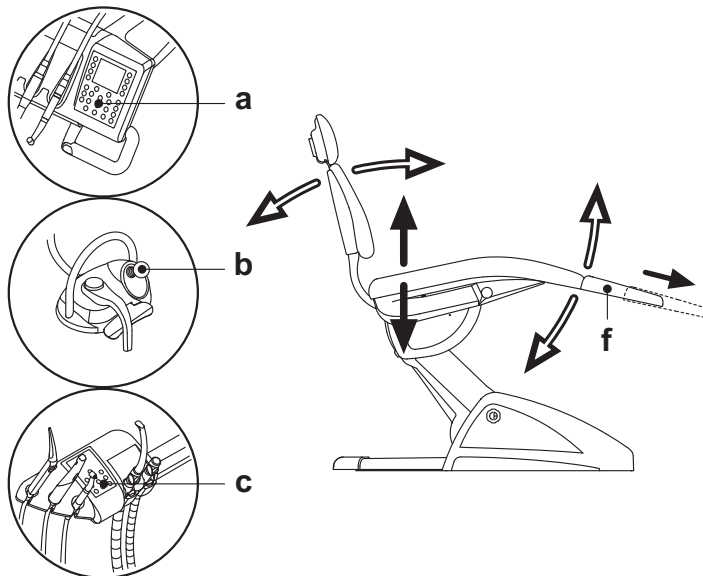
4.1. Safety devices

The equipment is supplied with the following safety devices:

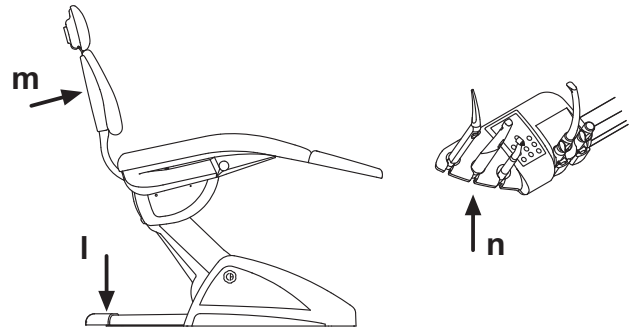- The floor box is equipped with a device [l] that immediately stops the dental chair from moving down in the presence of an obstacle and automatically lifts it up to free the obstacle.
- The backrest is equipped with a device [m] that immediately stops the backrest from moving down in the presence of an obstacle and automatically lifts it up to free the obstacle.
- The arms of the assistant's board are equipped with a safety device [n] that immediately stops the dental chair from moving down if an obstacle is encountered and automatically moves it up to clear the obstacle.


4.2. Emergency devices

ATTENTION!

Use the devices below When movement of the equipment needs to be blocked:

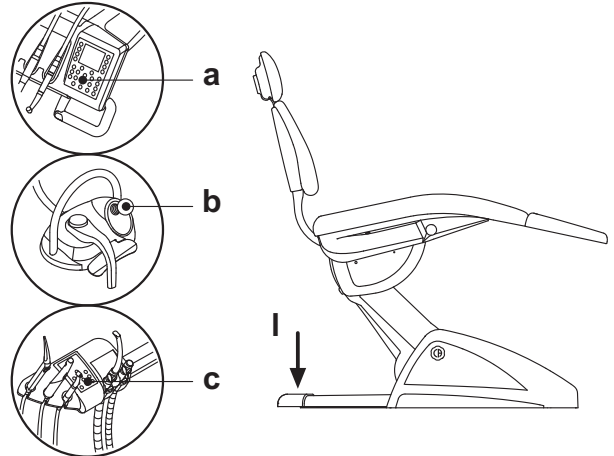- **Dental chair control buttons [a] or [c].**
Pressing any dental chair button blocks all movements.
- **Foot control [b].**
When the foot control is actuated, all movements of the equipment are blocked.
- **Foot board [i].**
Foot board activated: all movements are blocked.

4.3 Adjustable headrest

Two types of headrests are available:

- 1 Cushion secured in place by hand
- 2 Cushion secured pneumatically

Adjusting the height of the headrest

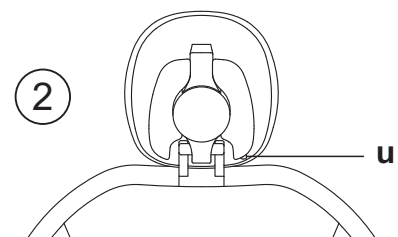
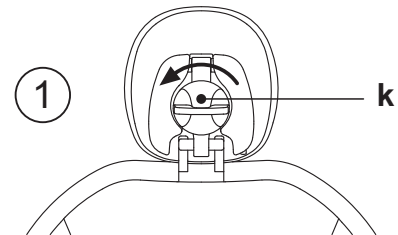
The head rest blade is positioned through a magnetic clutch. The operator should pull up and/or push down the headrest until it is in the desired position.

Adjusting the cushion

- secured in place by hand (1)
Turn the knob [k] counter-clockwise and place the cushion in the desired position. Once adjusted, tighten the knob.
- secured in place pneumatically (2)
Press and hold down the lock button [u] and place the cushion in the desired position. Once adjusted, simply release the button [u] to block the cushion in place.

NOTE: the locking device is active only When the air circuit is pressurized.

Safety guidelines

ATTENTION!

- **Maximum weight supported by headrest: 30 kg**
- **Do not change the position of the cushion without deactivating the locking device.**
- **In order to prevent the headrest from moving in an uncontrolled manner, hold it with both hands before disengaging the locking device.**

4.4. Adjustable armrest (optional)

Pushing down the adjustable armrest

Turn the adjustable armrest clockwise to move it down so that the patient can more easily get on and off the chair.

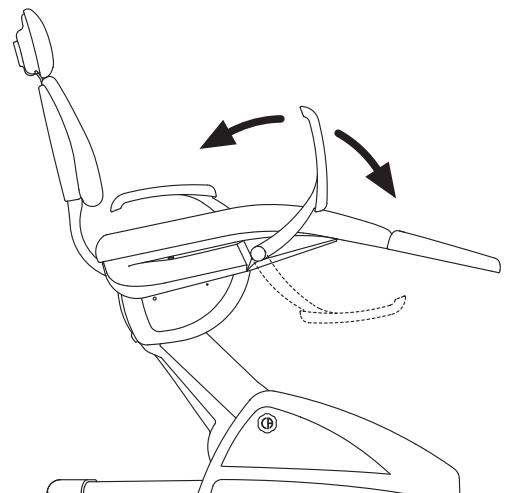
Taking off the adjustable armrest.

Put the armrest in a vertical position and pull it off.

ATTENTION!

- **Maximum weight supported by armrest: 35 kg**

5. Instrument board operation

Layout of instruments

The positions the instruments are placed in on the board are determined by the customer at the time of order.

Starting the instruments.

- The syringe is always on (see paragraph 5.3.).
- The curing light is turned on with the key When the instrument is withdrawn (see paragraph 5.7.).
- Intraoral camera and C-U2 PRO turn on When the instrument is extracted (see paragraph 5.8.).
- Once taken out, the L810 integrated laser is operated with the foot control (see paragraph 5.9).
- Once picked up, all the instruments are operated with the foot control. (see paragraph 5.2.).

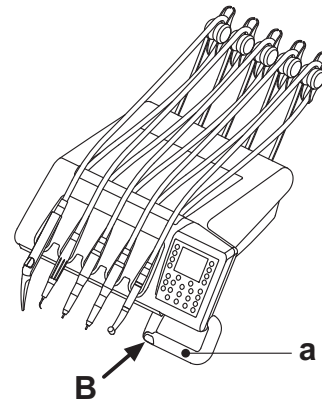
Simultaneous use of the instruments

A device sees that the instruments cannot be used simultaneously. The first instrument removed is operative while those removed there after are deactivated by this device. This device allows the chuck to be replaced in one instrument while another is used on the patient.

Positioning the instrument board.


The instrument board can move in all directions. To adjust the height of the board and/or direct it horizontally, simply grasp the handle [a].

 **NOTE:** Grasp the handle placing your thumb on point (B) to release the air lock for the pantograph arm.

Instrument return stop arms (only in SPRIDO versions)

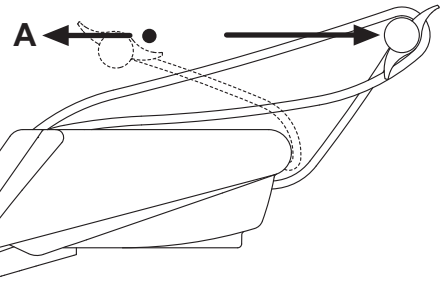If this device is provided, the instrument return arm can be locked in the instrument extracted position. When the device is used a click is heard about 2/3 of the total arm travel. To go back to the original condition, simply move the arm to the end of its travel [a].


Tray holder attached to instrument board.

The instrument holder [q] is removable and can be sterilized in autoclave at 135°.

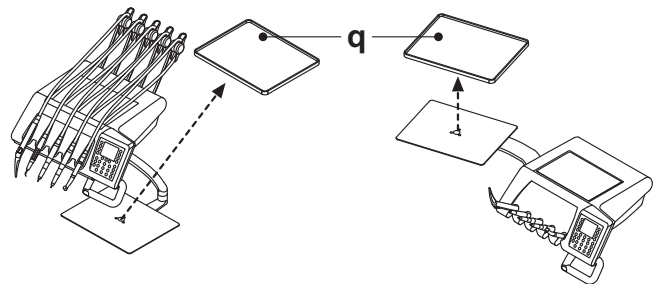ATTENTION!
Maximum weight allowed on instrument tray: 2 Kg evenly distributed.

Cleaning the instrument board handle.

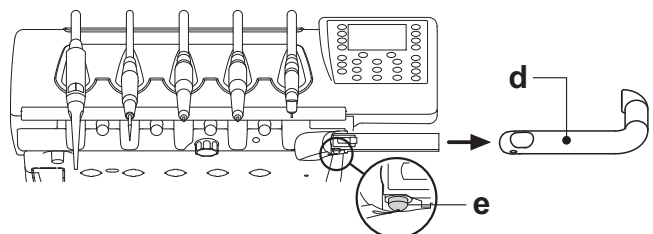
Press the release button [e] and then pull the handle [d] outwards and remove it.

Clean the instrument board handle with a suitable product approved against bacteria and fungi and effective against viruses (including HBV and HIV). Pay strict attention to all the cleaning product manufacturer's instructions, warnings and cautions; see paragraph 1.4.

 **NOTE:** The product should be sprayed on soft disposable paper towel or sterile gauze.

ATTENTION!
Do not use acids or harsh products.

Cleaning the instrument board

Clean the instrument board; use a suitable product approved against bacteria and fungi and effective against viruses (including HBV and HIV); see paragraph 1.4.

ATTENTION!

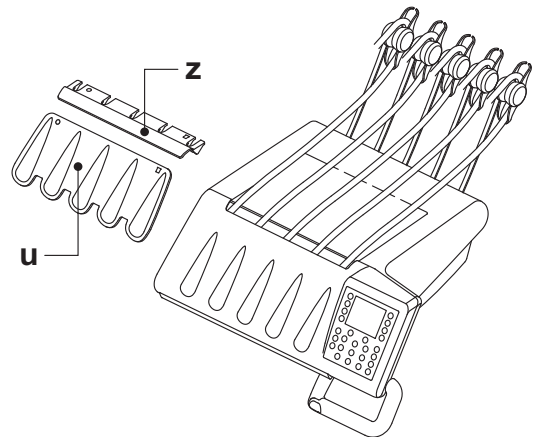
Do not use acids or harsh cleaning products.

NOTE RS version boards: The instrument holder [u] and cord cover [z] are autoclavable at 135°.

NOTE CP version boards: the board's protection [v] can be sterilized in autoclave at 135°.

Removable instrument cords

All the instruments have removable cords to ease cleaning.

Follow the directions given below to remove the instrument cords:

RS version boards

- Open the front of the board, after releasing it with the button [k] found behind the console as shown in the figure.
- Loosen the plastic ring nuts and remove the cords.

CP version boards

Loosen the plastic ring nuts found underneath the board and then remove the cords.

ATTENTION!

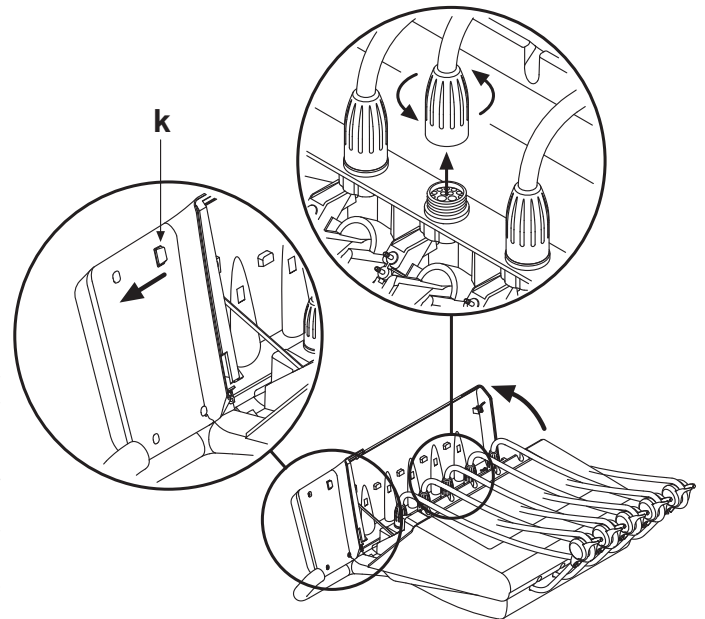
- Shut off the operatory unit before attempting to take off the cords.
- After shutting off the operatory unit, empty the syringe's ducts by pressing and holding down the relative air and water buttons directly on the bowl until water spray is no longer present.
- The cords of the TURBINE, MICROMOTOR and SCALER contain water, therefore hold the end of the cord on the handpiece side over the bowl! When removing the cord.
- When putting a cord back on, make certain the contacts are perfectly dry and the plastic ring nut is tight.
- Each cord may be remounted only in the position for the corresponding instrument.

Clean the instrument board with a suitable product approved against bacteria and fungi and effective against viruses (including HBV and HIV); see paragraph 1.4.

ATTENTION!

The instrument cords are NOT suitable for autoclave or cold sterilization.

5.1. Doctor's control console

Description of the buttons:

- | | | | |
|--|---|--|---|
| | Menu change/call up button. | | Dental chair functions save button |
| | Multifunction button. | | Automatic return position button |
| | Fast selection button, 100% of settable values | | Rinse position button |
| | Fast selection button, 50% of settable values | | Emergency position button |
| | Fast selection button, 1% of settable values | | Programmed position 1 retrieval and chair base up button. |
| | INCREASE button: increases the values | | Programmed position 2 retrieval and chair back up button. |
| | DECREASE button: decreases the values | | Programmed position 3 retrieval and chair base down button. |
| | Fiber optic lighting on/off button (on/off) | | Programmed position 4 retrieval and chair back down button |
| | Scaler ENDO function/micromotor reverse button | | |
| | Dental operatory light on/off button | | |
| | Cup fill button | | |
| | Bowl rinse function | | |
| | Only water to spray to instruments | | |
| | S.S.S. system on/off | | |
| | Assistant call button | | |
| | Bowl counter-clockwise rotation button (operative only with powered bowl) | | |
| | Bowl clockwise rotation button (operative only with powered bowl) | | |
| | Peristaltic pump start button (operative only if the peristaltic pump is present) | | |

NOTE: functioning of buttons that move the chair

- Briefly pressed: programmed position automatically reached.
- Held down: manual positioning.

Description of the icons.

- 1 S.S.S. system on.
- 1 Distilled water tank NOT pressurized.
- 1 (Black) Distilled water tank in reserve and active
- 1 (Red) Distilled water tank in reserve but NOT active
- 2 M.W.B. system on
- 2 (Black) Hydrogen peroxide tank in reserve and active
- 2 (Red) Hydrogen peroxide tank in reserve but NOT active
- H₂O COLD water to cup.
- H₂O WARM water to cup.
- H₂O HOT water to cup.

- OFF Pantograph arm brake engaged
- Peristaltic pump on
- ? Dental chair storage activated.
- 1 Dental chair automatic program 1.
- 2 Dental chair automatic program 2.
- 3 Dental chair automatic program 3.
- 4 Dental chair automatic program 4.
- * Dental chair rinse position automatic program.
- Dental chair reset position automatic program.
- Dental chair motion blocked.

5.1.1. User interface

When turned on, the dental unit performs a brief autodiagnosis cycle that ends when the main screen containing the name of the operator last set is displayed.

As of this moment a number of settings can be edited from user-friendly menus (see diagram).

How to navigate the various menus.

- Press button MENU to go from one icon to another.
- After selecting an icon, press multifunction button OK to enter the relative submenu.
- Press multifunction button ESC to go back to the previous screen.
- Press buttons INCREASE and/or DECREASE to edit the values.

Layout of the user interface menu.

The user interface menu is structured as shown in the diagram and includes the following menus:

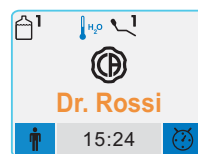
- 5.1.1.1. Operator settings.
- 5.1.1.2. GENERAL SETTINGS.
 - 5.1.1.2.1. Adjusting display contrast.
 - 5.1.1.2.2. Language selection.
 - 5.1.1.2.3. Adjusting/engaging the pantograph arm's brake.
 - 5.1.1.2.4. Entering the operator's personal identification.
 - 5.1.1.2.5. Chronometer.
 - 5.1.1.2.6. Setting the time.
 - 5.1.1.2.7. Selecting your favourite button.
 - 5.1.1.2.8. LAEC setting (only with LAEC).
- 5.1.1.3. SETTING THE DISINFECTING SYSTEMS.
 - 5.1.1.3.1. Setting the disinfecting cycle (only with AUTOSTERIL system).
 - 5.1.1.3.2. Setting the Time Flushing cycle (only with Time Flushing).
 - 5.1.1.3.3. Emptying the M.W.B. tank (only with M.W.B. system).
- 5.1.1.4. SETTING THE HYDROGROUP.
 - 5.1.1.4.1. Setting water deliver to bowl.
 - 5.1.1.4.2. Setting water deliver to cup.
 - 5.1.1.4.3. Movement control (only with powered bowl).

Error message.

During the autodiagnosis cycle the dental unit may detect malfunctions inside the system.

If an error message appears, it remains on the display until the operator presses the OK button. (see paragraph 11.).

If the malfunction is not dangerous, the dental unit remains operative.

5.1.1.1.

5.1.1.2.

5.1.1.3.

5.1.1.4.

5.1.1.4.1.

5.1.1.4.2.

5.1.1.4.3.

5.1.1.3.1.

5.1.1.3.2.

5.1.1.3.3.

5.1.1.2.1.

5.1.1.2.2.

5.1.1.2.3.

5.1.1.2.4.

5.1.1.2.5.

5.1.1.2.6.

5.1.1.2.7.

5.1.1.2.8.

5.1.1.1. Operator selection

The console of SKEMA 6 model dental units allows 3 different operators to be set.

The information given below can be set for each user:

- User name
- Adjusting the power of the turbine and scaler.
- 4 operating modes for the electric micromotor.
- Turn on and adjustment of the fiber optics of each instrument.
- Turbine and scaler power ON/OFF or incremental control
- Automatic dental chair programs.
- Alarm threshold of LAEC electronic apex locator

User selection.

The console allows 3 different users to be selected. Each user is identified by a specific name that consists of up to 20 characters (see paragraph 5.1.1.2.4.)

Press the multifunction button for icon (👤) on the main screen several times until the desired operator is reached.

NOTE:

The user is changed whenever multifunction button is pressed.

5.1.1.2. General settings

Perform the following steps from the main screen:

- Press button MENU to enter the main menu that contains the following icons:

GENERAL SETTINGS

DISINFECTING SYSTEM SETTINGS

HYDROGROUP SETTINGS

- Press button MENU to scroll the icons, selecting the one for " GENERAL SETTINGS" .
- Press multifunction button OK to enter the selected submenu that contains the icons for the following additional submenus:

Adjusting the display brightness

Language selection.

Adjusting/engaging the pantograph arm's brake.

Entering the operator's personal identification.

Chronometer.

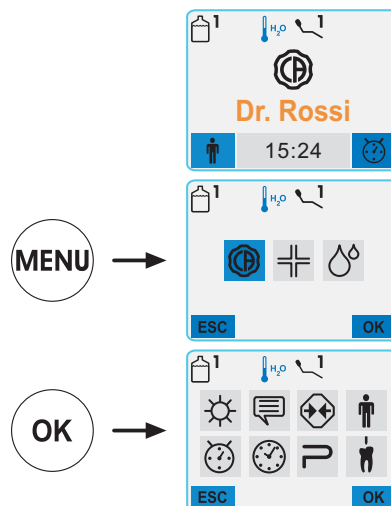
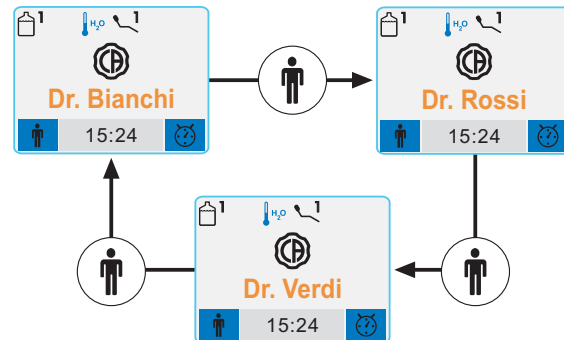
Setting the date and time

Selecting your favourite button.

LAEC settings

5.1.1.2.1. Adjusting the display brightness

This setting is used for all the operators.

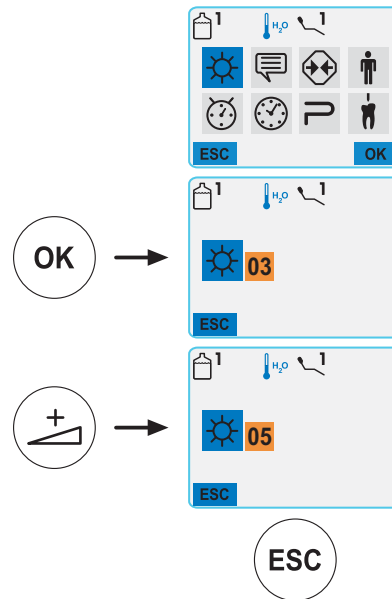
Perform the following steps from the GENERAL SETTINGS menu:

- Press button MENU to scroll the icons and select the one for "Display brightness adjustment" .
- Press multifunction button OK to enter the selected submenu.
- Adjust the display brightness with the increase and/or decrease buttons.

 **NOTE:** : the value which can be entered ranges from 1 to 10.

Use the fast selection buttons to automatically set the following values:

- MAX button: level 10,
 - MED button: level 5,
 - MIN button: level 1.
- Once the brightness has been adjusted, press function key ESC to quit the submenu automatically saving the brightness setting.


5.1.1.2.2. Language selection

This setting is used for all the operators.

Perform the following steps from the GENERAL SETTINGS menu:

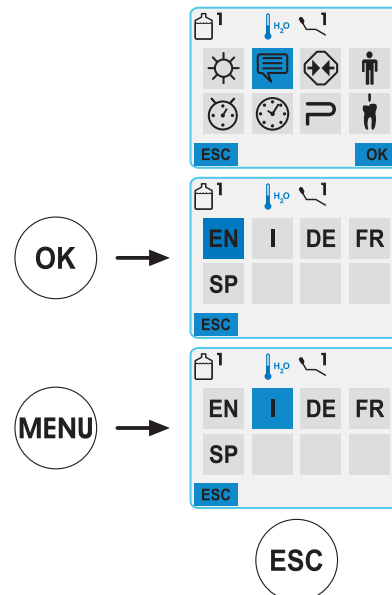- Press button MENU to scroll the icons, selecting the one for " Language selection" .
- Press multifunction button OK to enter the selected submenu.
- Use button MENU to scroll the language icons.
- Once the desired language has been found, press multifunction button ESC to quit the submenu automatically saving the data entered.

5.1.1.2.3. Adjusting/engaging the pantograph arm's brake

This setting is used for all the operators.

Perform the following steps from the GENERAL SETTINGS menu:

- Press button MENU to scroll the icons, selecting the one for "Adjusting the arm's brake".
- Press multifunction button OK to enter the selected submenu which contains the icons for the adjustments that can be made:

Engaging the pantograph arm's brake.

Engaging the pantograph arm's brake.

- Once adjusted as required, press multifunction button ESC to quit the submenu automatically saving the settings entered.

Engaging the pantograph arm's brake.

- Press button MENU to scroll the icons, selecting the one for "Engaging the pantograph arm's brake".
- Press buttons INCREASE and/or DECREASE to enable or disable release of the pantograph arm's brake.

 **NOTE:** Icon (A) found at the top of the console display signals if the brake cannot be disengaged. .

IMPORTANT INFORMATION!

As an additional safety measure, this operation is mandatory if an external electric scalpel needs to be used.

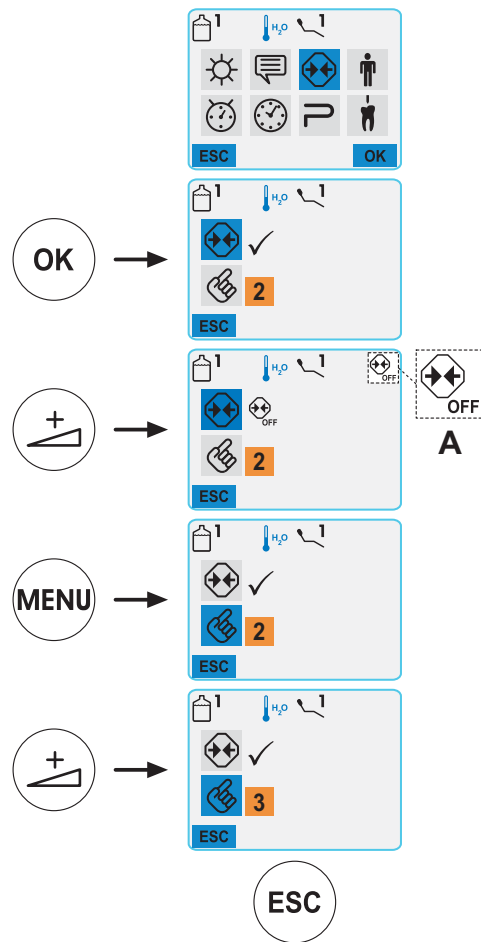
Adjusting brake sensitivity.

- Press button MENU to scroll the icons, selecting the one for "Adjusting brake sensitivity".
- Press buttons INCREASE and/or DECREASE to adjust sensitivity for brake activation.

 **NOTE:** the value which can be entered ranges from 1 to 5.

Use the fast selection buttons to automatically set the following values:

- MAX button: level 5,
- MED button: level 3,
- MIN button: level 1.

5.1.1.2.4. Entering the operator's personal identification

Perform the following steps from the GENERAL SETTINGS menu:

- Press button MENU to scroll the icons, selecting the one for "Entering the operator's personal identification".
- Press multifunction button OK to enter the selected submenu.
- To confirm the personal identification entered (max. 20 letters) simply quit the submenu by pressing multifunction button ESC.

NOTE: any edited personal identification always refers to the operator set on the main screen.

How to enter personal identification.

- Press multifunction button (▶) to move the cursor.
- Press button MENU to go from one group of letters to another.
- Press buttons INCREASE and/or DECREASE to enter a letter.
NOTE: Uppercase letters are present in each group of letters.
- To delete any errors, use icon (←C) canceling from right to left.

5.1.1.2.5. Chronometer

This setting is used for all the operators.

Perform the following steps from the GENERAL SETTINGS menu:

- Press button MENU to scroll the icons, selecting the one for "Chronometer".
- Press multifunction button OK to enter the selected submenu.
- Once the time has been set, press multifunction button (▶) to start the countdown.

NOTE: at this point, it is possible to quit this menu by pressing multifunction ESC without interrupting countdown.

- Press multifunction button (■) to interrupt countdown.
- Once the set time has elapsed, the dental unit beeps and the "Chronometer" menu appears on the console display again.

Press multifunction button ESC to shut off the beep.

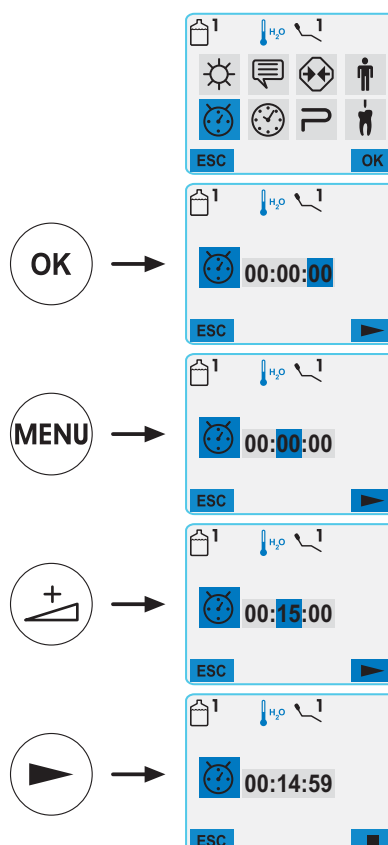
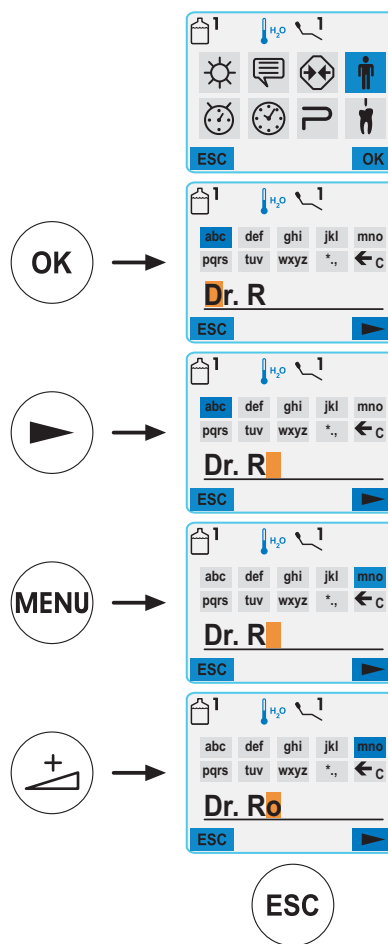
NOTE: the last time set is saved.

How to set the time.

- To change over from "seconds" to "minutes" or "hours" press button MENU.
- Press buttons INCREASE and/or DECREASE to change the time.

NOTE: the time which can be set ranges from 00.00.00 to 10.59.59.

- To reset the chronometer, go to "hours" and hold down button DECREASE until the value reaches 00.

5.1.1.2.6. Setting the date and time

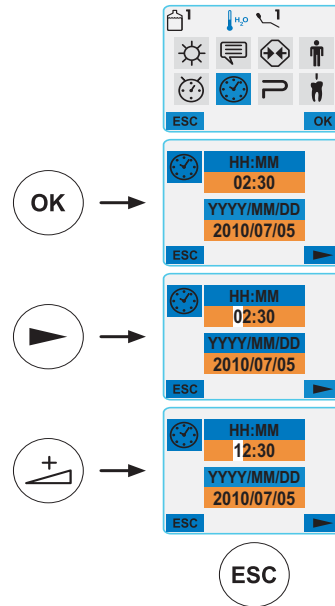
This setting is used for all the operators.

Perform the following steps from the GENERAL SETTINGS menu:

- Press button MENU to scroll the icons and select the one for "Date and time setting".
- Press multifunction button OK to enter the selected submenu.
- Once the date and time have been set, press function key ESC to quit the submenu automatically saving the changes made.

How to set the time.

- Press the function key (▶) to select the number to be modified.
- Press buttons INCREASE and/or DECREASE to modify the number selected.

5.1.1.2.7. Selecting your favourite button

This submenu is used to select the function to be assigned to the right multifunction button when the main page is displayed.

This setting is used for all the operators.

Perform the following steps from the GENERAL SETTINGS menu:

- Press button MENU to scroll the icons, selecting the one for "Selecting your favourite button".
- Press multifunction button OK to enter the selected submenu that contains the icons for the functions that can be selected:

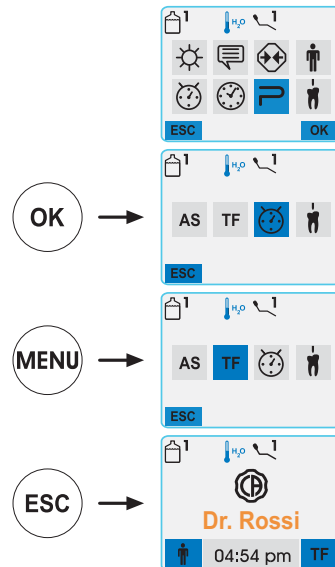
AS AS cycle setting

TF TF cycle setting

Chronometer

LAEC setting

- Press button MENU to scroll the icons, selecting the desired function.
- To confirm the selected function, press multifunction button ESC to quit this submenu.

5.1.1.2.8. LAEC setting

This submenu is used to set the alarm threshold of the LAEC electronic apex locator (see paragraph 5.11.).

Perform the following operations from the MAIN SETUP menu:


- Press button MENU to scroll the icons and select the one for "LAEC settings".
- Press function key OK to open the selected submenu.
- Use the increase and/or decrease buttons to set the alarm threshold. An orange dash on the left-hand bar graph will display the selected value.

NOTE: a value ranging from 0 to +2 can be used.

- Press function key ESC to quit this submenu and confirm the selected value.
- By pressing the MENU button, the alarm signal can be activated/deactivated when the set alarm threshold is reached.

The icon of the selected mode is shown next to the MENU button:

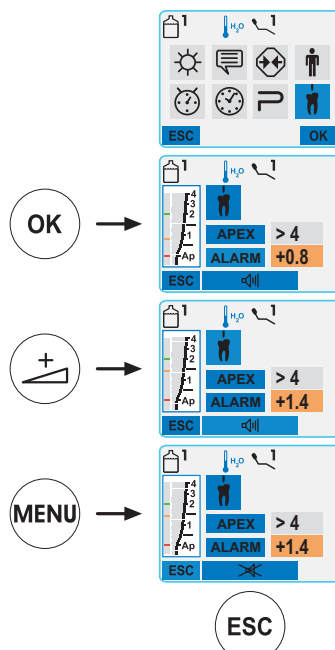- icon : alarm active.
- icon : alarm off.

5.1.1.3. Setting the disinfecting systems

Perform the following operations from the main menu:

- Press button MENU to scroll the icons, selecting the one for "DISINFECTING SYSTEM SETTINGS".
- Press multifunction button OK to enter the selected submenu that contains the following icons:

AS AS cycle setting

TF TF cycle setting

MWB Emptying the M.W.B. system tank

5.1.1.3.1. Setting the AUTOSTERIL cycle (AS)

 **NOTE:** this menu can be accessed only if the AUTOSTERIL system is present.

This setting is used for all the operators.

Perform the following steps from the DISINFECTING SYSTEM SETTINGS menu:

- Press button MENU to scroll the icons, selecting the one for "AS cycle setting".
- Press multifunction button OK to enter the selected submenu.

 **NOTE:** this submenu can also be accessed by holding down button AS found on the assistant's board for at least 2 seconds.

 **NOTE:** This submenu cannot be accessed if the reserve level is present in the liquid disinfectant tank (see paragraph 7.4.), if an instrument is withdrawn or if there is a fault in the MWB system. A beep will signal that the submenu cannot be accessed.

- Set the time liquid disinfectant has to be left in with the INCREASE and/or DECREASE buttons.

 **NOTE:** the time may range from at least 5 minutes to at most 30 minutes, with 30 second intervals.

Use the fast selection buttons to automatically set the following values:

- MAX button: 30 minutes
- MED button: 15 minutes
- MIN button: 5 minutes

ATTENTION!

Recommended time 3% hydrogen peroxide (10 volumes) should be left in: 10 minutes.

- Withdraw the instruments to be treated (the corresponding icon will appear on the display):

S1: syringe on instrument board.

A: instrument in position A

B: instrument in position B

C: instrument in position C

D: instrument in position D

S2: syringe on assistant's board.

F: instrument on assistant's board.

CA: suction tubes.

BC: water to cup duct.

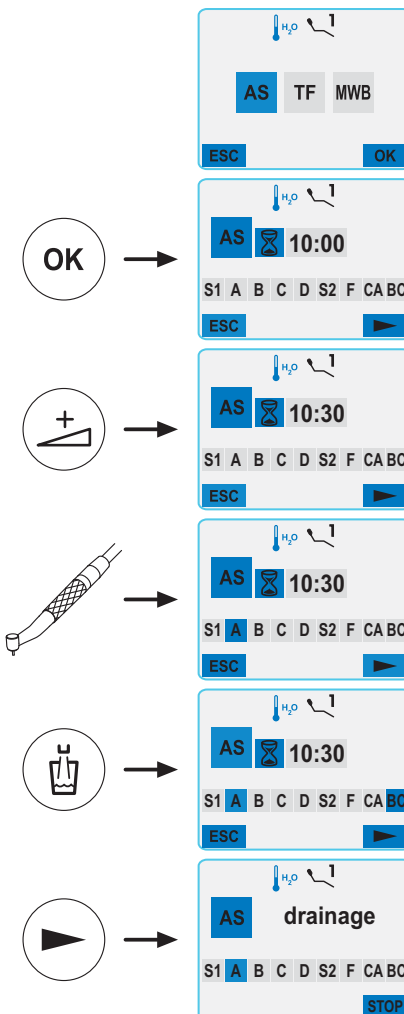
 **NOTE:** Use the "Water delivery to cup" button to select/deselect disinfection of the water to cup duct.

 **NOTE:** if the S.H.D. flushing system for the suction tubes is provided, it is possible to flush the latter by simply placing them in the fittings (see paragraph 7.5.).

- Press multifunction button (▶) to start the disinfection cycle (see paragraph 7.4.).

 **NOTE:** the disinfecting cycle can also be started by briefly pressing button AS found on the assistant's board.

 **NOTE:** While the disinfecting cycle is being set, hold down button AS found on the assistant's board for at least 2 seconds to immediately quit the submenu.

5.1.1.3.2. Setting the Time Flushing cycle (TF)

NOTE: this menu can be accessed only if the Time Flushing system is present.

This setting is used for all the operators.

Perform the following steps from the DISINFECTING SYSTEM SETTINGS menu:

- Press button MENU to scroll the icons, selecting the one for " TF cycle setting " .
- Press multifunction button OK to enter the selected submenu.

NOTE: this submenu can also be accessed by briefly pressing button AS found on the assistant's board.

NOTE: This submenu cannot be accessed if the reserve level is present in the S.S.S. system tank (see paragraph 7.2.). A message will appear on the console display and a beep will be given to signal that the submenu cannot be accessed.

- Set the washing time with the INCREASE and/or DECREASE buttons.

NOTE: the time may range from at least 1 minute to at most 5 minutes, with 1 minute intervals.

Use the fast selection buttons to automatically set the following values:

- MAX button: 5 minutes
- MED button: 3 minutes
- MIN button: 1 minutes

NOTE: it is advisable not to set a time longer than 2 minutes with the S.S.S. system tank.

- Withdraw the instruments to be treated (the corresponding icon will appear on the display):

S1: syringe on instrument board.

A: instrument in position A

B: instrument in position B

C: instrument in position C

D: instrument in position D

S2: syringe on assistant's board.

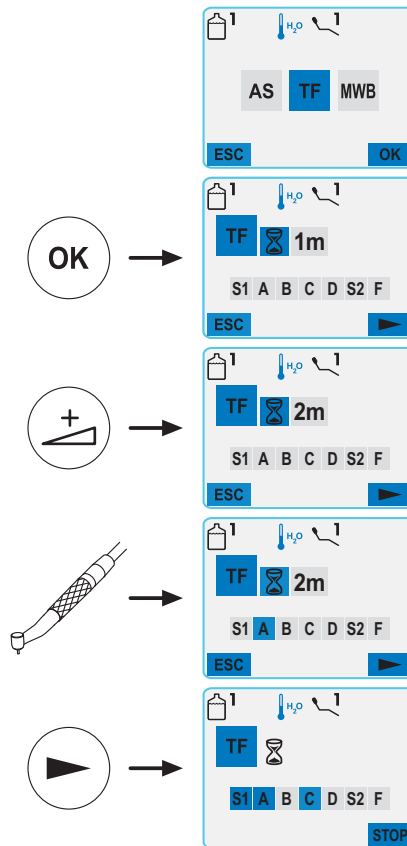
F: instrument on assistant's board.

NOTE: the Time Flushing cycle will not start if at least one instrument is not selected.

- Press multifunction button (►) to start the Time Flushing cycle (see paragraph 7.6.).

NOTE: the Time Flushing cycle can also be started by briefly pressing button AS found on the assistant's board.

NOTE: While the disinfecting cycle is being set, hold down button AS found on the assistant's board for at least 2 seconds to immediately quit the submenu.

5.1.1.3.3. Emptying the M.W.B. system tank

NOTE: this menu can be accessed only if the M.W.B. system is present.

This function is used to empty the M.W.B. system's water circuit (see paragraph 7.3.) if the dental unit is not going to be used for several days or if you want to drain the water present in the system.

Perform the following steps from the DISINFECTING SYSTEM SETTINGS menu:

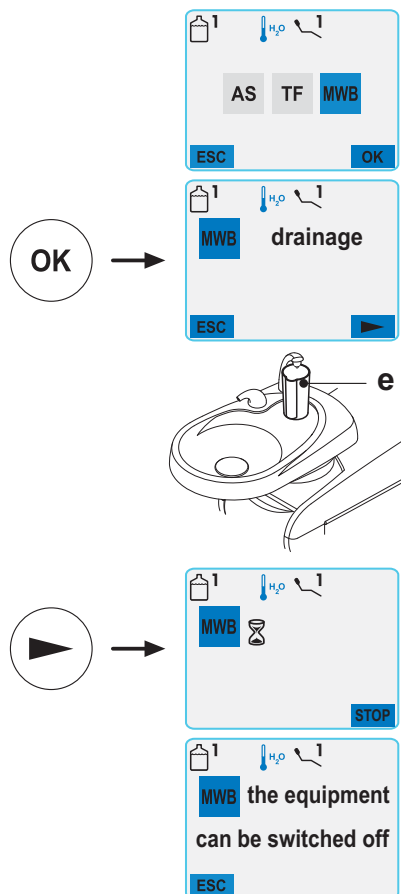
- Press button MENU to scroll the icons, selecting the one for "Emptying the M.W.B. system tank".
- Press multifunction button OK to enter the selected submenu.
- Put the cup [e] provided under the spout.
- Press multifunction button (▶) to start the the drain cycle.

NOTE: the drain cycle will not start if the S.S.S. system is turned on or a fault is present in the M.W.B. system.

- Press multifunction button STOP to stop the drain cycle.

NOTE: the system automatically resets itself.

- Once the drain cycle has been completed, it is possible to shut off the dental unit or press multifunction button ESC to reset the system to restart work.

5.1.1.4. Hydrogroup settings

Perform the following steps from the main menu:

- Press button MENU to scroll the icons, selecting the one for "HYDRO-GROUP SETTINGS" .
- Press multifunction button OK to enter the selected submenu that contains the following icons:

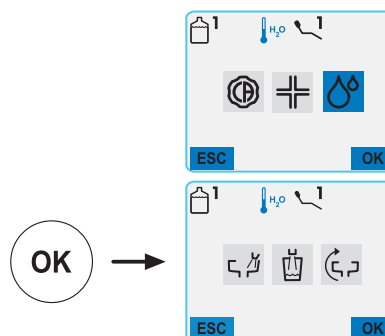
Water to bowl settings

Water to cup settings

Automatic bowl motion settings

5.1.1.4.1. Setting water to bowl

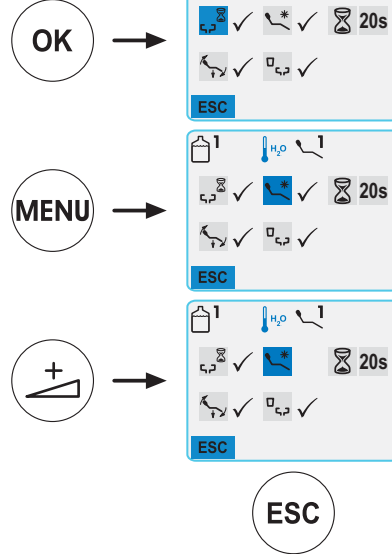
This setting is used for all the operators.

Perform the following steps from the HYDROGROUP SETTINGS menu:

- Press button MENU to scroll the icons, selecting the one for " Setting water to bowl"
- Press multifunction button OK to enter the selected submenu which contains the following icons:

-  Setting of timed or ON/OFF bowl flushing
-  Bowl flushing controller with dental chair brought to rinse position
-  Bowl flushing time setting
-  Bowl flushing controller with dental chair brought to reset position
-  A Bowl flushing controller with cup call

- Press button MENU to go from one icon to another.
- Press buttons INCREASE and/or DECREASE to select/deselect a controller or change the bowl flushing time.
- To confirm the settings, simply quit the submenu by pressing multifunction button ESC.

5.1.1.4.2. Setting water to cup

This setting is used for all the operators.

Perform the following steps from the HYDROGROUP SETTINGS menu:

- Press button MENU to scroll the icons, selecting the one for " Setting water to cup"
- Press multifunction button OK to enter the selected submenu which contains the following icons:

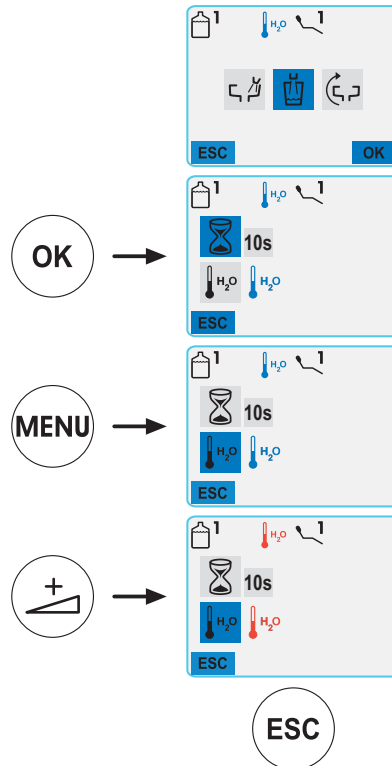
-  Cup filling level setting
-  Temperature of water to cup setting

- Press button MENU to go from one icon to another.
- Press buttons INCREASE and/or DECREASE to change the water delivery time to the cup or the temperature.

 **NOTE:** the temperature of the water to the cup can be set as follows:

-  COLD water (BLUE)
-  WARM water (ORANGE)
-  HOT water (RED)

- To confirm the settings, simply quit the submenu by pressing multifunction button ESC.

5.1.1.4.3. Automatic bowl motion settings

 **NOTE:** this menu can be accessed only if the powered bowl is present.


This setting is used for all the operators.

Perform the following steps from the HYDROGROUP SETTINGS menu:

- Press button MENU to scroll the icons, selecting the one for "Automatic bowl motion settings".
- Press multifunction button OK to enter the selected submenu that contains the following icons:

- | | |
|---|---|
|  | Bowl rotation controller with rinse position retrieval for the dental chair |
|  | Bowl rotation controller with automatic dental chair program retrieval |
|  | Bowl flushing controller with reset position retrieval for dental chair |

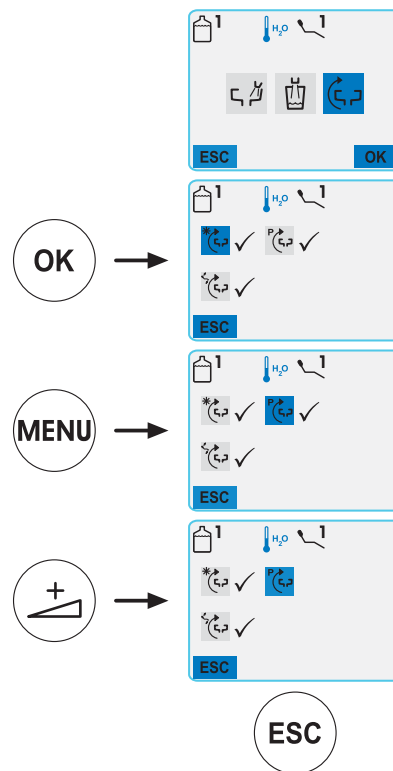- Press button MENU to go from one icon to another.
- Use buttons INCREASE and/or DECREASE to select/deselect a controller.
- To confirm the selections made, simply quit this submenu by pressing multifunction button

5.1.2. Setting the dental chair's "rinse" and "home" positions

This setting is specific for each operator.

Perform the following operations from the main screen:

- Bring the dental chair into the desired position with the manual movements buttons.

 **NOTE:** in the case of the "Rinse position" the position of the bowl can be saved, if powered.

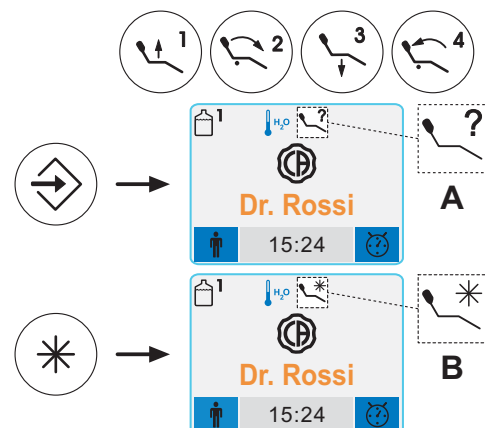
- Hold button "SAVE" for at least 2 seconds to activate memory mode. The system beeps and the relevant icon (A) appears on the console display to signal save mode has been activated.

 **NOTE:** Hold down button "MEMORIA" "SAVE" for at least 2 seconds to quit without saving the changes made.

- Press buttons "Automatic return" or "Rinse position" to assign the position to the button. The icon (B) for the selected program will appear on the console display to confirm the data have been saved.

 **NOTE:** The seat height cannot be changed when in the "Rinse position".

 **NOTE:** The "Rinse position" button brings the backrest and seat into the rinse position. When button "Rinse position" is pressed again, the backrest and seat return to the previous position.

5.1.3. Setting dental chair positions 1, 2, 3 and 4

This setting is specific for each operator.

Perform the following operations from the main screen:

- Bring the dental chair into the desired position with the manual movements buttons.

 **NOTE:** the position of the bowl can be saved, if powered. .

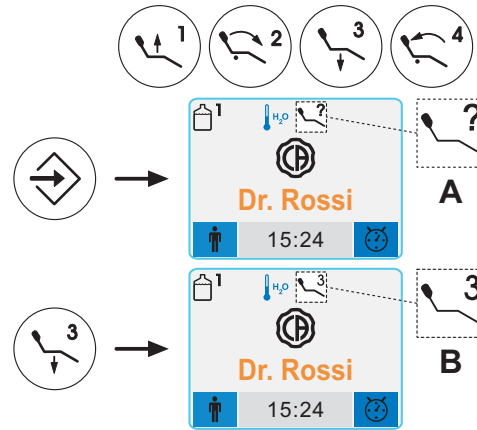
- Hold button "SAVE" for at least 2 seconds to activate memory mode.

 **NOTE:** A BEEP is given and an icon (A) appears on the top right-hand side of the screen when memory mode is activated.

- Press button "1" or "2" or "3" or "4" to assign the position to the button (e.g. 3).

 **NOTE:** The icon (B) that refers to the selected program (e.g. 3) appears on the display console to confirm the position has been saved.

 **NOTE:** To call up a set position simply briefly press the button assigned to the relative position.

5.1.4. Emergency stop button

This button can be used in the event of an emergency to bring the patient into the Trendelenburg position.

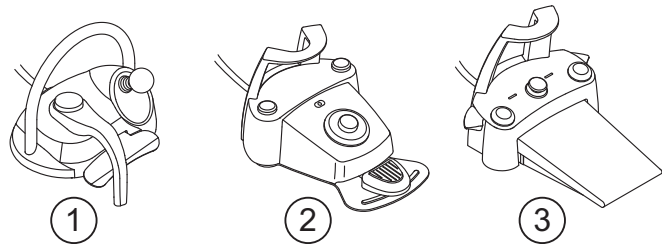
 **NOTE:** The Trendelenburg position is already set and cannot be changed.

5.2. Foot control

3 types of foot controls are available:

- 1 "Multifunction" foot control
- 2 "Push-pedal" foot control
- 3 "Power Pedal" foot control.

5.2.1. "Multifunction" foot control

Description of the parts

- [1] Handle
- [2] Horizontally actuated control lever
- [3] Joystick for dental chair movements
- [4] Vertically actuated control

Control lever [2]

When the instruments are in place, the control lever [2] has 2 different operating modes, depending on whether or not the pedal [4] is pushed before hand.

 **NOTE:** once the pedal [4] is pushed, the lever [2] must be actuated within 5 seconds. The function is automatically cancelled once this time has elapsed.

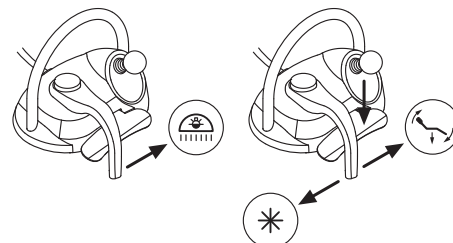
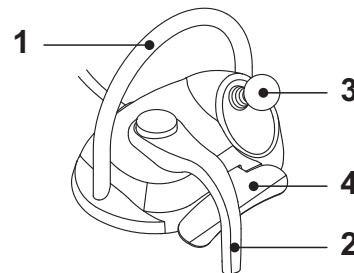
Operation with the instruments in place on the holder

- Lever [2] all the way to the right: dental operatory light on/off

Operation with the instruments in place on the holder after pushing the pedal [4]

- Lever [2] all the way to the right: dental chair automatic return
- Lever [2] all the way to the left: patient rinse position retrieval. If the lever [2] is pushed all the way to the left, the dental chair reaches the work position.

 **NOTE:** these dental chair functions are activated by holding down the pedal for at least 2 seconds.

Operation with the instruments withdrawn

- Lever [2] to the right: starts the instrument and adjusts its speed/power.

NOTE: if the camera is withdrawn, the image is frozen. .

- Lever [2] all the way to the left: CHIP-AIR for the turbine or the micromotor is started.

The air jet is stopped when the lever [2] is released.

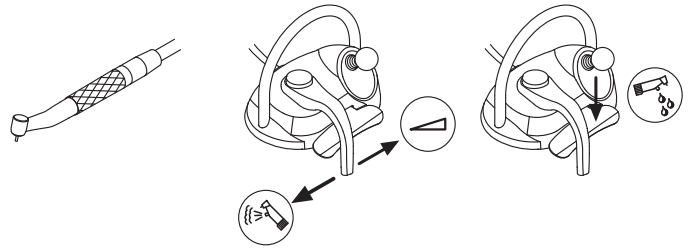
NOTE: if the camera is withdrawn, the saved images can be selected one after the other.

- Hold down the pedal [4] to start the CHIP-WATER for the turbine, micro-motor or scaler.

The water jet is stopped when the pedal [4] is released.

NOTE: if the camera is withdrawn, full-screen display of the selected image is activated/deactivated.

Water spray to instruments on/off

The spray is delivered by holding down the pedal [4] when the instrument is turned on (lever [2] activated).

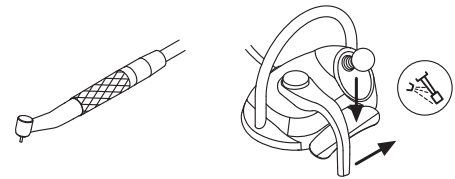
NOTE: the type of spray delivered can be set from the doctor's console.

NOTE: it is possible to set the foot control so that control [4] turns the water spray to the instrument on or off before actuating the lever [2] (bistable operation).

Only a CASTELLINI authorized technician can set the foot control in this manner.

NOTE: air is automatically blown at the end of work to eliminate any drops of liquid remaining in the spray ducts.

Joystick for dental chair movements [3]

The joystick for dental chair movements [3] has 2 different operating modes, depending on whether or not the pedal [4] is pushed before hand.

NOTE: once the pedal [4] is pushed, the joystick [3] must be actuated within 5 seconds. The function is automatically cancelled once this time has elapsed.

Operation with the instruments in place on the holder

The joystick [3] actuates the 4 manual dental chair movements.

NOTE: release the joystick to stop manual movement.

Operation with the instruments in place on the holder after pushing the pedal [4]

The joystick [3] retrieves the 4 dental chair programs.

NOTE: actuate the joystick again to stop automatic movement.

With instrument withdrawn but not turned on

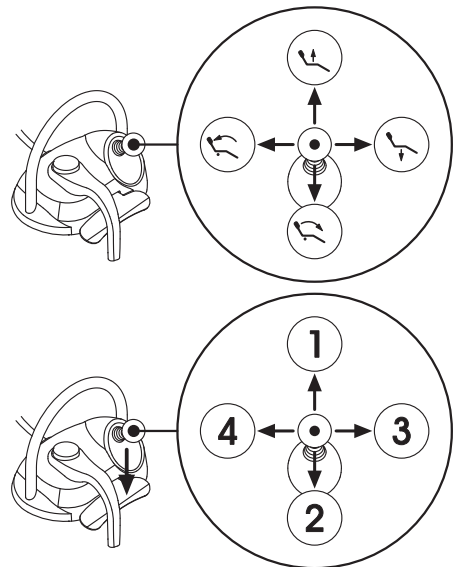
The joystick [3] actuates the 4 manual dental chair movements.

NOTE: the 4 programs are blocked under this condition.

With instrument turned on (lever [2] actuated)

Joystick disabled, all the commands to move the dental chair are blocked.

Foot control setup.

The foot control has 4 different configurations identified as 0, 1, 2 and 3:

0 = Joystick for only direction operation of the dental chair

1 = Auxiliary equipment and dental chair controllers (default configuration)

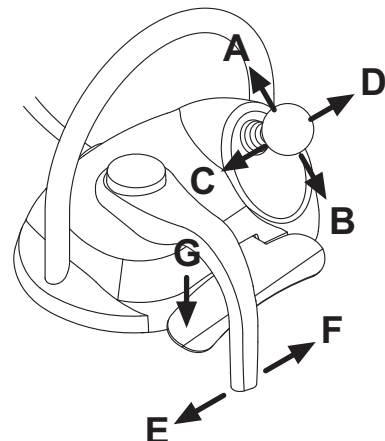
2 = ENDO mode and micromotor reverse rotation

3 = Surgical mode

Operation in the various configurations is explained below.

NOTE: the configuration of the foot control can be changed only by a CASTELLINI authorized technician.

SKEMA 6 - OPERATING INSTRUCTION

Dental chair movements (instruments in place)		CONFIGURATIONS			
ACTIVATION	FUNCTION	0	1	2	3
A	Up	X	X	X	X
B	Down	X	X	X	X
C	Back up	X	X	X	X
D	Back down	X	X	X	X

Automatic dental chair movements (instruments in place)		CONFIGURATIONS			
ACTIVATION	FUNCTION	0	1	2	3
G + A	Position 1 retrieval		X	X	X
G + B	Position 3 call up		X	X	X
G + C	Position 2 call up		X	X	X
G + D	Position 4 retrieval		X	X	X
G + E	Rinse position retrieval		X	X	X
G + F (End stop)	Reset position retrieval		X	X	X

Auxiliary equipment (instruments in place)		CONFIGURATIONS			
ACTIVATION	FUNCTION	0	1	2	3
F (End stop)	Dental operating light: on / off		X	X	X

Instruments (instruments withdrawn)		CONFIGURATIONS			
ACTIVATION	FUNCTION	0	1	2	3
F	Activation and/or speed/power adjustment	X	X	X	X
F e G	Startup with Spray	X	X	X	X
G	Scaler: start at max power Camera: image freeze	X	X	X	X
C	Micromotor: reverse rotation			X	X
D	Micromotor and Scaler: ENDO function			X	
A B D	Instrument menu navigation (surgical mode)				X
E	Activation with Chip-Air	X	X	X	X

Instrument menu navigation (surgical mode)		CONFIGURATIONS			
ACTIVATION	FUNCTION	0	1	2	3
A	INCREASE value (+)				X
B	DECREASE value (-)				X
C	Micromotor: reverse rotation				X
D	Micromotor: work programs M1, M2, M3 and M4 selection				X

Protection against liquid penetration.

The foot control is protected against liquid penetration.

Degree of protection: IPX1.

Cleaning

Clean the foot control with a suitable product approved against bacteria and fungi and effective against viruses (including HBV and HIV): see paragraph 1.4.

Pay strict attention to all the cleaning product manufacturer's instructions, warnings and cautions.

 **NOTE:** The product should be sprayed on soft disposable paper towel or sterile gauze.

 **ATTENTION!**

Do not use acids or harsh products.

 **NOTE:** If the foot control slips on the floor, dust the slip-proof rubber found under the base with a dry cloth.

5.2.2. "Push-pedal" foot control

Description of the parts.

- 1 Handle.
- 2 Control pedal.
- 3 Dental chair movements.
- 4 Water clean system/automatic dental chair return control.
- 5 Chip-air/patient rinse position control.
- 6 Operation with water spray LED.

Control pedal (2).

Operation:

- Remove the instrument
- Push the foot pedal to start the instrument (a).
- Adjust the rpm/power of the instrument with the control pedal:
 - to right: to increase;
 - to left: to decrease.

 **NOTE:** The control pedal adjusts the speed/power of the instrument from the minimum to the maximum value set from the instrument board.

- To stop the instrument, simply release the control pedal (a).

 **NOTE:** When water spray is used, air is automatically blown at the end of work to eliminate any drops of liquid remaining in the spray ducts.

ATTENTION!

Water spray to the instruments is turned on and off by holding down buttons (4) or (5) for at least 2 seconds. A beep sounds to signal the operating status has been changed.

LED (6) lights up to signal operation with water spray.

Dental chair buttons (3)

These buttons move the dental chair as follows:

Chair seat up

Chair back up

Chair seat down

Chair back down

To stop movement, release the button.

 **NOTE:** All the buttons used to move the dental chair are inoperative when an instrument is removed and the foot control pedal is actuated.

Water Clean System control (4).

Operation

- Key held down (at least 2 seconds) with the instrument removed:
Water Clean System control.

Running water is sent to the instruments such as the turbine, micromotor and scaler to flush the spray ducts. Water is delivered by pressing the button (4). Water is no longer delivered when the button is released and air is automatically blown to eliminate any drops of liquid remaining in the spray ducts.

- Briefly press key with instrument extracted:
Water spray to the instruments is turned on or off.

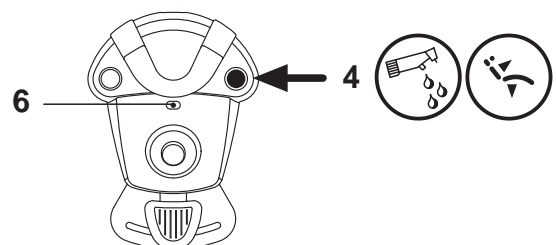
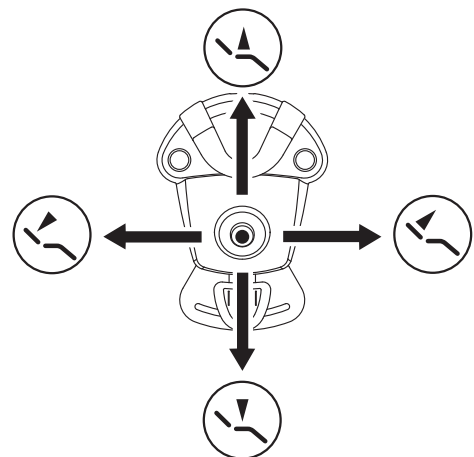
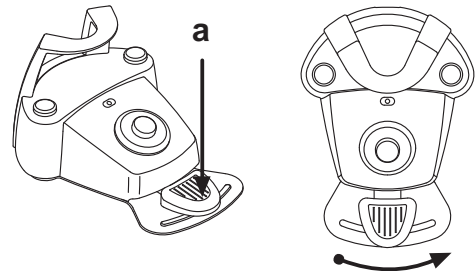
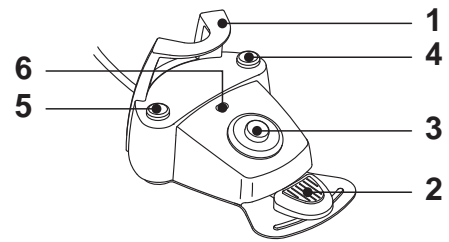
ATTENTION!

A beep sounds to signal the operating status has been changed.

LED (6) lights up to signal operation with water spray.

- Key held down (at least 2 seconds) with the instruments in place:
"Dental chair automatic return" function activated.

SKEMA 6 - OPERATING INSTRUCTION

Chip-air control (5)

Operation

- Key held down (at least 2 seconds) with the instrument removed:
Chip-air control.
Delivers air to the turbine or micromotor.
Air is delivered by pressing the button. Air is no longer blown when the button is released.
- Briefly press key with instrument extracted:
Water spray to the instruments is turned on or off.

ATTENTION!

A beep sounds to signal the operating status has been changed.
LED (6) lights up to signal operation with water spray.

- Key held down (at least 2 seconds) with the instruments in place:
"Dental chair automatic return" function activated.

NOTE: If button (5) is pressed and held down another time, the dental chair goes into its work position.

Protection against liquid penetration.

The foot control is protected against liquid penetration.

Degree of protection: IPX1.

Cleaning

Clean the foot control with a suitable product approved against bacteria and fungi and effective against viruses (including HBV and HIV): see paragraph 1.4.

NOTE: The product should be sprayed on soft disposable paper towel or sterile gauze.

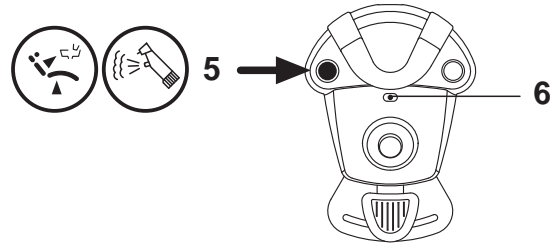
ATTENTION!

Do not use acids or harsh products.


NOTE: If the foot control slips on the floor, dust the slip-proof rubber found under the base with a dry cloth.

5.2.3. "Power Pedal" foot control

Description of the parts.

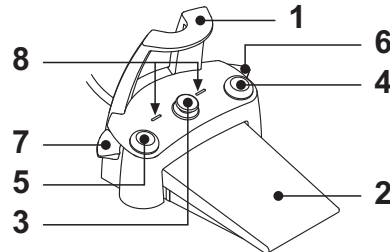- 1 Handle.
- 2 Foot control
- 3 Dental chair movements
- 4 Water Clean System control or spray to instruments function on/off
- 5 Chip-air control or spray to instruments function on/off
- 6 Program "B" retrieval or patient rinse position activation.
- 7 Program "A" retrieval or automatic dental chair return activation.
- 8 Operation with water spray LED

Foot control operation (2).

- With instrument removed
 - Push the foot control to start the instrument.
The instrument's rpm (or power) can be adjusted by varying the pressure exerted on the foot control.

NOTE: the foot control adjusts the speed/power of the instrument from the minimum to maximum value set from the instrument board.

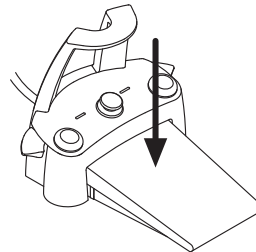
- Release the foot control to stop instrument operation.

NOTE: with water spray, air is automatically blown to eliminate any drops of liquid left in the spray ducts.

- With instrument in place

When the foot control is pressed, all automatic dental chair movements are automatically blocked.

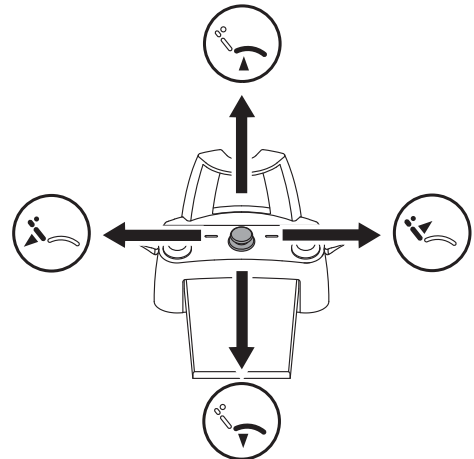
Operation of command for dental chair movements (3).

These buttons move the dental chair as follows:

- Chair seat up
- Chair back up
- Chair seat down
- Chair back down

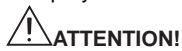
To stop movement, release the button.

 **NOTE:** all dental chair movements are blocked when an instrument is being used or the AUTOSTERIL system is running.


Right key (4) operation.

- Key held down (at least 2 seconds) with the instrument removed:
Water Clean System: sends a jet of running water to the instruments such as the Turbine, Micromotor and Scaler to flush the spray ducts. Water is delivered by pressing key (4); water is shut off when the key is released and air is automatically blown to eliminate any drops of liquid present in the spray ducts.
- Briefly press key with instrument extracted:
Spray water to instrument on or off.

ATTENTION!

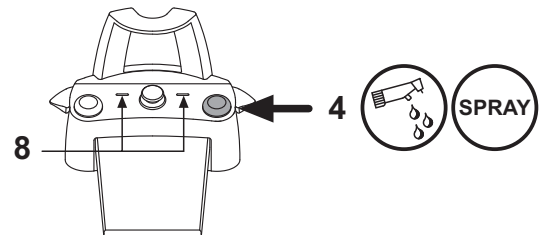A beep sounds to signal the operating status has been changed.
LED (8) lights up to signal operation with water spray.


Left key (5) operation.

- Key held down (at least 2 seconds) with the instrument removed:
Chip-air function: sends an air jet to the turbine or micromotor. Air is delivered by pressing the button. Air is no longer blown when the button is released.
- Briefly press key with instrument extracted:
Spray water to instrument on or off.

ATTENTION!

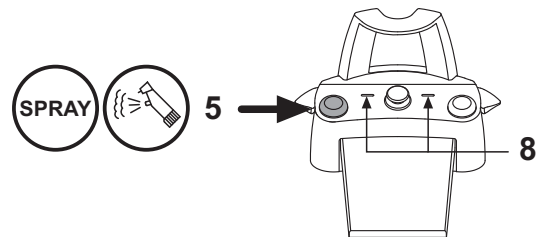A beep sounds to signal the operating status has been changed.
LED (8) lights up to signal operation with water spray.

Right lever operation (6).

 **NOTE:** the lever functions only with the instruments in their rest position.

For safety reasons, the selected function starts only after the switch has been briefly actuated and then released.

- Lever pushed down:
"Dental chair automatic return" program activated.
- Lever pulled up:
Dental chair program "2" start.

Left lever operation (7).

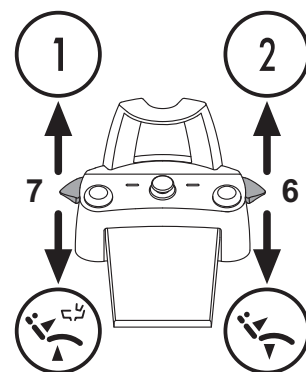
 **NOTE:** the lever functions only with the instruments in their rest position.

For safety reasons, the selected function starts only after the switch has been briefly actuated and then released.

- Lever pushed down:
"Rinse position" (PR) program activated

 **NOTE:** when the switch is actuated the second time, the dental chair reaches its work position.

- Lever pulled up:
Dental chair program "1" start.

Protection against liquid penetration.

The foot control is protected against liquid penetration.
Degree of protection: IPX1.

Cleaning

Clean the foot control with a suitable product approved against bacteria and fungi and effective against viruses (including HBV and HIV): see paragraph 1.4.

NOTE: The product should be sprayed on soft disposable paper towel or sterile gauze.

ATTENTION!

Do not use acids or harsh products.

NOTE: If the foot control slips on the floor, dust the slip-proof rubber found under the base with a dry cloth.

5.3. Syringe

Description of the instrument.

- [a] Nozzle.
- [b] Handpiece.
- [c] Syringe release button.
- [d] Air button.
- [e] Water button.
- [f] Hot/cold selector.
- [g] Hot/cold indicator light.

ATTENTION!

The instrument is supplied non-sterile.

Operation.

- **Operating times: work 5 sec., rest 10 sec.**
- Place the instrument in its work position.

NOTE: The relative icon is shown of the display screen to signal the instrument has been turned on.

- Operation with hot spray, air and water:
Turn the selector switch (f) clockwise (LED g on).
- Operation with cold spray, air and water:
Turn the selector switch (f) counter-clockwise (LED g off).
- Button (e) = water;
- Button (d) = air;
- Buttons (e + d) = spray.

Turning the fiber optic on/off.

- Set the instrument to its work position.

- Press the button to turn on the fiber optic.

NOTE: Icon (A) appears on the display to signal the fiber optic has been turned on.

- Adjust brightness with the “increase” (+) and/or “decrease” (-) buttons.
The set level is automatically saved.

Adjusting brightness of the fiber optic.

- Hold down the button (at least 2 seconds) to adjust brightness of the fiber optic .
- Use buttons INCREASE and/or DECREASE to adjust the brightness level.

NOTE: A value ranging from 1 to 16 can be set.

The following values can be set using the fast selection buttons:

- button “MAX”: level 16,
- button “MED”: level 8,
- button “MIN”: level 1.
- To confirm the selected brightness level, press multifunction button ESC to quit this submenu.

NOTE: the fiber optics automatically switch off after 30 seconds.

Removing the handpiece.

- The nozzle (a) is screwed onto the grip (b).
- Turn the selector switch counter-clockwise (LED g off) and press the button (c) to take the grip off the syringe casing.

Removable syringe cord

The syringe has a removable cord to ease cleaning (see chapter 5).

Cleaning.

Use soft disposable paper towel dampened with detergents/disinfectants.

ATTENTION!

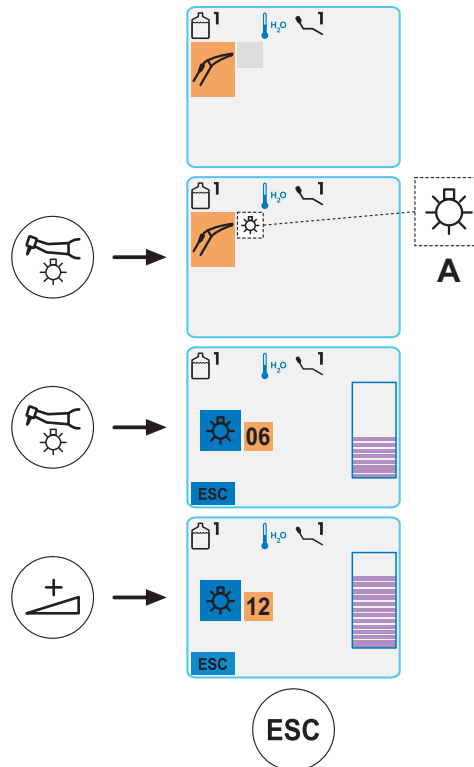
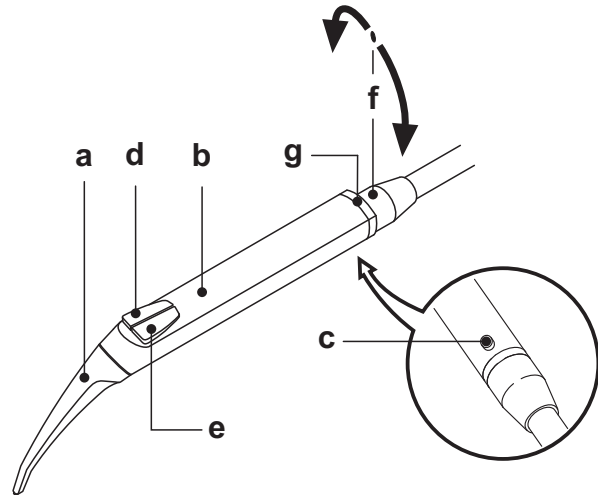
- **Do not soak the syringe in liquid disinfectants or detergents.**
- **Products not recommended: harsh products and/or products containing acetone, chlorine and sodium hypochlorites.**

Disinfection.

Syringe nozzle and handpiece: up to 135 °C (2 bar) in autoclave.
Minimum time: 15 min.

NOTE: Bag before sterilizing.

5.4. Turbine

Connecting the handpiece and changing the chuck.

Refer to the specific instructions furnished with the handpiece.

Use.

- **Operating time: work 5 min., rest 5 min.**
- Place the instrument in its work position.
- Use the foot control pedal to start the instrument (see paragraph 5.2).
- The cock (f) adjusts the water flow to the spray.
- The cock (e) adjusts the amount of air spray for all the instruments.

ATTENTION!

The instrument is supplied non-sterile.

Adjusting the turbine's operating speed.

- Place the instrument in its work position.
- Adjust the turbine's max. rpm with buttons INCREASE and/or DECREASE.

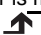
NOTE:

The following values can be set using the fast selection buttons:

- button "MAX": 100% power,
- button "MED": 50% power,
- button "MIN": 1% power.
- Use the right multifunction button to select the mode for varying the turbine's rpm.
The icon for the selected mode appears at the multifunction button.

a) Icon : Linear variation, proportional to the degree the foot control lever is moved.

b) Icon : ON/OFF mode which sees that the set maximum power is delivered when the foot control lever is actuated.

NOTE:

The set data are automatically saved.

Turning the fiber optic on/off.

- Set the instrument to its work position.
- Press the button  to turn on the fiber optic.

NOTE: Icon (A) appears on the display to signal the fiber optic has been turned on.

- Adjust brightness with the "increase" (+) and/or "decrease" (-) buttons.
The set level is automatically saved.

Adjusting brightness of the fiber optic.

- Hold down the button (at least 2 seconds) to adjust brightness of the fiber optic .
- Use buttons INCREASE and/or DECREASE to adjust the brightness level.

NOTE: A value ranging from 1 to 16 can be set.

The following values can be set using the fast selection buttons:

- button "MAX": level 16,
- button "MED": level 8,
- button "MIN": level 1.
- To confirm the selected brightness level, press multifunction button ESC to quit this submenu.

NOTE: The fiber optics shut off when the turbine is not used for 30 seconds (rheostat switch off).

Spray to instrument button.

It is possible to select the type of spray delivered to the instrument by pressing the relative button .

The relative icon (B) appears on the console display to signal the type of system operation.

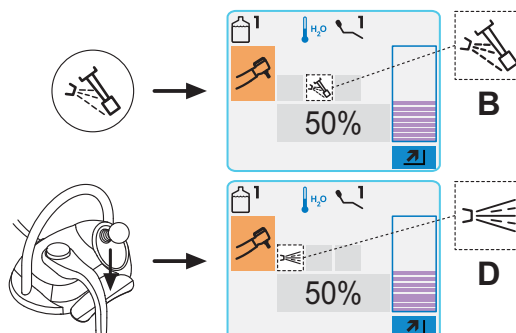
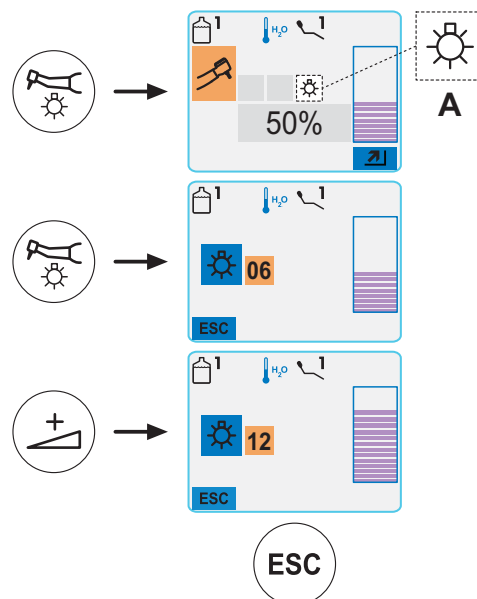
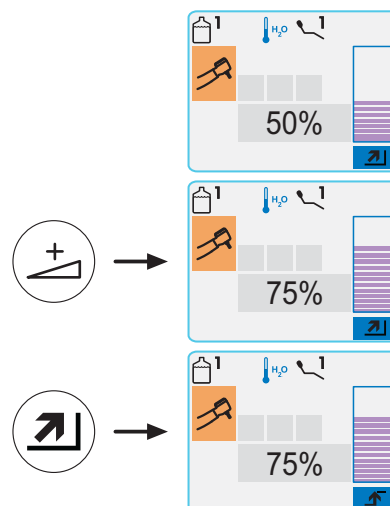
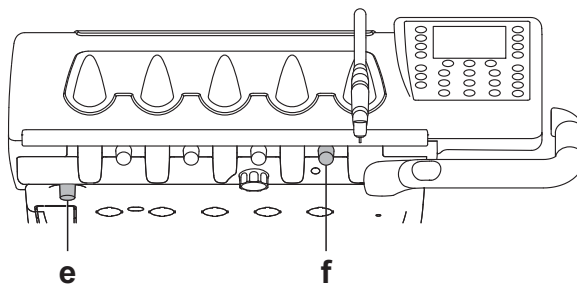
- Icon OFF: operation with water spray + air.
- Icon ON: operation with only water spray.

NOTE: if the foot control is set to bistable operation (see paragraph 5.2.1.), it is possible to set the instrument so that water spray is always delivered during operation.

In this case, the type of operation is indicated by the icon (D) on the console display:

- Icon OFF: operation without water spray.
- Icon ON: operation with water spray.

Removable cord

The turbine has a removable cord to ease cleaning (see paragraph 5.).

Cleaning and care.

Refer to the specific instructions furnished with the handpiece.

Disinfection.

Autoclave with steam at 135°C (2 bar). Minimum time: 15 min.

! ATTENTION!
Carefully read the operating instructions supplied with the handpiece before attempting to sterilize.

Safety guidelines.

! ATTENTION!

- The turbine must never be started without attaching the chuck or false chuck.
- The chuck release button must be held down during operation!
Friction between the button and micromotor rotor overheats the head and may cause burns.
- The patient's internal tissues (tongue, cheeks, lips, etc...) must be protected against contact with the button by using suitable instruments (mirror, etc...).
- The chucks and various instruments attached to the handpieces must comply to the standard ISO 10993-1 Biological evaluation of medical devices.

5.5. Micromotor

Coupling the handpieces and changing the chuck.

Refer to the specific instructions furnished with the micromotor and various handpieces.

Use.

- **Operating time: work 5 min., rest 5 min.**
- Place the instrument in its work position.
- Use the foot control pedal to start the instrument (see paragraph 5.2).
- The cock (f) adjusts the water flow to the spray.
- The cock (e) adjusts the amount of air spray for all the instruments.
- Press the left multifunction button (M1) with the instrument in its work position to select one of the 4 programs available. Each program stores the following data:
 - operating mode
 - torque value/maximum speed
 - fiber optic on/off
 - reverse rotation on/off
 - water spray on/off (only with "push-pedal" foot control)
 - peristaltic pump on/off (if present),
 - handpiece reduction ratio
- Press the MENU button with the instrument in its work position to select 3 different operating modes:

1) CONSERVATIVE (see paragraph 5.5.1)

- possibility to set a speed ranging from 100 to 50000 rpm regardless of the reduction/multiplication ratio.
- possibility to set torque values ranging from 1 to 100%
- customized list of reduction ratios
- rpm variation mode can be changed over from variable to fixed and vice versa
- alarm signaled when maximum torque is reached
- maximum speed quickly reached during motor operation

2) ENDO (see paragraph 5.5.2)

- possibility to set a speed ranging from 100 to 600 rpm with value always in reference to the burrs regardless of the reduction ratio
- possibility to set torque values ranging from 0.1 to 5.0 Ncm, excluding the gearbox 1:1 (4.5 Ncm)
- customized list of reduction ratios
- rpm variation mode can be changed over from variable to fixed and vice versa
- progressive alarm signal starting from 60% of the maximum torque
- button for calibration during motor operation

3) SURGERY (see paragraph 5.5.3)

- possibility to set a speed ranging from 5 to 2500 rpm with value always in reference to the burrs regardless of the reduction ratio (1:20 to 1:1000 gearboxes)
- possibility to set torque values ranging from 0.5 to 55.0 Ncm for certified gearboxes or 1 to 100%
- customized list of reduction ratios
- alarm signaled when maximum torque is reached
- button for calibration during motor operation

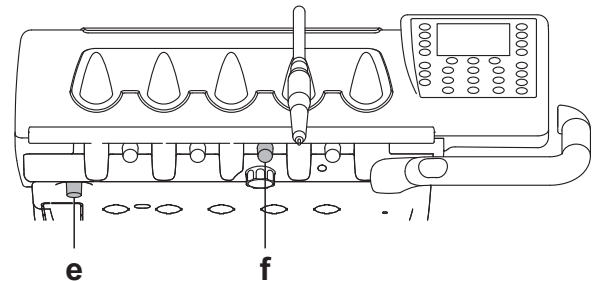
- Press the right multifunction button (M2) with the instrument in its work position to go to the setup submenu regarding the operating mode currently used.

- The following values are displayed:
 - torque value (set or actual) expressed in % or Ncm for certified gearboxes.
 - speed of operation of burrs (set or actual) expressed in rpm.

NOTE : the bar graph always shows the torque value expressed as a percentage in relation to the set/maximum allowable one.




! ATTENTION!

The instrument is supplied non-sterile.

MENU

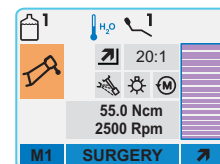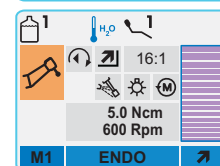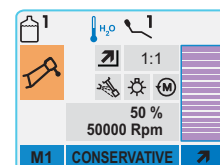MENU

**Turning the fiber optic on/off.**

- Set the instrument to its work position.
- Press the button  to turn on the fiber optic.

 **NOTE:** Icon (A) appears on the display to signal the fiber optic has been turned on.

- Adjust brightness with the “increase” (+) and/or “decrease” (-) buttons. The set level is automatically saved.

Adjusting brightness of the fiber optic.

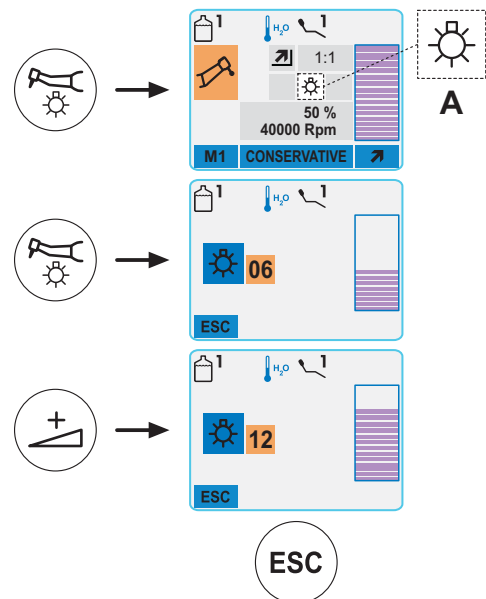
- Hold down the button (at least 2 seconds) to adjust brightness of the fiber optic .
- Use buttons INCREASE and/or DECREASE to adjust the brightness level.

 **NOTE:** A value ranging from 1 to 16 can be set. The following values can be set using the fast selection buttons:

- button “MAX”: level 16,
- button “MED”: level 8,
- button “MIN”: level 1.

- To confirm the selected brightness level, press multifunction button ESC to quit this submenu.

 **NOTE:** The fiber optics shut off when the turbine is not used for 30 seconds (rheostat switch off).

**Spray to instrument button**

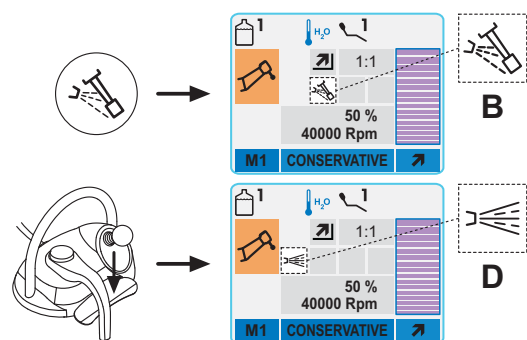
Use the foot control (see paragraph 5.2) to start the water spray to the instrument. The type of spray delivered to the instrument can be selected with the relative button . The type of system operation is signaled by the relative icon (B) shown on the display console.

- Icon OFF: operation with spray water and air
- Icon ON: operation with spray water only

 **NOTE:** if the foot control is set to bistable operation (see paragraph 5.2.1.), it is possible to set the instrument so that water spray is always delivered during operation.

In this case, the type of operation is indicated by the icon (D) on the console display:

- Icon OFF: operation without water spray.
- Icon ON: operation with water spray.

**Reversing the micromotor.**

Press REVERSE button  to reverse the direction in which the micromotor runs.


A beep is heard and the icon (C) appears on the display to signal the direction of rotation has been reversed.

**ATTENTION!**

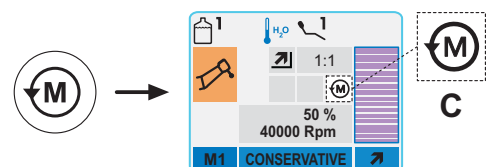The beeper BEEPS to signal the direction is reversed after the micromotor has been removed.

 **NOTE:**

When the rheostat's lever is actuated, the micromotor cannot reverse.

**Removable cord**

The turbine has a removable cord to ease cleaning (see paragraph 5.).

Cleaning and care.

Refer to the specific instructions furnished with the instrument.

**ATTENTION!**

- Do not soak the syringe in liquid disinfectants or detergents.
- Products not recommended: harsh products and/or products containing acetone, chlorine and sodium hypochlorites.

Disinfection.

Handpieces only: autoclave with steam at 135 °C (2 bar). Minimum time: 15 min.

**ATTENTION!**

Carefully read the operating instructions supplied with the instrument before attempting to sterilize.

Safety guidelines.**ATTENTION!**

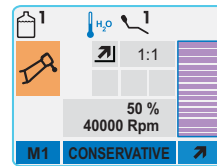
- Never put the contra angle on the micromotor while it is running.
- The chuck release button must be held down during operation! Friction between the button and micromotor rotor overheats the head and may cause burns.
- The patient's internal tissues (tongue, cheeks, lips, etc...) must be protected against contact with the button by using suitable instruments (mirror, etc...).
- The chucks and various instruments attached to the handpieces must comply to the standard ISO 10993-1 Biological evaluation of medical devices.

5.5.1 CONSERVATIVE operating mode

Menu with micromotor withdrawn but not turned on.

- Press the left multifunction button (M1) to select one of the 4 programs available.
- Press the MENU button to change the operating mode:
- Press the right multifunction button (↗) to go to the submenu used to edit the operating data.
- Use buttons INCREASE and/or DECREASE to change the maximum speed of the burrs.
The new value entered is automatically saved in the selected work program (e.g. M1).

 **NOTE:** : Use the fast selection buttons to automatically set the following values:


- MAX button: 100% of the maximum allowable speed
- MED button: 50% of the maximum allowable speed
- MIN button: 1% of the maximum allowable speed

Menu with micromotor withdrawn and turned on.

- Press the left multifunction button (▶|◀◀) to set the current speed to the maximum allowable one simultaneously selecting the ON/OFF variation mode for the foot control switch.

 **NOTE:** : in order to quit this data item, go to the submenu used to set the operating data and change the foot control switch over from ON/OFF to linear speed variation mode.

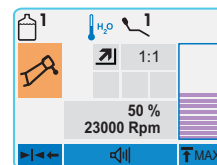- Press the MENU button to activate/deactivate an alarm signal once the set maximum torque is reached.
The icon indicating the selected mode appears near the MENU button.
 - icon : alarm active
 - icon : alarm not active
- Press the right multifunction button (↑ MAX) to set the current speed to the maximum allowable one.
- Use buttons INCREASE and/or DECREASE to change the maximum speed of the burrs.

Set up menu

Press the right multifunction button (↗) with the micromotor withdrawn but not turned on to go to the submenu containing the following icons which indicate the settable functions:

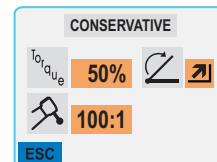
Reduction ratio setting

Maximum torque variation setting

Speed variation mode

- Press the MENU button to scroll the icons until the desired function is reached.
- Setting the reduction ratio
See paragraph 5.5.4.
- Setting maximum torque variation
Use buttons INCREASE and/or DECREASE to change the maximum torque value that can be reached.

 **NOTE:** Use the fast selection buttons to automatically set the following values:

- MAX button: 100% of the maximum allowable speed
- MED button: 50% of the maximum allowable speed
- MIN button: 1% of the maximum allowable speed

Speed variation mode

- Use buttons INCREASE and/or DECREASE to select the speed variation mode:
 - icon : linear variation, proportional to movement of the foot control switch
 - icon : ON/OFF variation. The set maximum power is delivered when the foot control switch is actuated.
- To confirm the selections made, simply quit this submenu by pressing multifunction button ESC.

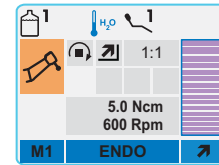
5.5.2 ENDO operating mode

Menu with micromotor withdrawn but not turned on.

- Press the left multifunction button (M1) to select one of the 4 programs available.
- Press the MENU button to change the operating mode:
- Press the right multifunction button (↗) to go to the submenu used to edit the operating data.
- Use buttons INCREASE and/or DECREASE to change the maximum allowable torque value..
The new value entered is automatically saved in the selected work program (e.g. M1).

 **NOTE:** Use the fast selection buttons to automatically set the following values:

- MAX button: 100% of the maximum allowable torque value
- MED button: 50% of the maximum allowable torque value
- MIN button: 1% of the maximum allowable torque value

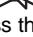
Menu with micromotor withdrawn and turned on.


- Press the left multifunction button (.....) to set the current torque to the maximum allowable one simultaneously selecting the ON/OFF variation mode for the foot control switch.

 **NOTE:** in order to quit this data item, go to the submenu used to set the operating data and change the foot control switch over from ON/OFF to linear variation mode.

- Press the MENU button to activate/deactivate a progressive alarm signal starting from 60% of the set maximum torque.

The icon indicating the selected mode appears near the MENU button.

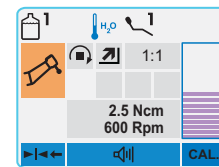- icon  : alarm active
- icon  : alarm not active
- Press the right multifunction button (CAL..). While the handpiece is idling at maximum, the current torque value is set to value 0 (handpiece calibration).

Set up menu

Press the right multifunction button (↗) with the micromotor withdrawn but not turned on, to go to the submenu containing the following icons which indicate the settable functions:

Reduction ratio setting

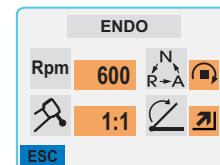
Maximum speed setting

Operation with maximum torque reached

Speed variation mode

- Press the MENU button to scroll the icons until the desired function is reached.
- Setting the reduction ratio
See paragraph 5.5.4.
- Setting the maximum speed
Use buttons INCREASE and/or DECREASE to change the maximum speed that can be reached.

 **NOTE:** Use the fast selection buttons to automatically set the following values:

- MAX button: 600 rpm
- MED button: 300 rpm
- MIN button: 100 rpm

- Setting operation with maximum torque reached

Use buttons INCREASE and/or DECREASE to change how the micromotor reacts when the set maximum torque is reached:

Micromotor stops operating

Micromotor stops operating and then reverses direction of rotation

Micromotor stops operating, reverses the direction of rotation and then goes back to the normal direction of rotation.

- Speed variation mode

Use buttons INCREASE and/or DECREASE to select the speed variation mode:

- icon  : linear variation, proportional to movement of the foot control switch
- icon  : ON/OFF variation. The set maximum power is delivered when the foot control switch is actuated.

- To confirm the selections made, simply quit this submenu by pressing multifunction button ESC.

5.5.3 SURGERY operating mode

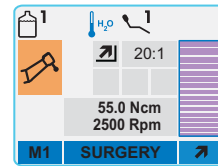
Menu with micromotor withdrawn but not turned on.

- Press the left multifunction button (M1) to select one of the 4 programs available.
- Press the MENU button to change the operating mode:
- Press the right multifunction button (↗) to go to the submenu used to edit the operating data.
- Use the increase and/or decrease buttons to modify the torque of the burs.

The new value entered is automatically saved in the selected work program (e.g. M1).

NOTE: Use the fast selection buttons to automatically set the following values:

- MAX button: 100% of the maximum settable torque,
- MED button: 50% of the maximum settable torque,
- MIN button: 1% of the maximum settable torque.

Menu with micromotor withdrawn and turned on.

- Press the left multifunction button (↔) to set the current speed to the maximum allowable one simultaneously selecting the ON/OFF variation mode for the foot control switch.

NOTE: in order to quit this data item, go to the submenu used to set the operating data and change the foot control switch over from ON/OFF to linear speed variation mode.

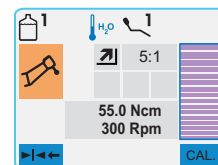
- Press the MENU button to activate/deactivate an alarm signal once the set maximum torque is reached.

The icon indicating the selected mode appears near the MENU button.

- icon : alarm active
- icon : alarm not active

- Press the right multifunction button (CAL.). While the handpiece is idling at maximum, the current torque value is set to value 0 (handpiece calibration).

- Use buttons INCREASE and/or DECREASE to change the maximum speed of the burs.

Set up menu

Press the right multifunction button (↗) with the micromotor withdrawn but not turned on, to go to the submenu containing the following icons which indicate the settable functions:

Reduction ratio setting

Maximum speed setting

Speed variation mode

- Press the MENU button to scroll the icons until the desired function is reached.

- Setting the reduction ratio
See paragraph 5.5.4.

- Maximum speed setting

Use the increase and/or decrease buttons to modify the highest speed that can be reached.

NOTE: use the fast selection buttons to automatically set the following values:

- MAX button: 100% of the maximum settable speed,
- MED button: 50% of the maximum settable speed,
- MIN button: 1% of the maximum settable speed.

- Speed variation mode

Use buttons INCREASE and/or DECREASE to select the speed variation mode:

- icon : linear variation, proportional to movement of the foot control switch
- icon : ON/OFF variation. The set maximum power is delivered when the foot control switch is actuated.

- To confirm the selections made, simply quit this submenu by pressing multifunction button ESC.

5.5.4 Reduction ratio setup menu

Select the icon for "Reduction ratio setting" on the menu used to set the operating data and press the right multifunction button (↗) to go to the submenu containing the following icons:

Current reduction ratio

Maximum speed that can be reached with selected reduction ratio

NOTE: icon Rpm cannot be edited as it only shows the maximum speed that can be reached with the reduction ration selected above.

How to change the current reduction ratio:

- Use buttons INCREASE and/or DECREASE to scroll the saved reduction ratios.

NOTE: the fast selection buttons can be used to set the preset values that vary according to the current operating mode.

- To confirm the selections made, simply quit this submenu by pressing multifunction button ESC.

How to create customized reduction ratios:

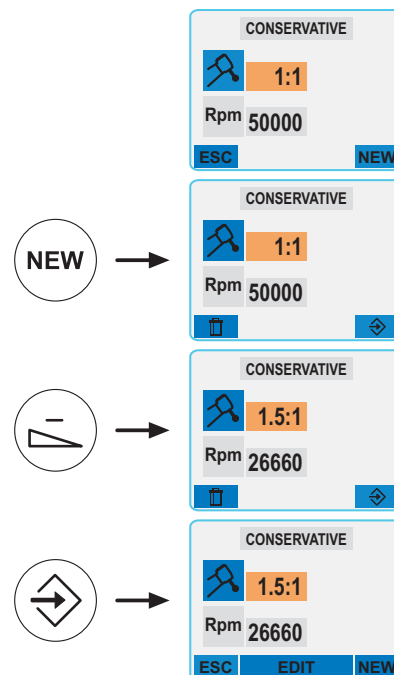
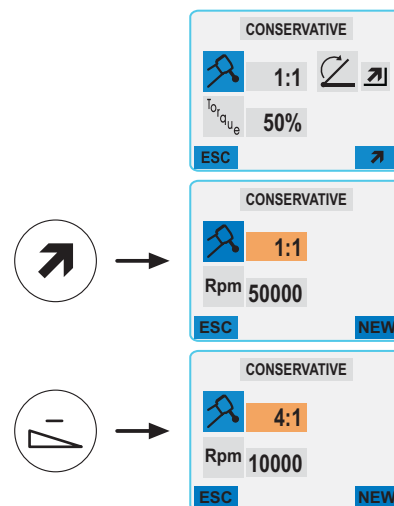
Follow the directions given below to create and save customized reduction ratios.

- Press the right multifunction button (NEW).
- Press buttons INCREASE and/or DECREASE to modify the tenths.
- Press the MAX and/or MIN fast selection buttons to change the units.
- Once the desired reduction ratio has been set, press the right multifunction button (↗) to save it.

How to edit and/or delete customized reduction ratios:

NOTE: Only customized reduction ratios can be edited and/or deleted.

- Press buttons INCREASE and/or DECREASE to scroll the saved reduction ratios.
- Press the MENU button (EDIT) to change over to editing mode.
- Press buttons INCREASE and/or DECREASE to modify the tenths.
- Press the MAX and/or MIN fast selection buttons to change the units.
- Press the right multifunction button (↗) to save it again.
- Press the left multifunction button (↖) to delete.

5.6. Scaler

Connecting the handpiece and inserts.

Refer to the specific instructions furnished with the handpiece.

ATTENTION!

Before attempting to connect the handpiece, make certain the contacts are perfectly dry. Blow air from the syringe, if necessary, to dry.

Use.

- **Operating time:**
 - PIEZOSTERIL 5 scalers: continuous operation,
 - PIEZOSTERIL 6 scalers: continuous operation,
 - PIEZOLIGHT 6 scalers: continuous operation.
 - SURGISON 2: scalers: continuous operation.
- Place the instrument in its work position.
- Use the foot control pedal to start the instrument (see paragraph 5.2).
- The cock [f] adjusts the cooling water flow.

ATTENZIONE!

The instrument is supplied non-sterile.

Safety guidelines

ATTENTION!

- Make sure the threaded sections of the inserts and handpiece are perfectly clean.
- Do not change the shape of the inserts.
- Check wear and tear of the inserts on a regular basis, replacing them in the following cases:
 - obvious wear
 - drop in performance
 - out of shape or banged.
- **Note for PIEZOLIGHT 6 scalers:**
 - Class 1 LED apparatus;
 - Do not direct the light beam in anyone's eyes When cleaning or servicing the device (it is recommended to keep the fiber optics shut off).

Adjusting the scaler's power.

- Put the instrument in its work position.
- Use the buttons INCREASE and/or DECREASE. to adjust the maximum power delivered.

NOTE: The following values can be set using the fast selection buttons:





- button "MAX": 100% power,
- button "MED": 50% power,
- button "MIN": 25% power,
- Use the right multifunction button to select the mode used to vary the scaler's power.

The icon indicating the selected mode appears on the display beside the set maximum power percentage.

- Icon : Linear variation, proportional to the degree the foot control lever is moved.
- Icon : ON/OFF mode which sees that the set maximum power is delivered when the foot control lever is actuated.

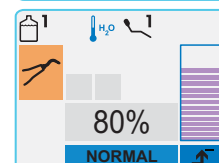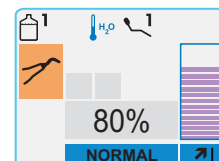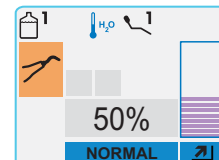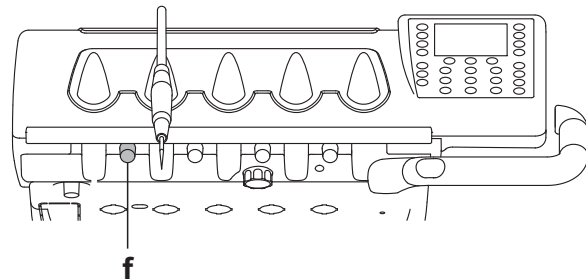NOTE: The set data are automatically saved



ENDO function.

The scaler functions at 1/2 the set maximum power.

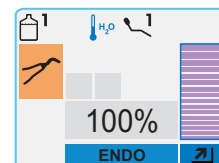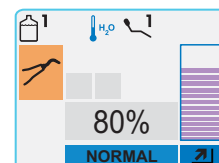- Put the instrument in its work position.
- Press the MENU button to activate the ENDO function. ENDO will appear in place of NORMAL on the display to confirm the mode of operation has been changed.
- Press the MENU button again to go back to normal operation.

NOTE: the ENDO function is deactivated when the foot control switch is actuated.

Turning the fiber optic on/off.

- Set the instrument to its work position.
- Press the button  to turn on the fiber optic.

 **NOTE:** Icon (A) appears on the display to signal the fiber optic has been turned on. If the handpiece used does not have the fiber optic, icon (E) will appear on the display.

Adjusting brightness of the fiber optic.

- Hold down the button (at least 2 seconds) to adjust brightness of the fiber optic .
- Use buttons INCREASE and/or DECREASE to adjust the brightness level.

 **NOTE:** A value ranging from 1 to 16 can be set.
The following values can be set using the fast selection buttons:

- button "MAX": level 16,
- button "MED": level 8,
- button "MIN": level 1.

- To confirm the selected brightness level, press multifunction button ESC to quit this submenu.

 **NOTE:** The fiber optics shut off when the turbine is not used for 30 seconds (rheostat switch off).

Instrument cooling water.

If the foot control is set to bistable operation (see paragraph 5.2.1.), it is possible to set the instrument so that cooling water is always delivered during operation.

In this case, the type of operation is indicated by the icon (D) on the console display:

- Icon OFF: operation without cooling water.
- Icon ON: operation with cooling water.

Removable cord

The scaler has a removable cord to ease cleaning (see paragraph 5.).

Cleaning and care.

Refer to the specific instructions furnished with the instrument.

ATTENTION!

- Do not soak the handpiece in liquid disinfectants or detergents.
- The torque wrench should not be sterilized, cold disinfect with isopropyl alcohol.

Sterilization.

- Inserts: Scaler tips: autoclave at 135°C (2 bar), minimum time: 15 min.
- Handpiece: autoclave at 135°C (2 bar), minimum time: 15 min.

ATTENTION!

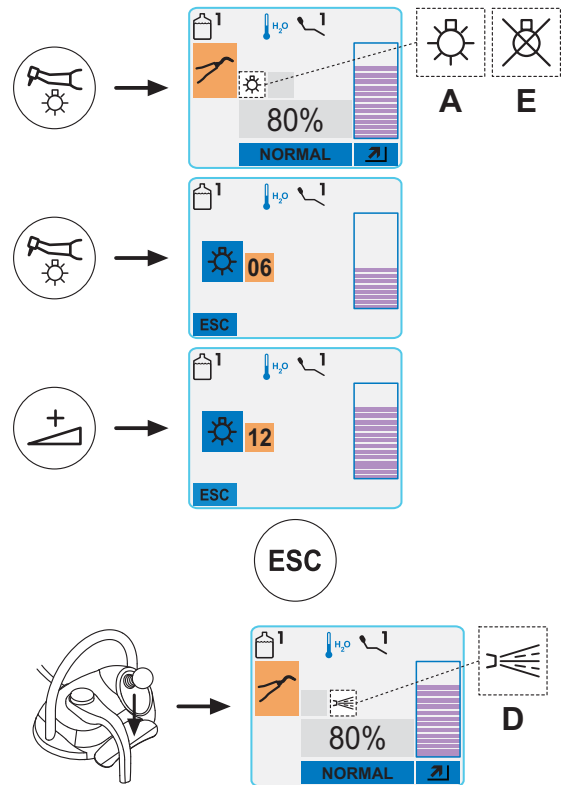
Carefully read the operating instructions supplied with the instrument before attempting to sterilize.

Safety standards.

ATTENTION!

- To avoid hazards or malfunctions When connecting the board, do not reverse the positions of the cords for scalers that are different brands.
- The inserts attached to the handpiece must comply to Biocompatibility standard ISO 10993-1.

5.6.1. SURGISON 2 surgical detartarer device

The Surgison 2 handpiece is intended exclusively for surgical bone operations in the odontoiatric field. Any other use can cause serious injury to the patient and the damage or breakage of the instrument.

Usage instructions.

CAUTION!

- Ordinary water from the mains must **NOT** be used for cooling the operational field. Make sure the cooling liquid used is of a suitable sterile type.
- The SURGISON 2 handpiece must only be used with the bits supplied with the instrument or other bits supplied by CASTELLINI S.p.A.
- Do not change the shape of the bit
- The instrument is supplied in a non-sterile state.

Use.

- **Operational times:** continuous operation.
- Put the instrument in the working position. The system automatically detects the presence of a SURGISON 2 handpiece and displays the relative menu.
- The maximum power generated can be changed by pressing the INCREASE and/or DECREASE buttons

NOTE: with the rapid selection buttons it is possible to automatically set the following values:

- MAX button: 100% power
- MED button: 50% power
- MIN button: 25% power

- With the instrument in the working position, it is possible to choose two different operational modes by pressing the MENU button:

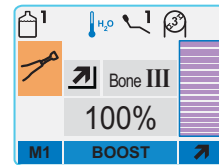
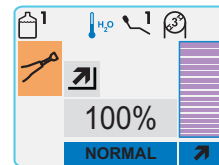
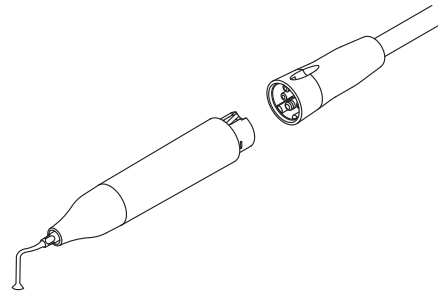
1) NORMAL (see paragraph 5.6.1.1.):

- The bit vibrates subjected to a continuous (not modulated) power the amplitude of which depends on the set value of the power.

2) BOOST (see paragraph 5.6.1.2):


- in addition to the ultrasonic vibration of the bit, a modulation in the amplitude of the generated power (approx. +40% of the set power) is performed. The modulation creates a percussion effect on the bone thereby making incision easier.
- three different modulation times can be set depending on the thickness of the bone tissue to be treated.
- There are 3 separate working programs.

- By pressing the  button for a longer time (at least two seconds), it is possible to modify the amount of physiological solution supplied by the peristaltic pump (see paragraph 5.10).
- To start up the instrument press the foot control (see paragraph 5.2).


5.6.1.1. Normal operational mode

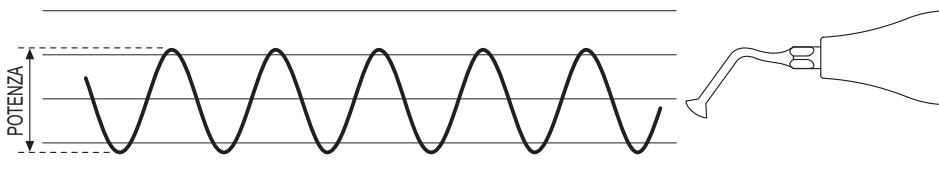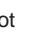The bit vibrates subject to a continuous power, the amplitude of which depends on the power value set.

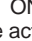
Menu with surgical detartarer device extracted but not active.

- When the right-hand multifunction button () is pressed, the submenu is accessed to change the following operational settings:



Power variation mode

Press the INCREASE and/or DECREASE buttons to select the power variation mode:

- icon : linear variation, proportional to the movement of the foot control lever,
- icon : ON/OFF variation that entails the supply of maximum power set on the activation of the foot control lever.

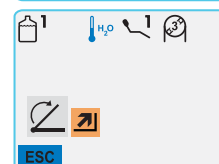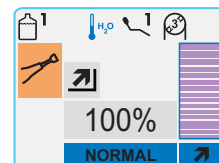
To confirm the settings selected you just need to leave this sub menu pressing the ESC multifunction button.

5.6.1.2. BOOST functioning mode

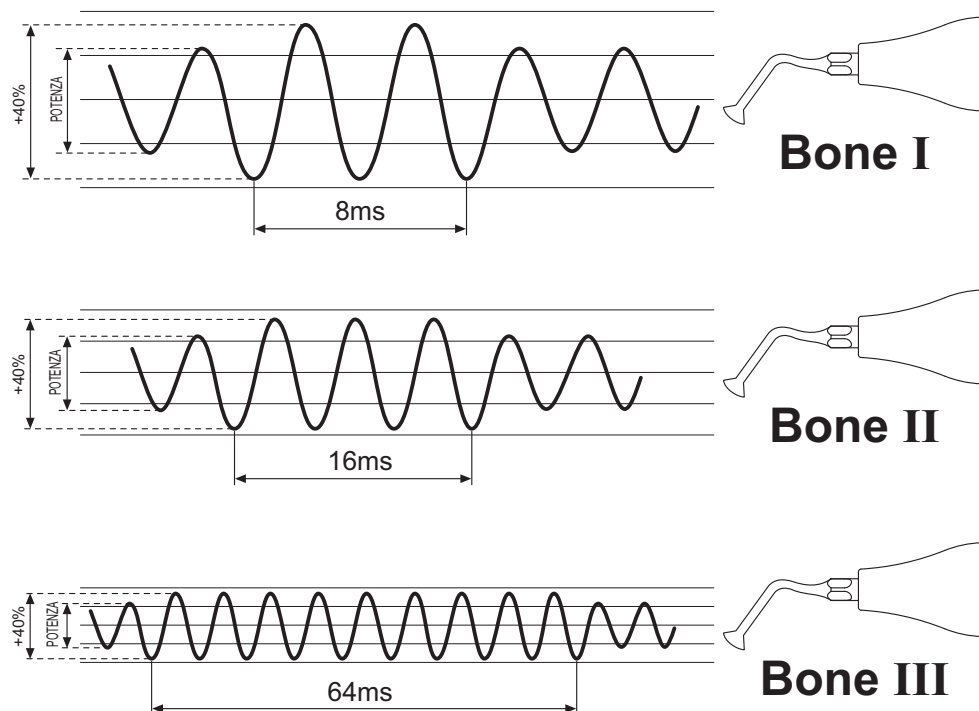In addition to the ultrasonic vibration of the bit, a modulation with the amplitude of the generated power (approx. +40% of the set power) is performed. The modulation creates a percussion effect on the bone thereby making incision easier. Three different modulation times can be set depending on the thickness of the bone tissue to be treated.

Menu with surgical detaratarer device extracted but not active.

- When the left-hand multifunction button (M1) is pressed, one of the three available working programmes is selected. Each working programme stores the following data:
 - maximum power generated,
 - the quantity of physiological solution supplied by the peristaltic pump,
 - set modulation (Bone)
 - the power variation mode during functioning.
- By pressing the right-hand multifunction button (↗) the submenu for changing the operational settings is accessed.

Operational settings menu.

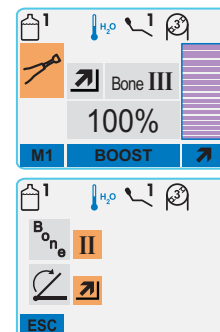
With the instrument extracted but not active, when the right-hand multifunction button (↗) the submenu is accessed where the following icons regarding the settable functions are to be found:

For setting the type of modulation

Method for varying the supplied power

- Press the MENU button to scroll through the icons until you select the function required.
- Setting the type of modulation.
Press the INCREASE and/or DECREASE buttons to modify the modulation depending on the thickness of the bone tissue to be treated:
 - Bone I : hard bone tissue (8ms);
 - Bone II : Medium bone tissue (16ms);
 - Bone III : soft bone tissue (64ms).
- Supplied power variation mode.
Press the INCREASE and/or DECREASE buttons to select the rotation speed variation mode:
 - icon ↗: linear variation, proportional to the movement of the foot control lever,
 - icon ↗: ON/OFF variation that entails the supply of maximum power set on the activation of the foot control lever.
- To confirm the settings selected you just need to leave this submenu pressing the ESC multifunction button.

Safety guidelines

ATTENTION!

- The LED is a class 2 source ($P_o < 1\text{mW}$) in compliance with international standard IEC 60825-1. NEVER CONTINUOUSLY STARE INTO THE LASER SOURCE.

The light emitted may cause eye damage in the event of direct radiation without eye protection.

Always use the light with the protective screen (provided) or eyewear that protects against radiation emitted by the blue light up to a wave length of 520 nm.

Do not aim the light beam in anyone's eyes.

- **The light emitted may damage soft tissues (oral cavity mucous, gums, skin).**
Be extremely careful to direct the light precisely on the material to be cured.
- **Do not use the device on patients taking photosensitising drugs.**
- **People with eye diseases, such as those who have had cataracts removed or retina diseases must be adequately protected When the curing lamp is used, for example with suitable protective eyewear.**
- **Do not pull the power cord.**
- **Do not expose the handpiece to excessive vibrations.**
- **Do not drop the handpiece and in particular the fiber optic. The lamp may break if accidentally banged.**
Check the condition of the handpiece if it has been banged or dropped before using the curing light. Try to turn on the light and check operation first without using it on the patient.

If cracked, broken or if there are any other faults, do not use the curing light on the patient and contact the technical service department.

The fiber optic is rather delicate and may crack or break if banged, affecting the final amount of light emitted. If dropped, carefully inspect the fiber optic to verify if it is cracked or broken. If cracked, a strong light appears in the spot in which the fiber is cracked. In all these cases, the fiber optic must be replaced.

- **The curing light handpiece (sold separately) can be connected only to dental units with connections for this curing lamp. Connection to any other equipment may damage the circuits inside the lamp and seriously injure the user and patient.**
- **The curing lamp handpiece is not protected against liquid penetration (IP20).**
While cleaning and disinfecting, be careful not to let the detergent get inside the device.
- **The curing lamp handpiece is not suitable for use in the presence of a mix of inflammable anaesthetic gas with oxygen or nitrous oxide.**
- **The equipment is designed for intermittent operation: (60 sec work – 13 min. pause).**

Cleaning

The curing lamp may be a vehicle for cross contamination between patients.

The most contaminated parts are the fiber optic and eye protection.

Only the fiber optic is autoclavable with water steam at 135 °C. Always clean the fiber optic before sterilizing.

ATTENTION!

- **The fiber optic is able to support 500 autoclave cycles after which it tends to become dull and may therefore emit less light. Contact the manufacturer to purchase original spare parts**
- **Periodically check the autoclave as directed by the manufacturer; temperatures that exceed the indicated limit may damage the fiber optic.**

The handpiece and eye protection cannot be put in autoclave; disinfect them with suitable products and cover the handpiece with disposable plastic wrap.

Use soft disposable paper towels to disinfect the handpiece and eye protection. Do not use harsh products or soak in liquids.

ATTENTION!

- **The curing light handpiece is NOT suitable for autoclave.**
- **The curing light handpiece is not protected against penetration of liquids therefore it CANNOT be soaked in solution to be sterilized.**
- **The outside of the lamp should be disinfected with the fiber optic on. Do not use any type of disinfectant on the exposed optical surface of the handpiece When the fiber is removed. The surface will become irreparably opaque if it comes into contact with disinfectant.**

Maintenance

This equipment does not require any particular type of maintenance. Only technicians authorized by the manufacturer can replace and/or repair the handpiece and dental unit. The handpiece has been purposely constructed in a manner that requires specific tools to open it and therefore it cannot be removed by the user. The warranty is automatically void if the handpiece is altered in any way.

Troubleshooting

- **Limited curing capacity**
Check the condition of the fiber optic and make sure it is clean
If the problem persists, call Technical Service
- **Light shuts off before the end of the cycle (LED indicator flashes)**
The overload cutout has tripped; let the light cool down.

Disposal at end of service life

- Never throw out the equipment in regular trash.
- Observe current local regulations regarding disposal of the equipment at the end of its service life.
- Due to the possibility of cross contamination, disinfect the equipment before disposing.


5.7. LEDA curing light

Technical specifications.

Supply voltage:	: 24-36 Vdc
Max. power absorbed:	: 10,5 VA
Light source	: 1 Led da 5 W
Wave length	: 440÷480 nm
Settable times	: 20 – 40 – 60 sec.
Audible signals:	: at start up every 5 sec and at end of cycle
Type of operation	: intermittent (60 sec. work – 13 min. pause).

General description of the light.

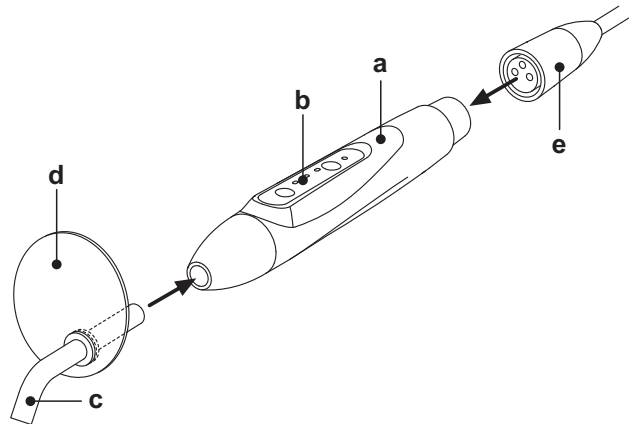
- [a] Light handpiece.
- [b] Control panel.
- [c] Fiber optic.
- [d] Eye protection.
- [e] Power cord.

 **NOTE:** The curing light is delivered in its original packing which should be kept for future shipment.

Description of the control pad.

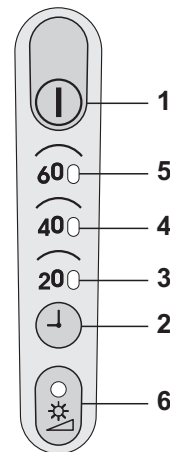- [1] START button.
- [2] Button to set work times.
- [3] LED indicator, time set to 20 seconds.
- [4] LED indicator, time set to 40 seconds.
- [5] LED indicator, time set to 60 seconds.
- [6] LED indicator, Brightness set on increasing scale

Operation.

WARNING!

The instrument is supplied non-sterile.

Before attempting to use the device, disinfect the light handpiece and eye protection. The fiber optic is autoclavable in water steam at 135 °C.

- Insert the fiber optic [c] in the handpiece.

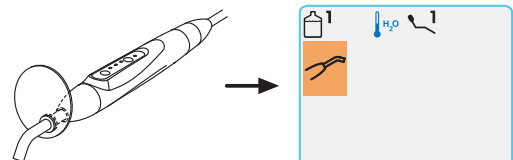
NOTE: be careful to put the fiber optic fully in place in its housing (end of travel) to assure top performance.

- Attach the LEDA light to its power cord by simply snapping it.

NOTA: the mating surfaces must be clean and dry.

- Put the instrument in its work position.
The relative image appears on the console display to signal the instrument has been withdrawn.
- Select the desired time by pressing and releasing button [2]. The led [3], [4] or [5] that lights up indicates the selected time.
- Holding down button [2] does not cause the times to elapse but changes the operating mode of the light. When the button [2] is held down the first time, brightness is gradually increased. The led [6] comes on to signal this condition. If the operation is repeated, the device operates at half power. The led [6] will flash to signal this condition. Press and hold down button [2] until led [6] goes out to restore normal operation.
- Press button [1] to turn on the light. Hold down the button to emit the pointing light. The work light is emitted when the button is released. The light automatically shuts off once the set time is over. Three fast BEEPS are provided in this case. The light can be shut off before the set time has elapsed by pressing button [1] again.

NOTA: a beep is given every 5 seconds during operation.

ATTENTION!

A safety device stops operation of the light if the temperatures established by current safety standards are exceeded. Wait until the light has cooled down before restarting work. When the instrument is not usable because the maximum temperature has been reached, the leds will start to flash.

5.8. C-U2 PRO intraoral camera system

C-U2 PRO is an intraoral camera system specially designed for simple and well-conceived usability in examining the oral cavity complete with extremely lightweight handpiece, auto-exposure and variable focus features. This system is designed to facilitate dentist-patient communication, to show and explain to patients all oral conditions and reasons for planned treatment so that they more clearly understand it. The C-U2 PRO system allows live images of the section in question to be taken through a touch of a fingertip on the touch-sensitive area of the handpiece. The live intraoral images are displayed on the monitor or workstation.

! ATTENTION!

The camera may be used as a tool to aid in diagnosis however, the result must always be compared to visual examination and/or other diagnostic indications. Evaluations and conclusions based only on the image taken by the camera may be poor as the colors and shapes, electronically processed, may not perfectly correspond to the actual ones.

Safety guidelines.

! ATTENTION!

- The disposable infection control sheaths must be used with the device. Change the sheath for each new patient.
- After putting on a new disposable infection control sheath, check it over before using the camera, making sure it is not torn anywhere. If it is, take it off and put on a new one.
- Do not place the handpiece in liquids or in autoclave under any circumstances.
- Store the handpiece in a clean dry area.
- Do not bend the connecting cable excessively.
- Be extremely careful not to drop the handpiece and do not expose it to excessive vibrations.
- Never use a damaged handpiece. Make certain the camera is in good condition and has no sharp edges before attempting to use it. If in doubt, do not use the handpiece, carefully put it away and contact technical assistance.
- Before starting the equipment, check the condition of the lens protection.
- Do not direct the light beam in the operator's or patient's eyes during operation.
- During continuous use (example, more than 10 minutes consecutively), the temperature of the camera's tip usually increases significantly; if this is uncomfortable, put the handpiece in its holder for a few minutes to allow the light source to cool down. When the camera needs to be used a prolonged time, reduce the strength of the light emitted with the cursor on the OSD control panel (see paragraph 5.8.1).
- If left running for extended periods, make sure the temperature of the tip is acceptable before attempting to use the camera. Briefly touch the clear plastic part with your fingertip being careful not to touch the lens in the middle.

Connecting the handpiece.

Attach the camera handpiece [a] to the end of its power cord and tighten the ring nut [b].

! ATTENTION!

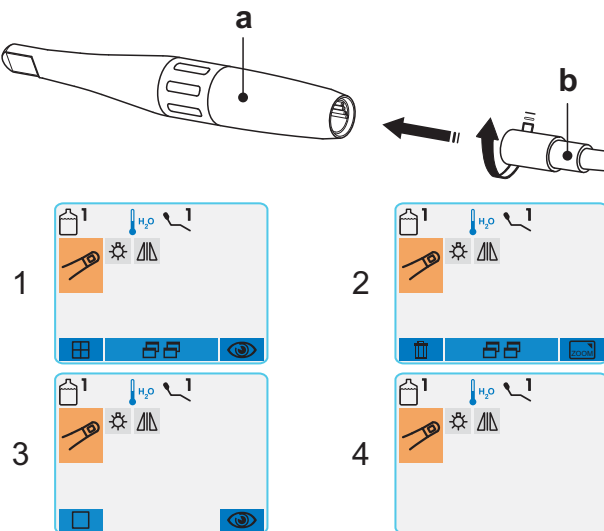
Make sure the cord is firmly screwed onto the handpiece.

! ATTENTION!

Make sure the cord is firmly screwed onto the handpiece.

Camera system usage.

- Put the instrument in its work position.
At this point the camera is on and may be in *LIVE* mode (the monitor shows "live" images) or *FREEZE* mode (the last images frozen appear on the monitor), more precisely:
1- *LIVE* status in multi image mode
2- *FREEZE* status in multi image mode
3- *LIVE* status in single image mode
4- *FREEZE* status in single image mode.

NOTE: If the camera is in *LIVE* mode, the main screen is displayed when the instrument is put back in place.

If the camera is in *FREEZE* mode, the relative menu remains on the monitor when the instrument is put back into place.

- Briefly press the foot control to stop from 1 to 16 images divided in 4 pages on the monitor.

NOTE: The images shown on the monitor by the camera are only temporarily saved. To permanently save the images, connect the camera to a PC that complies to standard IEC 60950 which has a USB 2.0 HIGH SPEED port and image software program.

Turning on the camera's light (led).

Press the button  to turn the camera's led lights on/off.

NOTE: Icon (A) appears on the display to signal the leds have been turned on.

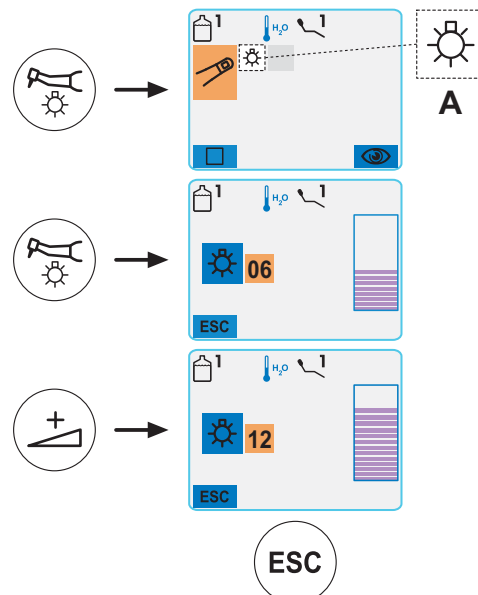
Use buttons INCREASE and/or DECREASE to adjust the brightness level, if necessary.

The set level is automatically saved.

NOTE: the value which can be entered ranges from 1 to 16.

Use the fast selection buttons to automatically set the following values:

- MAX button: level 16,
- MAX button: level 8,
- MIN button: level 1,

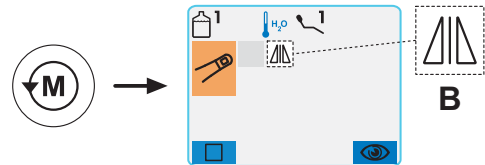
“MIRROR” function.

When real images are displayed, press the button  to change over to mirror images.

This function is possible only in *LIVE* mode.

 **NOTE:** The relevant icon (**B**) appears on the display to signal this function has been activated.

“FREEZE” function.

This camera allows images to be frozen on the monitor.

This function can be activated in different ways:


- With the foot control (see paragraph 5.2.).
- By tapping the touch button [g] on the camera's handpiece.

 **NOTE:** the OSD control panel (see paragraph 5.8.1.) can be used to select if the touch button [g] is to be activated or deactivated leaving only operation with the foot control operative.

In order to go back to the “live” image, simply tap the button again or actuate the foot control.

These images can be displayed in two different ways: single image or multi image.

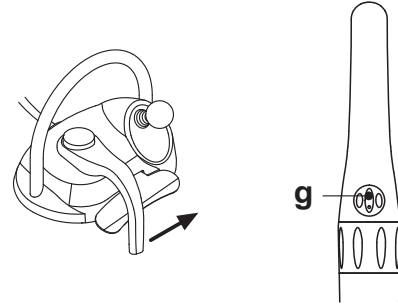In the latter mode, the screen is divided into four parts and 4 frozen images are displayed simultaneously.

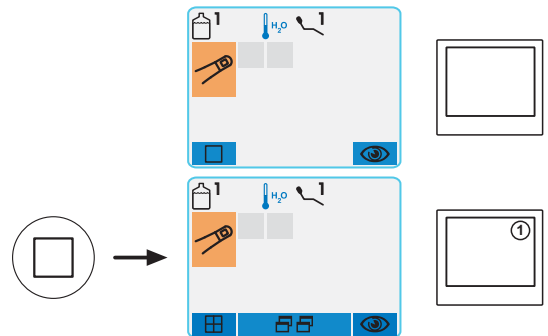
Setting single image or multi image mode.

With the camera turned on and set to *LIVE* mode, simply press the left multifunction button to change over from single image to multi image mode (and vice versa).

The icon for the selected mode appears at the multifunction button:

- icon : single image mode activated,
- icon : multi image mode activated.

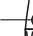
 **NOTE:** The relative icon appears in the top right-hand corner of the monitor to signal multi image mode has been selected.

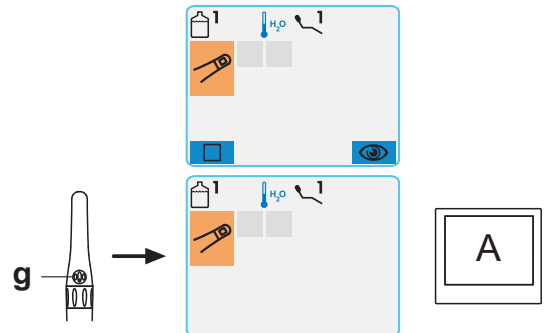
Operation in single image mode.

- Take out the camera in *LIVE* mode and set to single image mode, the “live” image appears on the monitor:

- Tap the touch button [g] on the camera's handpiece (or actuate the foot control) to freeze the image which is immediately displayed on the monitor, deleting any previous images.

- Press the right multifunction button () to view the last image frozen.

 **NOTE:** The last image frozen remains on the monitor even if the camera is put back in place.


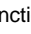
Operation in multi image mode.

- If the camera is withdrawn when set to *LIVE* and multi image mode, the “live” image is displayed and the icon with the number of the storage page active (e.g. 1) appears in the top right-hand corner of the monitor:

- press button MENU () to change the storage page rotating amongst the 4 pages available.
- press the right multifunction button () to view the 4 images already present on the selected storage page,

 **NOTE:** The number of the page selected appears in the middle of the monitor.

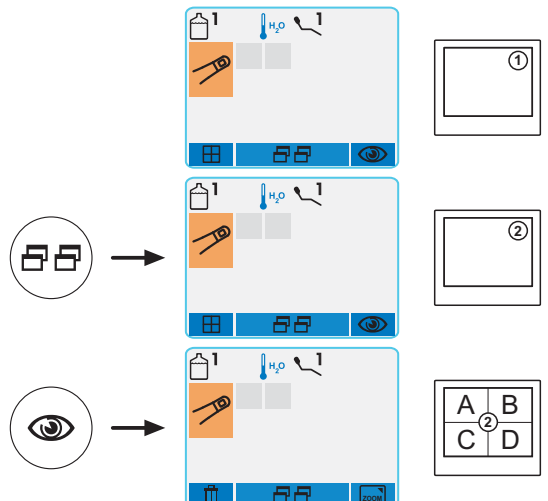- press the left multifunction button () to change over to single image mode.

SKEMA 6 - OPERATING INSTRUCTION

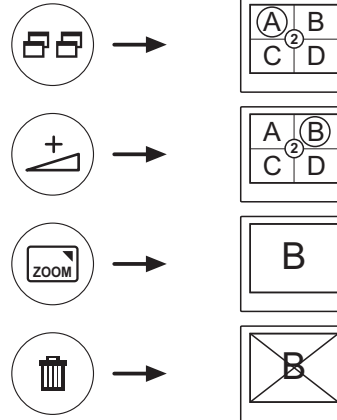
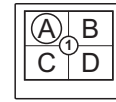
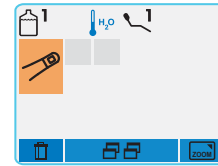
- While displaying (FREEZE mode):

- Press button MENU () to view the 4 pages available one after the other.
- Use buttons INCREASE and/or DECREASE to select one of the 4 images present on the page displayed.
- Press the right multifunction button () to see the full screen version of the selected image.

 **NOTE:** press the button again to display the 4 images.

- Press the left multifunction button () to delete the selected image.


 **NOTE:** : hold down the button () for 3 seconds to delete all 4 images present on the page displayed.

- Saving the image: the frozen image appears directly on the monitor in the first box available (e.g. 1) on the page currently displayed. Each subsequent frozen image is placed in the next box, turning clockwise.

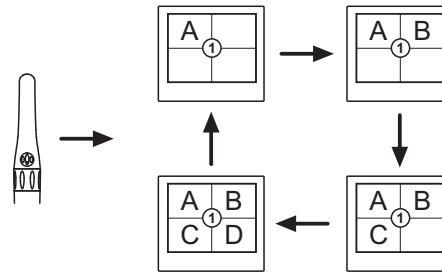
 **NOTE:** Once the 4 boxes available are full, the next frozen images will replace the existing ones, going clockwise.

 **NOTE:** in order to save other images without deleting the previous ones, simply press the MENU button () and change the storage page active.

Description of the handpiece knob

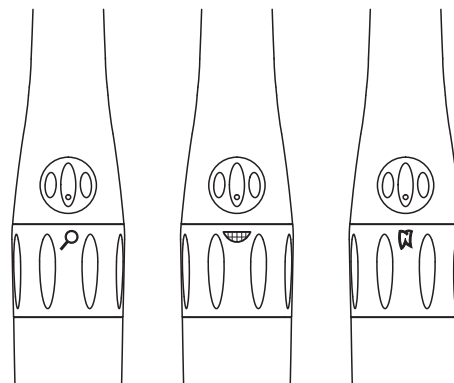- Position “Macro” ().
Captures images in MACRO mode, i.e. with subjects very close to one another.
Object distance: 1 - 3 mm.
- Position “Intraoral” ().
Captures images in INTRAORAL mode, i.e. to focus on areas inside the oral cavity.
Object distance: 13 - 23 mm.
- Position “Extraoral” ().
Captures images in EXTRAORAL mode, i.e. to focus on areas outside the oral cavity.
Object distance: more than 70 mm.

OSD control panel.

This menu allows a number of handpiece C-U2 (see paragraph 5.8.1) operating parameters to be edited. It appears directly on the screen. Press fast selection button MIN to access the OSD control panel. Press multifunction button ESC to quit the OSD control panel.

NOTE: The OSD control panel can also be accessed with the foot control: place the handpiece in the work position and actuate the foot control at least 2 – 3 seconds within the first 10 seconds.

MyRay Grabber

This program allows the C-U2 PRO camera to be setup when it is connected to a PC / WORKSTATION.

A full description of how the MyRay Grabber works is given in the instructions, provided in electronic format, enclosed with the C-U2 PRO camera.

Disposable infection control sheaths

The camera can be a vehicle of cross-contamination between patients. Always use disposable infection control sheaths (code 97901338) and disinfect it on the outside after use every day.

The sheath (with white paper backing) is enclosed in two protective layers: a transparent one with blue tab at the front and a paper one at the back. Follow the directions below to install a new infection control sheath:

1. Insert the camera handpiece tip between the layer with white tab and the rear paper backing. The lens, surrounded by the LEDS, must face downwards towards the rear paper layer. Gently push the handpiece to the end of the sheath.
2. Pull the blue tab removing the protective covers.
3. The dental camera is now protected and ready for use.

ATTENTION!

- Always make certain the handpiece is correctly inserted inside the infection control sheath.
- Always change the disposable infection control sheath before using the device on a new patient. This is a fundamental step in preventing cross-contamination.
- Disposal: the disposable infection control sheaths are to be treated as special waste materials (like surgical gloves).

Cleaning and disinfecting.

Clean the handpiece after each use (see paragraph 1.4.).

ATTENTION!

- The intraoral camera is not designed for cold sterilization by being soaked, for example in solutions such as glutaraldehyde or hydrogen peroxide.
- All products must be used as directed by the manufacturer.
- All material used to clean and disinfect must be thrown away.

Maintenance and repairs.

The C-U2 PRO camera system does not require any particular maintenance. In the event of malfunctions, please send back the complete handpiece.

ATTENTION!

There are no parts that can be repaired on site. In the event of a malfunction, please contact an authorized dealer.

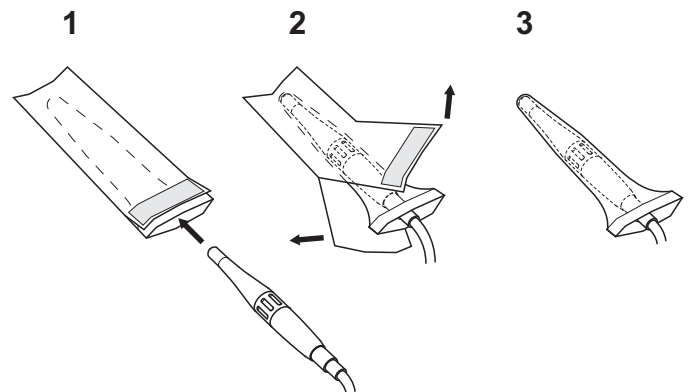
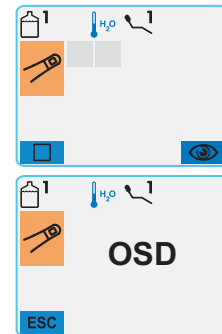
Returning parts.

- Please send back any defective devices in their original packaging. Do not reuse damaged boxes.
- The device must be disinfected before being shipped to prevent cross-contamination. Handpieces that have not been adequately cleaned and disinfected will not be accepted.

ATTENTION!

The sender shall be held responsible for any equipment damaged incurred during shipment regardless of Whether or not the devices are under warranty.

5.8.1. OSD control panel (On Screen Display)

This menu allows the user to change some operating parameters of the C-U2 PRO handpiece.

The menu displayed is automatically cleared if the foot control or touch key is not pressed for a few seconds or if the handpiece is put back in place.

The selected options are automatically saved.

All the parameters that regard the aspect of the image are saved in the handpiece therefore they accompany it if it is used in more than one dental unit.

 NOTE: *The factory settings are shown in the figures below.*

How to navigate inside the menus

- Press the foot control to scroll the items available (“vertical” scroll).

 **NOTE:** The active option is always the one highlighted in yellow.

- Press the touch key on the handpiece to select the possible options for a certain item ("horizontal" scroll).
- Use buttons INCREASE and/or DECREASE to change the parameters given on the adjustment menu.
- Press the foot control to go to the next data item. The selected option is automatically activated.

Menu layout.

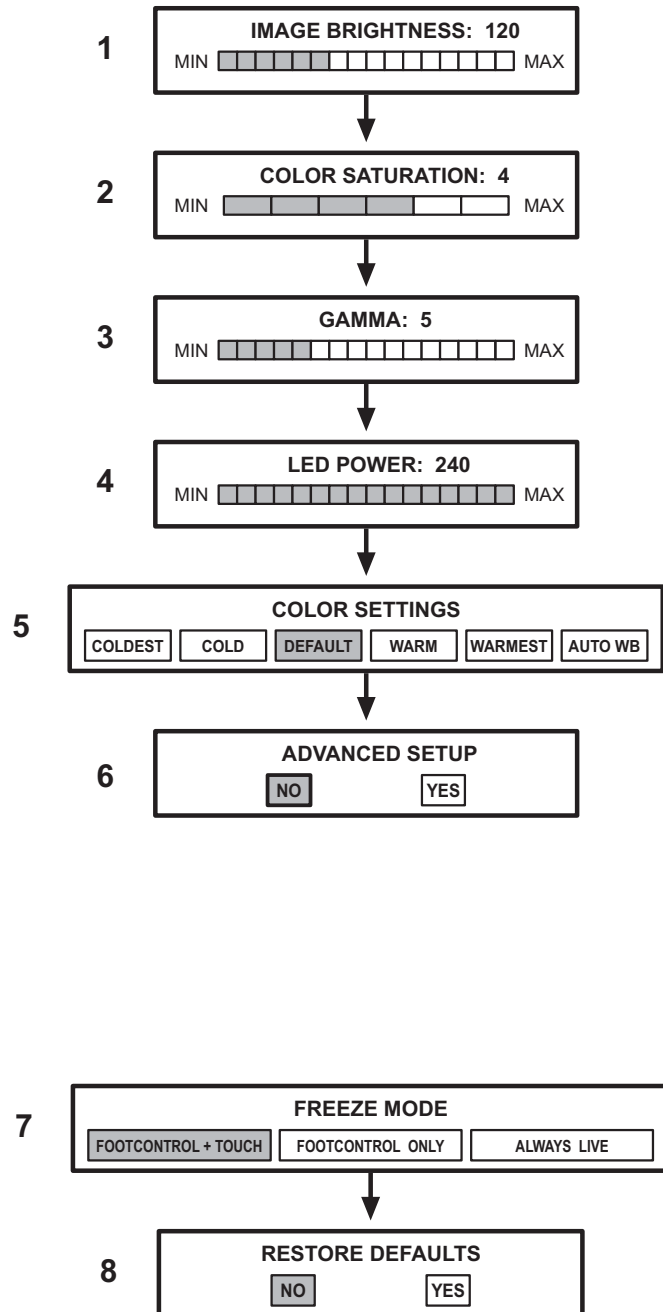
The data items given in the menu are shown below in the order in which they appear on the display screen:

- 1 **IMAGE BRIGHTNESS:** increases or decreases the brightness of the image displayed according to the power of the lighting unit. It affects the exposure time.
- 2 **COLOR SATURATION:** : increases or decreases color saturation. If set to minimum the image appears in black and white, at maximum the color is very strong.
- 3 **GAMMA:** similar to a contrast control but with a stronger effect. It changes how dark/bright areas render in the image; the image is "softer" at the lowest value; it is "harder" at the maximum value.
- 4 **LED POWER:** used to adjust the brightness of the lighting unit. It should be used if saturation cannot be eliminated (amalgam, etc..) even if brightness is decreased.
- 5 **COLOR SETTINGS:** adjusts the color profile.
It allows the user to gradually change from a "colder" profile (to the left) to a "warmer" one (to the right). "Auto WB" activates the automatic white balancing function which is to be used under particular conditions in which the previous profiles are not satisfactory.
- 6 **ADVANCED SETUP:** last data item on the menu it accesses options used less frequently. Actuate the foot control to go back to the beginning to continue modifying other parameters. On the other hand select "YES" (with the touch key) and actuate the foot control to go to the advanced menu.

Advanced menu.

This menu provides access to seldom used configuration options.

- 7 FREEZE MODE:** menu to set the “Freeze” function (image freeze).
- **Footcontrol + Touch:**Function “Freeze” can be activated both with the foot control and touch button for the C-U2 PRO camera.
 - **Footcontrol only:** function “Freeze” is turned on only with the foot control.
 - **Always Live:** function “Freeze” disabled, “live” mode is used.
- 8 RESTORE DEFAULTS:** restores the original factory settings. Press the foot control to go back to the basic level menu (“Image Brightness”). Select “YES” to restore the factory settings.

5.9. Integrated laser L810

This is a diode laser for medical use.

Use.

- **Operating times: 6 min. ON, 3 min. OFF.**
- Place the instrument in its work position (the screen regarding laser extraction appears on the display of the instrument board).

- 1 When the foot control is pressed the first time (see paragraph 5.2.) the pilot beam is turned on and a buzzer rings intermittently at the same time.

 **NOTE:** when the foot control is released, the pilot beam is emitted for about 1 min. before it automatically shuts off (time-out).

- 2 Press the foot control another time to turn on the laser beam. In this case the buzzer will ring continuously.
- 3 When the foot control is released, the LASER beam shuts off and the pilot one automatically comes back on (a 1 min. time-out is always used).

 **NOTE:** While the laser is functioning, the other instruments on the board along with all dental chair movements are disabled.

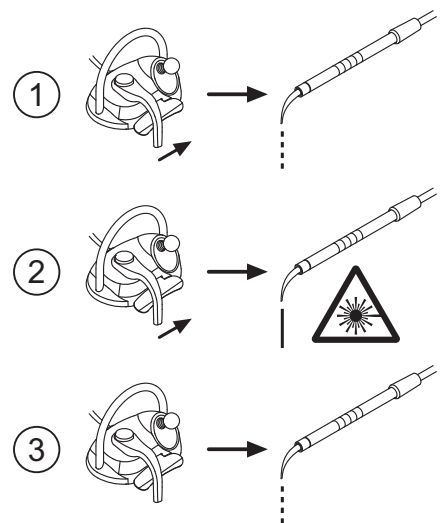
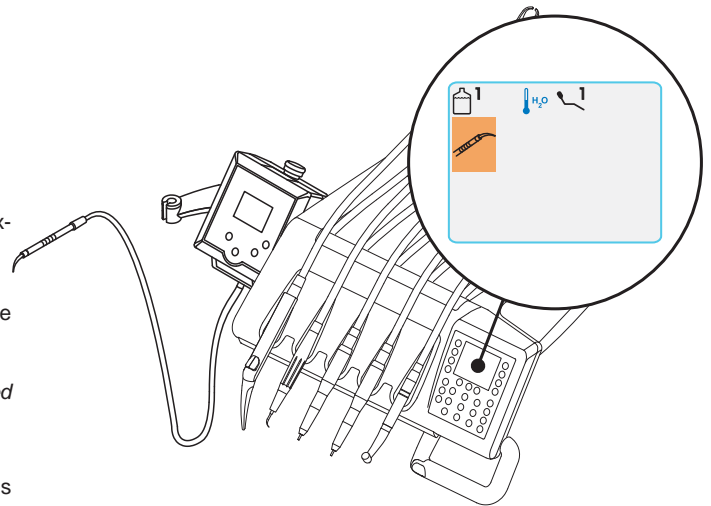
ATTENTION!

See the specific Operations and maintenance instructions enclosed with the equipment for information about all the other laser device functions.

Safety guidelines.

ATTENTION!

- This laser device is classified as a Class 4 laser. Only expert personnel adequately instructed on safe operation may operate this laser device.
- The laser beam should be directed on the patient taking all the required precautions set forth by current safety regulations. Any other use or use in any other way is considered erroneous use of the device.
- The manufacturer shall not be held responsible for any consequences deriving from incorrect use of the instrument, the specific use includes:
 - observe the instructions given in the specific Operation and maintenance manual.
 - observe all national and international standards (CEI EN 60825-1, 2003-02).

5.10. Peristaltic pump

This device allows saline solution to be supplied through a single-use administration line without any contact.

This device is provided together with the micromotor and scaler.

NOTE:

- In order to use with the micromotor, contrangles with external cooling or for hollow burs are required (Goldspeed R20-L type).
- In order to use with the scaler, handpiece SURGISON 2 is required.

Description of the symbols present on the device.

- 1) Material meets and exceeds the essential requirements of directive EEC 93/42 and subsequent changes.
- 2) CAUTION: PINCH HAZARD.
Do not put your fingers in rotating parts.
- 3) Material sterilized with ethylene oxide
- 4) Expiration date (yyyy-mm).
- 5) Single-use material.
- 6) Material identification code

Putting into service

- Direct and put the IV drip pole [a] in place and hang the flask or bottle [b] that contains the saline solution.
- Open the bag [c] and take out the sterile administration line.

WARNING!
Use sterile disposable gloves.

WARNING!
Check the condition of the packaging as well as the expiration date of the administration line. Only CASTELLINI administration kits guarantee proper trouble-free operation. These lines are sterile and disposable, reusing them may put patients at risk of microbiological contamination.

- Open the cover [d] of the peristaltic pump, turning it upwards.
- Attach the tube, being careful to place the part with the largest diameter inside the pump's V seats. The pump rotates clockwise. Place the tube so that the section that runs from the bag enters from the left side of the pump (see figure).
- Close the cover [d]. If it does not close, open the cover again and check the position of the tube.

WARNING!
Do not start the pump with the cover [d] open, finger pinch hazard.

- Pierce the cap [b] of the bottle of saline solution with the outflow tip of the administration line [c].
- Attach the administration line to the instrument cord using the plastic clips provided in the sterile kit.

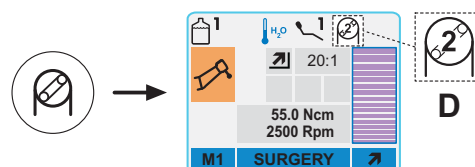
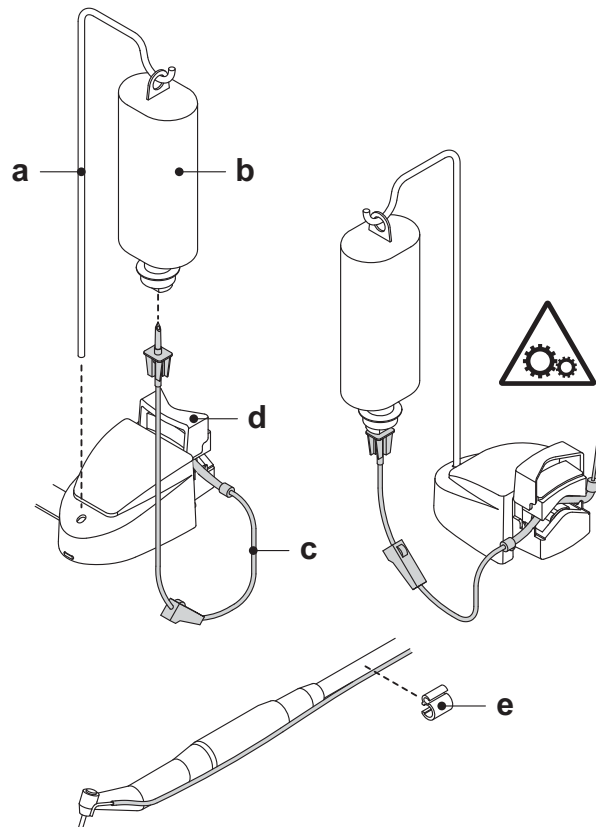
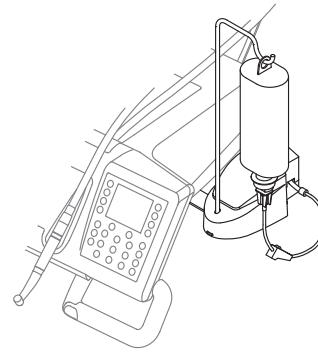
NOTE: use type A for the scaler's cord and type B for the micromotor's cord.

Operation

- In order to start the peristaltic pump, withdraw the desired instrument (micromotor or scaler) and press the relative button on the doctor's console.

NOTE: the icon (D) appears on the display and a BEEP is emitted to signal the peristaltic pump has been started.

- If necessary, set the amount of saline solution delivered by the pump (see paragraph below).
- At this point, the internal air and water supply is interrupted and the peristaltic pump will automatically deliver saline solution when the instrument is started.
- To shut off the peristaltic pump, press the button again with the instrument withdrawn. .


Setting the amount of saline solution delivered by the pump.

- Hold down the button (at least 2 seconds) to set the amount of saline solution delivered by the pump.
- Use the increase and/or decrease buttons to select the desired amount.

NOTE: a value ranging from 1 to 5 can be used.

The amount of saline solution delivered that is associated to the settable values is given below:

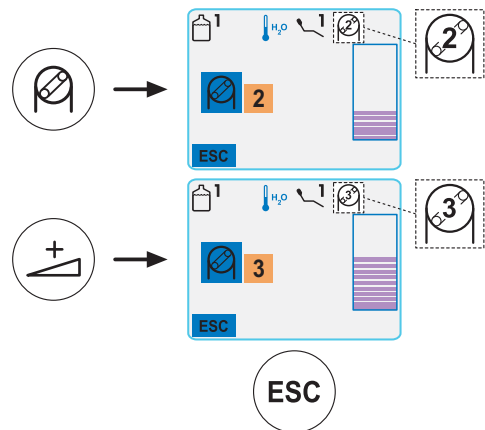- value 1: approximately 35 cc/min
 - value 2: approximately 50 cc/min
 - value 3: approximately 70 cc/min
 - value 4: approximately 90 cc/min
 - value 5: approximately 100 cc/min
- Press function key ESC to quit this submenu and confirm the selected value.

Faults.

If the peristaltic pump fails, the relative icon will appear in red on the display.

In this case, call Technical Service.

5.11. Electronic apex locator (LAEC).

By analyzing resistance to a small electric current, LAEC is able to help in locating the root apex. If used together with a "file" (not provided) for manual treatment, it is an important tool in measuring the length of the canal itself.

Along with regular manual use in this dental unit, it is possible to use the apex locator together with "ENDO" mode of the micromotors and scaler. The position of the instruments used on the handpieces can be monitored as the LAEC signals directly reach the files, through the instrument cords, thereby allowing the position of the canal to be monitored during treatment.

Description of the components.

- [1] LAEC outer cable.
- [1.1] LAEC outer cable - neutral pole.
- [1.2] LAEC outer cable - active pole.
- [2] Hook electrode.
- [3] Probe.
- [4] LAEC clip connecting forceps.
- [5] Outlet for LAEC outer cable .

Operation.

- In this dental unit, the LAEC automatically turns on when the outer cable [1] is inserted in the outlet [5] placed under the instrument board. When turned on, the menu for setting the alarm threshold appears on the display (see paragraph 5.1.1.2.8.).
- Applying the electrodes:
 - Connect the hook electrode [2] to the neutral pole [1.1] and place it on the patient's lip.
 - Connect the active pole [1.2] to the file (not provided) placed on the root canal; connection to the file can be done using the probe [3], forceps [4] or with the fittings provided for the handpieces.

IMPORTANT INFORMATION!

The electrodes are not sterile when supplied.



Indications given on the display.

- The bar graph on the left side of the display shows the position of the file in relation to the apex. Numbers "1, 2, 3" represent the distance between the instrument and apex.
- Icon APEX shows the distance of the instrument from the apex.

NOTE: indication "> 4" means that the file is too far away from the apex in order for it to be measured.

- Icon ALARM shows the set alarm threshold. The alarm threshold identifies the distance between the instrument and apex, beyond which a beep is emitted that gradually increases as it gets closer to the apex. See paragraph 5.1.1.2.8 for instructions on how to set the alarm threshold.

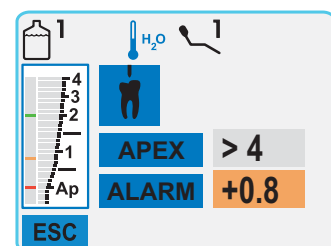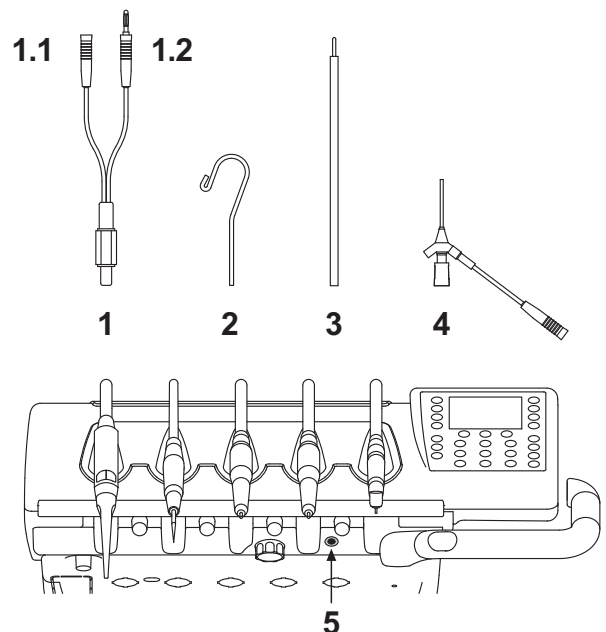The graphic and numeric readings are constantly updated while the file is being introduced in the canal.

SKEMA 6 - OPERATING INSTRUCTION

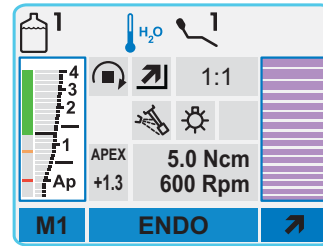
LAEC combined with the electric micromotor.

It is possible to use the LAEC together with the electric micromotor when operating in ENDO mode.

When the LAEC is turned on, information about the micromotor and the LAEC (bar graph and APEX value) is displayed when the electric micromotor is withdrawn in ENDO mode.

While the electric micromotor is running, the keys are associated to the instrument's functions and therefore it is not possible to modify the alarm threshold of the LAEC if the instrument is not put back in place.

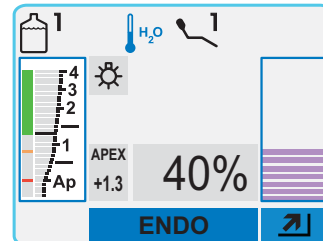
LAEC combined with scaler.

It is possible to use the LAEC combined with the scaler when it is in ENDO mode.

When the LAEC is turned on, information about the scaler and the LAEC (bar graph and APEX value) is displayed when the scaler is withdrawn in ENDO mode.

While the scaler is running, the keys are associated to the instrument's functions and therefore it is not possible to modify the alarm threshold of the LAEC if the instrument is not put back in place.

 **NOTE:** with the scaler, it is necessary to take the measurements with the instrument stopped as good electrical contact cannot be obtained if the file is vibrating. Simply release the pedal for a few seconds to correctly measure the position of the instrument in relation to the apex.

Measuring the length of the root canal.

- Use of the manual file is fundamental to measure the canal. The correct procedure consists of placing the file into the canal and advancing until 0.5 is reached.
- Move the file further forward slowly turning clockwise until APEX appears on the instrument.
- Once APEX appears, back up the file turning it counter-clockwise until 0.5 is obtained again. Place the rubber stop at the occlusal surface as the reference point for determining the length of the root canal.
- Take an x-ray to verify the file is in the correct position.
- Remove the file from the canal and measure the length with a ruler. Subtract 0.5-1 mm from the measured value.

IMPORTANT INFORMATION!

Always use the LAEC together with an x-ray to precisely determine the position of the apex.

Various morphological situations which are not always foreseeable may result in readings that are not always precise.

For example:

- excessively wide root canal;
- shrinkage;
- broken roots;
- presence of metal crowns,

6. Assistant's board operation

Main features:

- Two articulated arms secure the board (a) to the hydrogroup (b) allowing it to be placed in the most convenient work position.
- The pantograph arm (e) allows the assistant's board to be moved 335 mm vertically through 6 work positions.

NOTE: press button [h], found on the pantograph arm, to move the assistant's board fully down.

- The assistant's board (a) comes with a control console (d) with buttons used to operate the dental chair and hydrogroup.
- The assistant's board can be equipped with 2 suction tubes and 3 instruments.
- The assistant's board comes with sliding rollers (f) that guide and hold up the suction tubes.

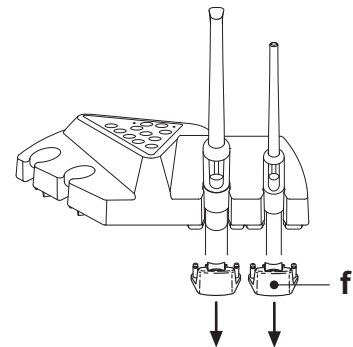
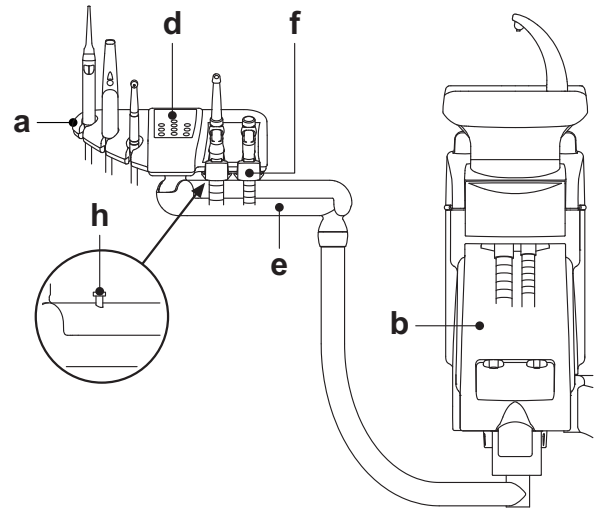
NOTE: The assistant's board is equipped with a safety device that locks out dental chair movement when the board's arms are obstructed.

Cleaning the sliding rollers.

Push down and take off the sliding rollers (f).
Clean the sliding rollers with a suitable product approved against bacteria and fungi and effective against viruses (including HBV and HIV): see paragraph 1.4.

NOTE: The product should be sprayed on soft disposable paper towel or sterile gauze.

ATTENTION!
Do not use acids or harsh products.


6.1. Assistant's touchpad

Description of the buttons:

- Dental operatory light on/off button
- Cup fill button
- Bowl rinse button
- Automatic return position retrieval button
- Rinse position button
- Programmed position 1 retrieval and chair base up button.
- Programmed position 2 retrieval and chair back up button.
- Programmed position 3 retrieval and chair base down button.
- Programmed position 4 retrieval and chair back down button
- Disinfecting procedure start button

NOTE: **functioning of buttons that move the chair**

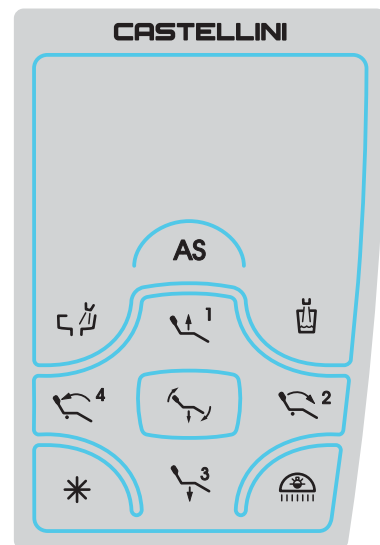- Briefly pressed: programmed position automatically reached.
- Held down: manual positioning.

6.2. Instruments on assistant's board

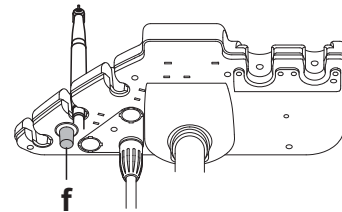
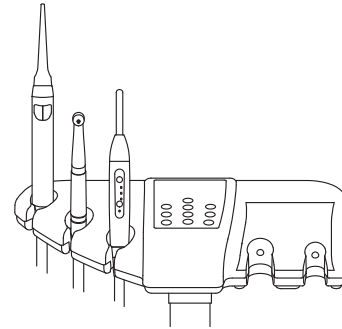
All the instruments provided on the assistant's board operate in the same manner as those on the instrument board.

More specifically:

- Syringe, see paragraph 5.3
- Turbine, see paragraph 5.4
- Scaler, see paragraph 5.6
- LEDA curing light, see paragraph 5.7
- C-U2 camera, see paragraph 5.8

NOTE: the turbine and scaler are equipped with cock [f] to adjust the water spray found under the assistant's board. The amount of air spray cannot be adjusted in the turbine.

6.3. Suction tubes

To adjust the suction power, use the slider [a] located on the tube handpiece.

The SKEMA 6 comes with the V.D.S. system that allows the suction line to be dried using an automatic stop delay (about 2 seconds).

Removing the suction tubes

ATTENTION!

Always wear gloves to prevent contact with infected material When removing the suction tubes.

Remove the suction tubes from the conveyor fittings by turning and twisting the tube fitting. Detach the suction tubes from the holders by turning and twisting the tube fitting.

ATTENTION!

Never directly grasp the suction tube.

Cleaning the suction tubes.

The suction tubes need to be regularly cleaned and disinfected (see paragraph 9.4). Use a product specially designed for this purpose.

ATTENTION!

Use only a diluted solution of STER 3 PLUS (Castellini) as a liquid disinfectant. As concerns the amount of disinfectant and percentage of solution follow the instructions provided with the product.

Sterilization.

- Suction tube holders: water steam autoclave up to 135°C (2 bar), minimum time: 15 min.
- Suction tubes: soak to cold sterilize.

ATTENTION!

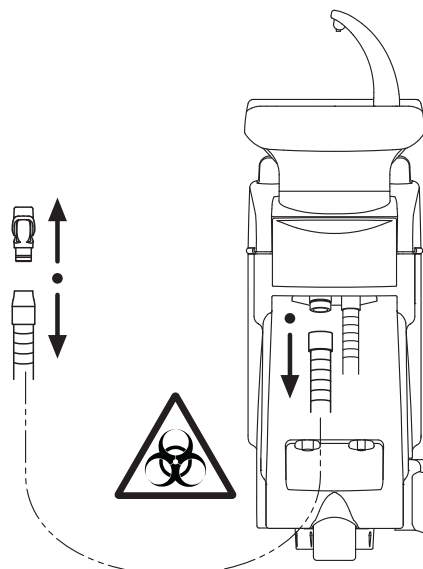
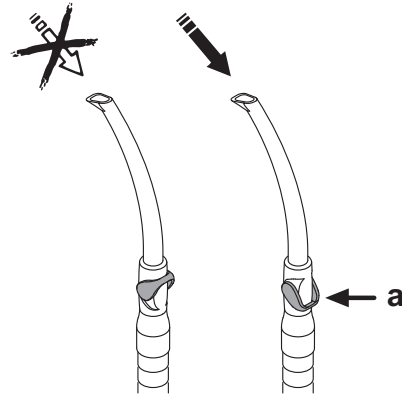
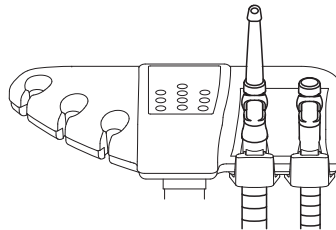
Never use procedures in which the temperature goes over 55 °C with the tubes.

Maintenance

Periodically lubricate the o-rings of the suction tube holders with S1 – Protection for o-rings lubricant (see paragraph 9.4).

Note about biocompatibility

Only suction tubes supplied with the dental unit and there after original replacement tubes can be used. The suction tubes must comply to the standard ISO 10993-1 Biological evaluation of medical devices.

6.4. Instrument tray on assistant's board.

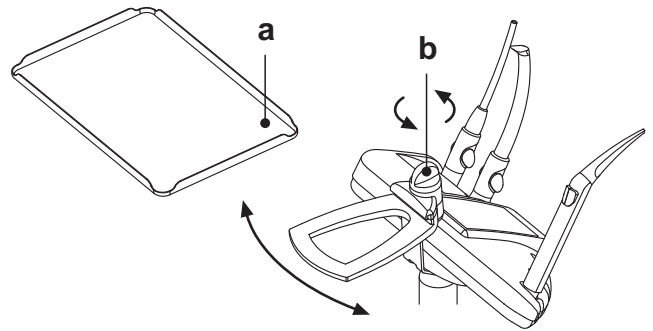
The instrument tray [e] is constructed in stainless steel and can be conveniently taken off its support. The tray holder can be turned either clockwise or counter-clockwise, allowing it to be placed in the most convenient position for the operator.

To lock/unlock the tray holder, simply using the clutch knob [b].

ATTENTION!

Maximum allowable load that can be applied on instrument tray: 1 Kg evenly distributed.

6.5. Hydraulic saliva ejector

The hydraulic saliva ejector starts running When the tube is removed from the support.

Cleaning after each use.

Draw in about half a liter of hydrogen peroxide with maximum 3% concentration (10 volumes).

Cleaning the saliva ejector filter.

NOTE: This operation must be carried out at the end of each work day.

ATTENTION!

Put on gloves before attempting to perform this operation!

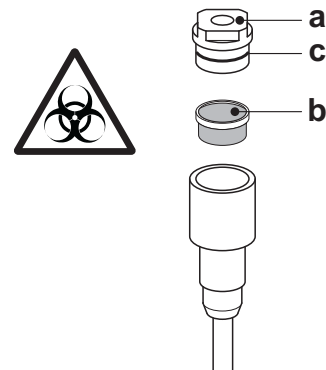
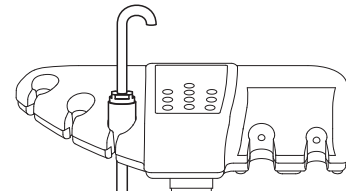
- Take off the cap [a] by turning and twisting at the same time.
- Remove the filter [b].
- Clean/replace the filter (code 97290060).
- Put the filter and cap back in place.

NOTE: To prevent liquids and matter from dripping from the filter taken out, suck air only for approximately 5 seconds before performing these operations.

Routine maintenance.

Lubricate the o-rings [c] with S1 – Protection for o-rings lubricant.

7. Hydrogroup operation

7.1. Fill cup and bowl

The bowl can be turned 305° on the hydrogroup either by hand or with the drive unit.

The bowl and water to cup spout can be removed to ease cleaning.

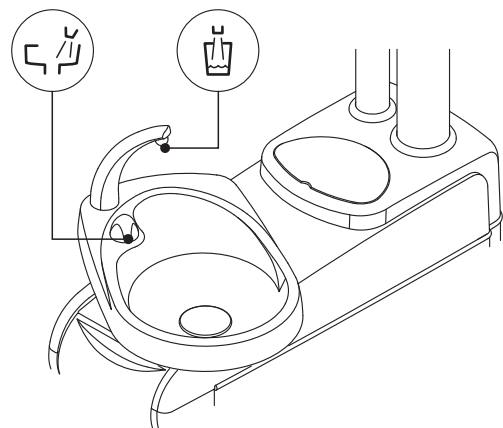
Control buttons

Water to cup button

Bowl flush button.

Adjusting the amount of water used to fill the cup.

See paragraph 5.1.1.4.1.

Adjusting the temperature of the water sent to the cup.

See paragraph 5.1.1.4.2.

Setting bowl flushing.

Water can be delivered to the bowl either in manual (with the ON/OFF button) or timed.

See paragraph 5.1.1.4.1 to set the desired mode of operation and water delivery time.

Setting automatic bowl flushing.

The bowl is automatically flushed in the following cases:

- when button "Water to cup" is pressed,
- when button "Dental chair automatic return" is pressed,
- when button "Rinse position" is pressed.

To change operation see paragraph 5.1.1.4.1.

Powered bowl movement.

Bowl counter-clockwise button.

Bowl clockwise button.

NOTE: The bowl can also be moved directly by hand.

Automatic bowl motion (only with powered bowl).

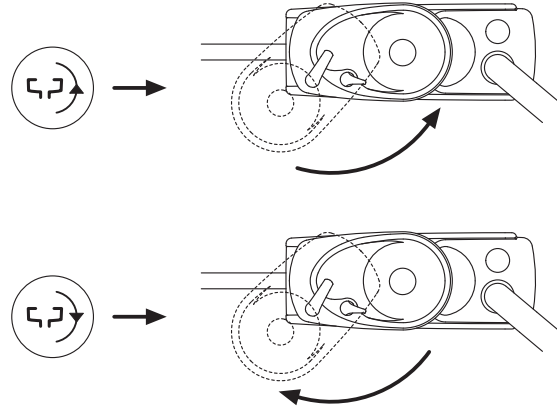
The bowl moves automatically in the following cases:

- press button "Dental chair rinse position",

NOTE: in this case, the position of the bowl can also be set (see paragraph 5.1.2.).

- by pressing button "Dental chair reset position",
 - by pressing the button to retrieve 1 of the 4 programs available.
- See paragraph 5.1.1.4.3 for information on how to modify operation.

Removing the spouts, bowl and bowl filter.

- Pull up the spouts [l] and [n] and remove them.
- Pull up the filter [q] and its cover [p] to remove them.
- Remove the bowl [m] by first turning it counter-clockwise.

Disinfecting and cleaning.

ATTENTION!

Always wear gloves to prevent contact with infected material When cleaning the bowl and bowl filter.

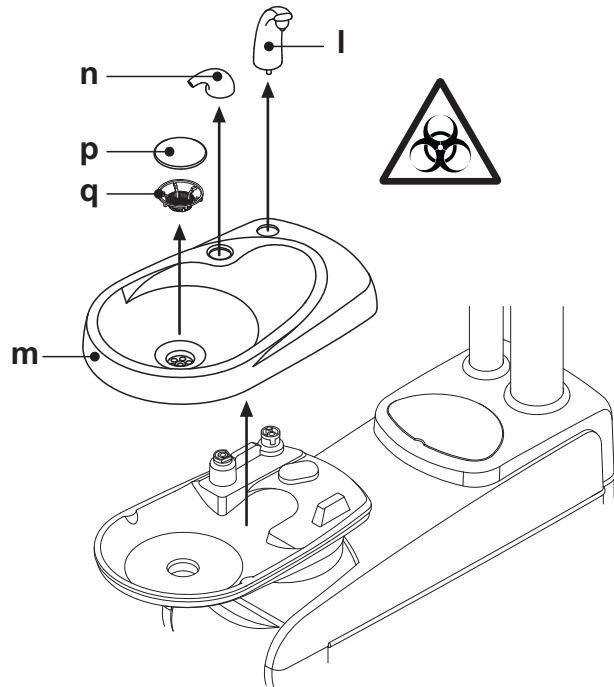
The parts are to be cleaned daily at the end of each work day.

- Spouts and bowl: thoroughly wash with a specially formulated cleaner.
- Bowl filter: clean with running water and commercially-available cleaning products.

ATTENTION!

Do not use acids or harsh products.

7.2 Distilled water tank system

Description of the system

The system is equipped with a tank (a), placed under the hydrogroup, that can contain distilled water.

The tank can hold 1.8 liters.

Distilled water is delivered to:

- the sprays of all the instruments found on the instrument and assistant's board.
- The syringe on the assistant's board.
- to fill the cup
- water quick-connect coupler (if present)

Button (b) present on the assistant's board's control panel is used to shut off the S.S.S.

NOTE: the relative icon (A) on the console display signals with the S.S.S. system is turned on.

A disinfecting cycle for the instrument spray ducts can be carried out with hydrogen peroxide with this system (see paragraph 7.2.1.).

Tank reserve level.

When the liquid in the tank falls below the reserve level, the relative icon (B) appears on the instrument board's console.

NOTE: the icon (B) is black if the S.S.S. system is turned on while it is red if the S.S.S. is turned off.

Filling the tank.

When the water level in the tank is low (about 500 cc), fill it as directed below:

- Push lever (b) in the direction of "CLOSE AIR PRESSURE".
- In this mode, the icon (A) appears in a different color on the display.

NOTE:

While performing this operation, the pressurized air in the tank is automatically released outside and the relative icon (C) appears on the console.

- Turn the tank counter-clockwise (a) and remove it.
- Pour distilled water into the tank until the maximum level is reached.

ATTENTION!

Use only distilled water which 600 parts per million (ppm) of hydrogen peroxide (20 ml of 3% hydrogen peroxide in 1 liter of distilled water) can be added to.

- Put the tank back in place turning it clockwise.

ATTENTION!

Assure the tank is firmly secured in place.

- Press the switch (b) again in direction "CLOSE AIR PRESSURE" to confirm the tank has been filed
- The icon (A) on the display goes back to its original color.

NOTE:

The tank will be pressurized only after an instrument has been extracted or when water is delivered to the cup.

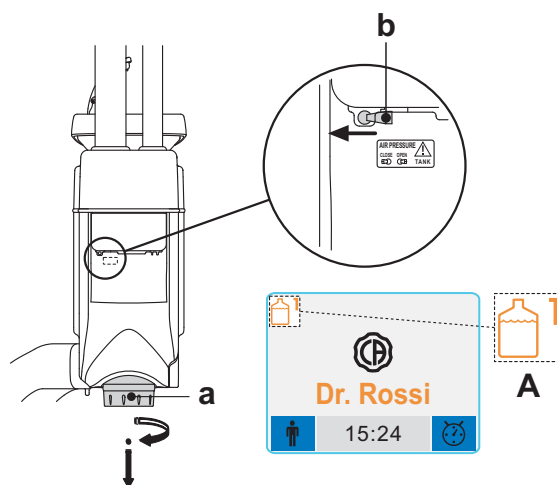
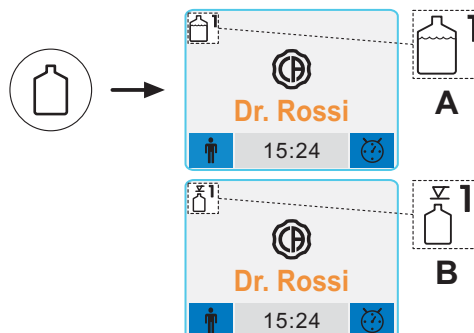
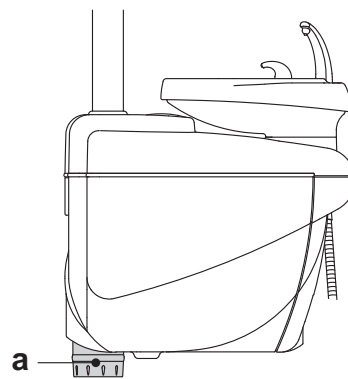
ATTENTION!

- Fill only with distilled water.
- When you are going to be away from the office for an extended period of time (holiday), always empty the hydrogen peroxide tank (a) before leaving.

Cleaning the tank.

Only the tank should be cold disinfected on a regular basis (at least once a month) using a peracetic acid based product proceeding as directed below:

- take the tank out of the dental unit and empty it;
- prepare the solution of peracetic acid based product following the instructions provided by the manufacturer;
- fill the tank up to the rim with the peracetic acid based product;
- let the peracetic acid solution soak in the tank for the time stated by the manufacturer.
- empty the peracetic acid solution from the tank;
- rinse the tank with distilled water;
- fill the tank with distilled water adding 600 ppm of hydrogen peroxide (20 ml of 3% hydrogen peroxide in one liter of distilled water),
- put the tank back in place in the dental unit.

7.2.1. Manual disinfecting cycle with hydrogen peroxide

With the S.S.S. cycle it is possible to perform a manual disinfecting cycle with hydrogen peroxide for the water ducts of all the instruments found on the instrument board and the syringe on the assistant's board.

To disinfect, proceed as directed below:

A) Prepare the disinfectant:

- Pour 3% pure hydrogen peroxide (10 volumes) without diluting it into the pink bottle.

B) Putting in the disinfectant:

- Replace the bottle of distilled water [a] with the one containing hydrogen peroxide (see paragraph 7.2).

NOTE: make sure there is at least 300 cc of liquid present.

- Make certain the spray cocks [d] found at the bottom of the board are opened (if they are not, either water does not come out or too little water flows out).
- Press the key to fill the cup 5 (five) consecutive times, thereby filling 5 cups with water. This step is very important as it assures that the distilled water in the ducts is replaced by hydrogen peroxide.
- Withdraw the instruments one at a time and let the water run for at least 2 minutes using the CHIP-WATER function of the foot control for the dynamic instruments (see paragraph 5.2) and the water key for the syringes.

NOTA: at this point, the ducts contain hydrogen peroxide.

- Put the instruments back in place.

C) Disinfectant contact time:

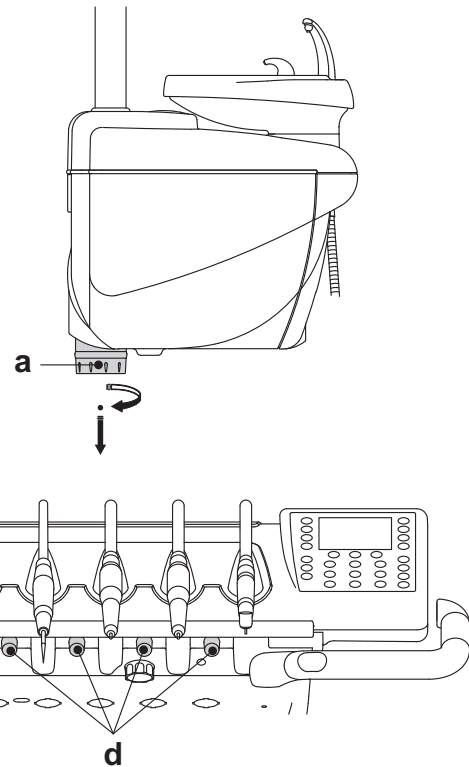
- Let the hydrogen peroxide act for at least 10 minutes but not more than 30 minutes.

D) Rinsing the ducts:

- Replace the bottle of distilled hydrogen peroxide with the one containing distilled water (see paragraph 7.2).
- As during the previous step, take out the instruments one at a time a let the water run at least 2 minutes using the CHIP-WATER function of the foot control for the dynamic instruments (see paragraph 5.2) and the water key for the syringes.

NOTA: at this point, the ducts contain distilled water again and the dental unit is ready to be used.

ATTENTION!

- Once disinfected, make sure the hydrogen peroxide tank is closed (it becomes less effective When exposed to air).
- It is good practice to perform a disinfection cycle at least once a day, preferably at the end of the day.
- Do not leave hydrogen peroxide in the ducts for more than 30 minutes under any circumstances.

Storing the hydrogen peroxide.

Follow the manufacturer's instructions given on the package to correctly store the hydrogen peroxide solution.

ATTENTION!

Never leave hydrogen peroxide in the tank [a] for more than one month.

When you are going to be away from the office for an extended period of time (holiday), always empty the tank [a] before leaving.

7.3. Continuous disinfecting system combined with the device for separation from potable water line to prevent backflow in compliance with European standard EN 1717 (M.W.B. system)

The M.W.B. system assures secure physical separation of the dental unit hydrogroup from the municipal water line thanks to a section with free-fall water. In addition, the system continuously lets hydrogen peroxide into the water circuit obtaining a final concentration of 0.06% in the ducts which is suitable for bacteriostasis.

Description of the system.

The M.W.B. system is located in the utility service center and is always active.

A specific icon [a] appears on the doctor's console to signal the M.W.B. system is turned on.

A tank [C], located in the hydrogroup, holds approximately 970 cc of hydrogen peroxide to feed the system.

The M.W.B. system automatically turns off when the S.S.S. system (if present) turns on.

Signaling the hydrogen peroxide level is low.

When the hydrogen peroxide present in the tank [a] is about to run out (at about 230 cc.) the icon [C] displayed changes into [D] and 1 BEEP is emitted.

NOTE: the icon (D) is black if the M.W.B. system is turned on while it is red if the M.W.B. system is turned off.

ATTENTION!

If the hydrogen peroxide runs out, the operatory continues running but uses municipal water that is NOT treated.

Filling the hydrogen peroxide tank.

Follow the instructions given below to fill the hydrogen peroxide tank:

- Remove the hydrogen peroxide tank [a].
- Take off the cap [k] and fill the tank with hydrogen peroxide.

NOTE: the shape of the cap allows it to be used as a funnel to more easily fill the tank.

- Put the cap and tank back into place.
- Close the cover on the side of the hydrogroup.

ATTENTION!

Fill only with a 3% hydrogen peroxide solution (10 volumes).

Emptying the M.W.B. system tank

This function permits the M.W.B. system's water circuit to be drained, if the dental unit is going to be left off for several days.

The procedure for emptying the tank is given in paragraph 5.1.1.3.3.

Fault messages shown on the console display.

If the system detects a malfunction, a fault message appears on the display screen (see paragraph 11.).

If the fault found is minor, the dental unit continues to operate. On the other hand, if a serious fault is detected, the dental unit shuts down and it is necessary to call technical support.

Hydrogen peroxide storage.

Follow the manufacturer's instructions given on the package to correctly store the solution of hydrogen peroxide.

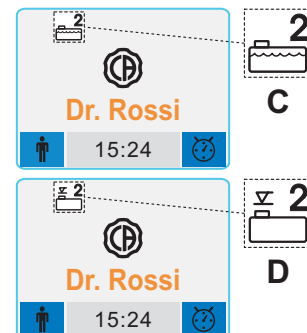
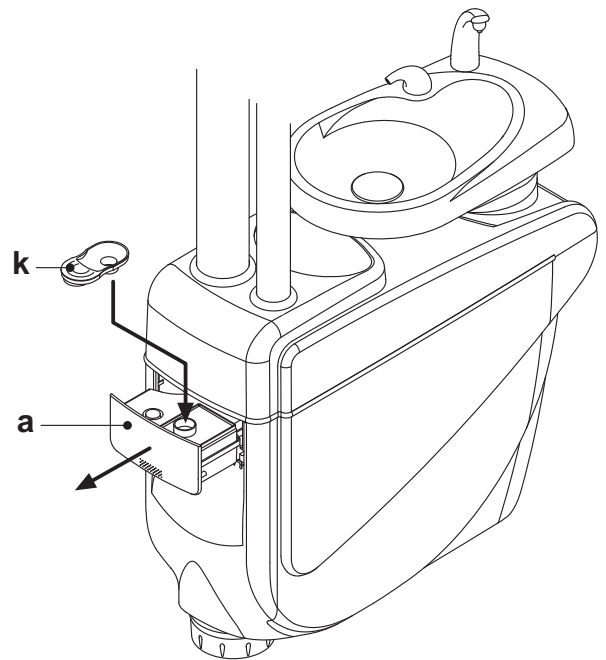
ATTENTION!

Never leave hydrogen peroxide in the tank (a) for more than one month.

When you are going to be away from the office for an extended period of time (holiday), always empty the tank (a) before leaving.

NOTE: use a suction tube to empty the tank.

7.4. AUTOSTERIL automatic disinfection system

Description of the system.

This system performs an automatic disinfection cycle for the water circuits of the following instruments:

- all instruments found on the instrument board
- the syringe and instrument placed on the assistant's board,
- the suction tubes (if the S.H.D. system is present),
- water to cup ducts

The system is equipped with a tank (a), located in the hydrogroup that holds approximately 970 cc of hydrogen peroxide with 230 cc as reserve.

The disinfection cycle can be set and is equipped with an electronic safety system in compliance with CEE 93/42, and subsequent changes, Medical Device Directive.

ATTENTION!

- Perform a disinfecting cycle at the end of each work day.

Signaling the hydrogen peroxide level is low.

When the hydrogen peroxide in the tank (a) is about to run out, a specific icon (D) appears on the doctor's console.

NOTE: the icon (D) is black if the AUTOSTERIL system is turned on while it is red if the AUTOSTERIL system is turned off.

Filling the hydrogen peroxide tank.

Follow the instructions given below to fill the hydrogen peroxide tank:

- Remove the hydrogen peroxide tank [a].
- Take off the cap [k] and pour hydrogen peroxide into the tank until it is full.

NOTE: the shape of the cap allows it to be used as a funnel to more easily fill the tank.

- Put the cap and tank back into place.
- Close the side cover of the hydrogroup.

ATTENTION!

To top up, use only 3% pure hydrogen peroxide (10 volumes) without diluting it.

Setting the disinfecting cycle.

- Check the level of liquid disinfectant in the tank, topping up if necessary.

NOTE: the disinfecting cycle will not start if the liquid in the tank is below the reserve level.

- Use the doctor's console or button AS found on the assistant's board to go to the "AUTOSTERIL cycle setup" menu and set the time the liquid disinfectant has to remain in the water ducts of the instruments (see paragraph 5.1.1.3.1.).
- Place the instruments to be disinfected in the container provided in the hydrogroup.

ATTENTION!

Use the adapter [f] provided for the syringe. The heater should be turned off.

The micromotor should be put in without the handpiece.

- In order to disinfect the suction tubes, insert the suction tube terminals in the fittings found underneath the manifold (see paragraph 7.5.).

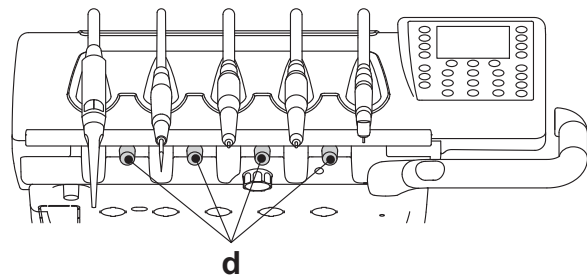
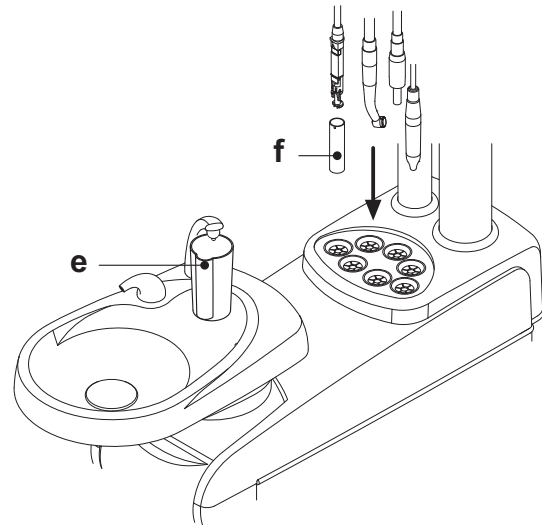
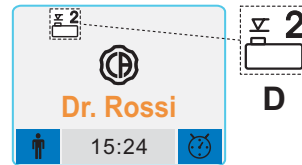
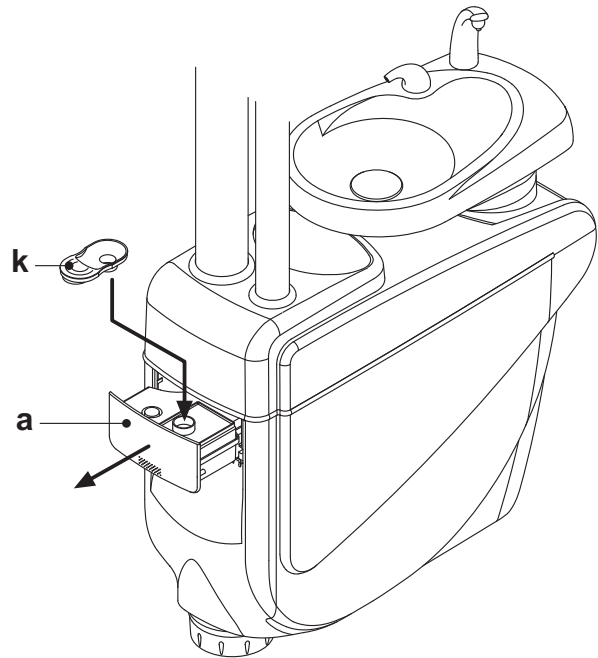
NOTE: make sure the suction tube terminals are opened.

The micromotor cord should be put in complete with the motor body.

- If the water to cup ducts need to be disinfected, place the container [e] provided under the cup spout.
- Make sure the spray cocks [d], found at the bottom of the board, are opened.

Performing the disinfecting cycle.

- Press the multifunction button (▶) on the doctor's console (see paragraph 5.1.1.3.1.) or AS button found on the assistant's board to start the disinfecting cycle.
- At this point, the system performs the following steps automatically:
 - water ducts for instruments emptied with air;
 - disinfectant let in and time it has to remain in the ducts previously set starts to clock down;
 - once this time is over, ducts are emptied with air again;
 - ducts flushed with water from the municipal line or distilled water (only if a S.S.S system is present and turned on).
- Once the disinfection cycle has been completed ("Cycle completed: put instruments back in place" appears on the display) simply put the instruments removed back into place to resume work.

**Stopping the disinfection cycle.**

- Press multifunction button STOP to interrupt the disinfecting cycle at any time.
- A confirmation message will appear on the console display:
 - Press multifunction button ESC so that the disinfecting cycle is no longer stopped. The cycle menu can now be viewed.
 - press multifunction button OK to interrupt the disinfecting cycle. An intermediate menu showing the set time and instruments withdrawn is displayed.

NOTE: At this point, the dental unit is shutdown.

- The following operations can now be performed:
 - press multifunction button ESC to go back to the initial menu used to set the times. Here it is possible to restart the disinfecting cycle from the beginning, if desired, changing the time the disinfectant has to be left in and/or adding other instruments to be disinfected.
 - press button MENU to go to the "Instrument flushing" menu to flush the ducts of the withdrawn instruments.
 - press multifunction button () to restart the disinfecting cycle from the point at which it was interrupted.
- press multifunction button () on the "Instrument flushing" menu to start the cycle to drain and flush the ducts of the instruments withdrawn using municipal water or distilled water (if the S.S.S. system is present).

NOTE: Once the disinfection cycle has been completed ("Cycle completed: put instruments back in place" appears on the display) simply put the instruments removed back into place to resume work.

How to store the hydrogen peroxide.

Follow the manufacturer's instructions given on the package to correctly store the hydrogen peroxide solution.

**ATTENTION!**

Never leave hydrogen peroxide in the tank (a) for more than one month. When you are going to be away from the office for an extended period of time (holiday), always empty the tank (a) before leaving.

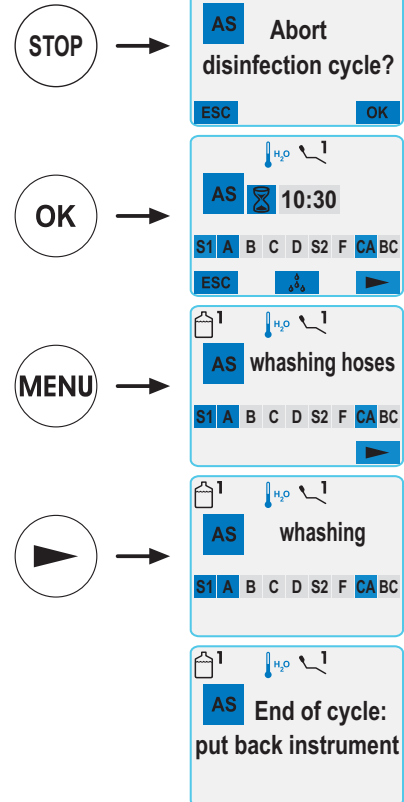
NOTE: a suction tube should be used to empty the tank. :

Error messages shown on the console display.

If the system detects a malfunction, a fault message appears on the display screen (see paragraph 11.).

**ATTENTION!**

If the disinfecting cycle is incorrectly interrupted, the equipment will shutdown until either the disinfecting cycle is performed again or the washing cycle is carried out.

7.5. S.H.D. system for flushing the suction tubes.

Description of the system.

This system flushes, washes and sterilizes the surgical suction system. The system consists of a tank [a] containing disinfectant liquid and two connectors [b] which are used when flushing the cannulae. The tank can hold 500 cc of disinfectant. The washing cycle is automatically carried out and should usually be performed at the end of each surgical procedure and whenever the dental unit is cleaned and disinfected.

ATTENTION!

Use only a diluted solution of **STER 3 PLUS (Castellini)** as the liquid disinfectant.

As regards the amount of disinfectant and percentage of solution, follow the instructions enclosed with the product.

How to start the washing cycle.

To start the washing cycle, follow the directions given below:

- Make sure there is disinfectant in the tank [a].
- Remove both suction tube terminals from the assistant's board, making sure the suction motor starts running.
- Open the mechanical closing of the suction tube terminals.
- Insert the suction cannulae terminals in the corresponding fittings [b] under the manifold.

The vacuum created by the Venutri meters triggers the washing cycle.

 **NOTE: an indication is also given on the console display to inform the user the flushing cycle has been started.**

- Washing cycle stages:
 - deliver municipal water for 50 sec. using intermittent operation (2 sec. ON - 1 sec. OFF);
 - stop the water flow and let in 10 cc of liquid disinfectant;
 - stop letting in liquid disinfectant and continue sucking for 10 sec.
 - The washing cycle ends when the suction flow is interrupted and the motor stops running.
- "Put the suction tubes back in place" appears on the display.
- At this point, put the ends of the suction tubes in the supports on the assistant's board to go back to the work conditions.

Filling the tank.

When the disinfectant [a] tank goes into reserve, proceed as follows:

- Raise the patient chair to maximum height.
- Remove the tank by rotating it anticlockwise.
- Pour disinfectant into the tank until it is full.
- Re-fit the tank; screw it on clockwise.

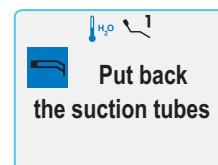
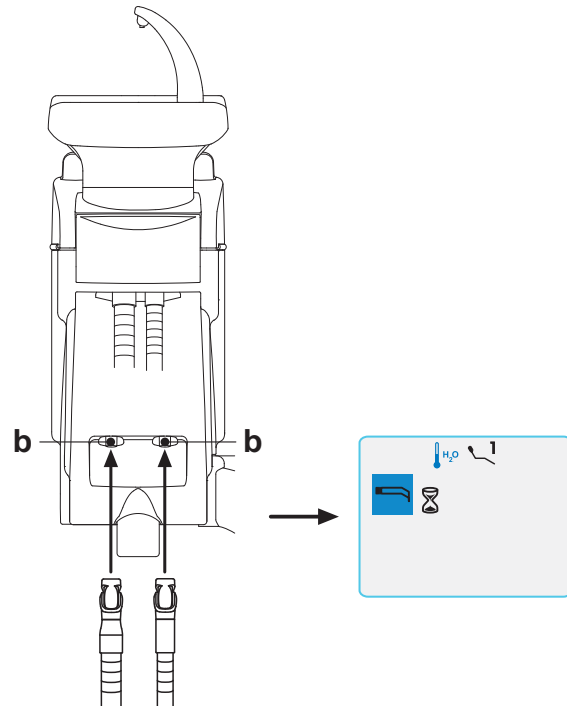
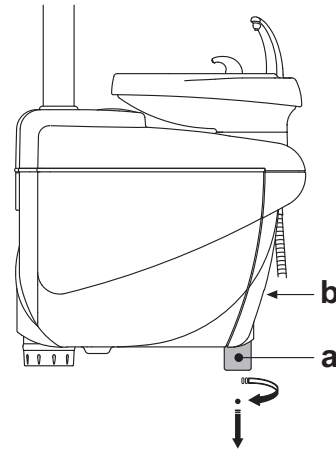
Stopping the washing cycle.

If the system detects a malfunction, a fault message appears on the display screen (see paragraph 11.).

NOTE:

Once the problem has been solved, the washing cycle automatically restarts.

7.6. TIME FLUSHING automatic cycle

Description of the system

The automatic TIME FLUSHING cycle allows an automatic flushing cycle to be carried out to renew the water present in the water lines of the instruments found on the instrument board and assistant's board.

Municipal water, treated water (if the M.W.B. system is present) or distilled water (if the S.S.S. system is present) can be used to flush the ducts. The cycle takes 2 minutes to complete.

ATTENTION!

A flushing cycle should be performed before starting work and in between patients.

Setting the TIME FLUSHING cycle.

- If the S.S.S. system is present. Check the level of distilled water present in the tank, topping up if necessary.

NOTA: : the TIME FLUSHING cycle will not start if the reserve level is present in the S.S.S. system tank.

- Use the doctor's console or button AS found on the assistant's board to go to the "TIME FLUSHING cycle setup" menu and set the time the cycle lasts (see paragraph 5.1.1.3.2.).
- Place the instruments to be treated in the container provided in the hydrogroup

ATTENTION!

Use the adapter [f] provided for the syringe. The heater should be turned off.

The micromotor should be put in without the handpiece.

The turbine and scaler should be put in without any burrs or tips.

- Make sure the spray cocks [d], found at the bottom of the board, are opened.

Performing the TIME FLUSHING cycle.

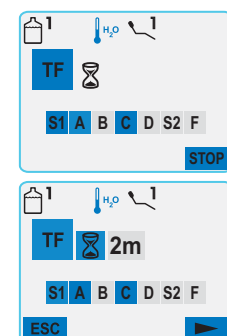
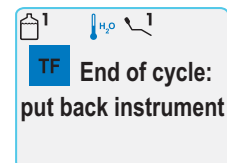
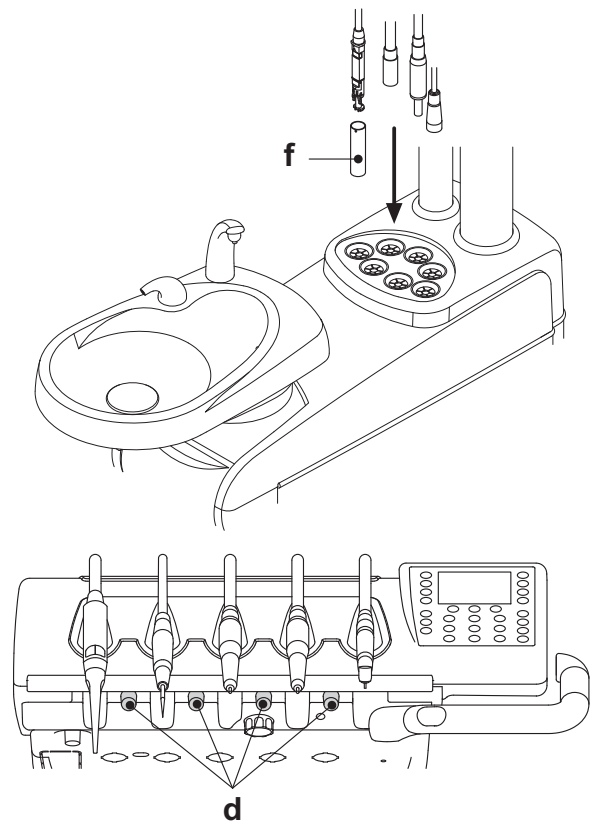
- Press the multifunction button (▶) on the doctor's console (see paragraph 5.1.1.3.2.) or AS button found on the assistant's board to start the disinfecting cycle.
- Once the TIME FLUSHING cycle has been completed ("Cycle finished: put instruments back in place" appears on the display) put the withdrawn instruments back into place to restore work conditions.

Interrupting the TIME FLUSHING cycle.

Press multifunction button STOP to interrupt TIME FLUSHING cycle at any time and go back to the initial cycle setup menu.

Error messages shown on the console display.

If the system detects a malfunction, a fault message appears on the display screen (see paragraph 11.).


7.7. Opening/closing the cover at the side of the hydrogroup

Opening the cover:

- Grasp the handle [h] and pull it up to unlock the cover [a] at the side of the hydrogroup.
- Turn the cover outwards to open it.

Closing the cover:

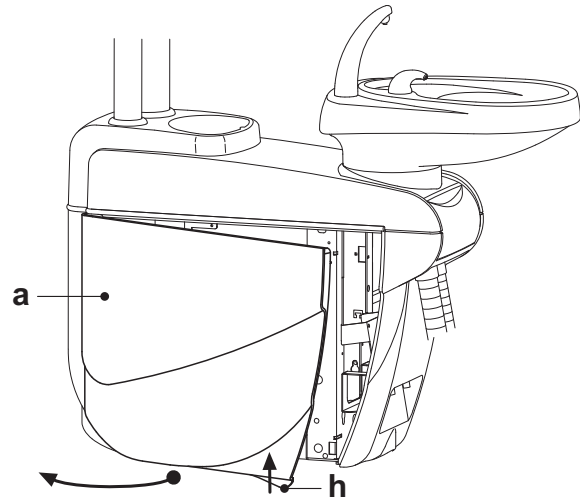
- Press the handle [h] to block the cover until a click is heard.

8. Accessories


8.1. Operating lamp

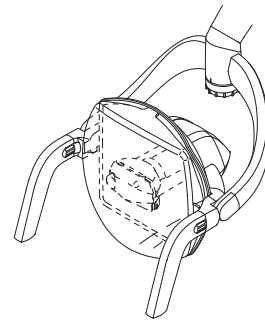Operation and care instructions for the lamp are enclosed with the equipment.

8.2. Monitor / Workstation with light pole

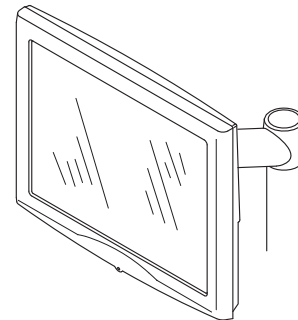Instructions on how to use the monitor/workstation are enclosed with the dental unit.

8.3. Air/water/230V quick-connect couplers

The air/water/230V quick-connect couplers are placed to the side of the electrical box.

ATTENTION!

Shut off the equipment before attempting to connect or disconnect the air/water outlets.

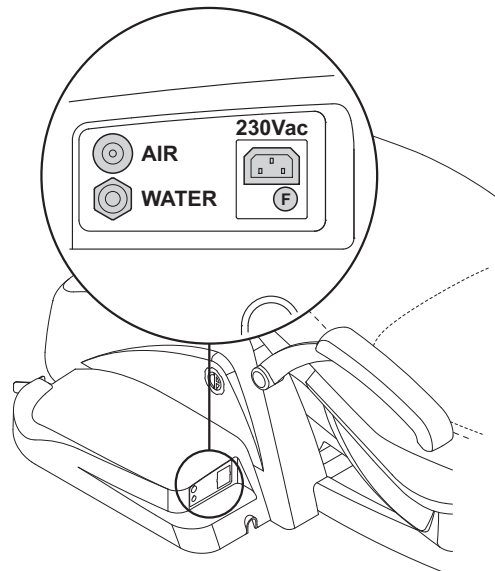
Technical specifications:

- Electrical outlet: 230Vac 2A, in accordance with standard IEC/EN 60320-2-2/F.
- Air quick-connect coupler pressure: 6 Bar.
- Water quick-connect coupler pressure:
 - municipal water, 2.5 Bar
 - with S.S.S. system, 3Bar
 - with M.W.B. system, 3 Bar
- Water quick-connect coupler delivery rate:
 - municipal water, 1800 ml/min
 - with S.S.S. system, 950 ml/min
 - with M.W.B. system, 400 ml/min

NOTE:

with the S.S.S. system: in order to use the quick-connect coupler with municipal water shut off the distilled water tank (see paragraph 7.2.).

8.4. Viewer for panoramic x-rays

A viewer for panoramic x-rays can be installed on all the CP version instrument boards.

The dimensions of the screen are: H=200mm, L=300mm.

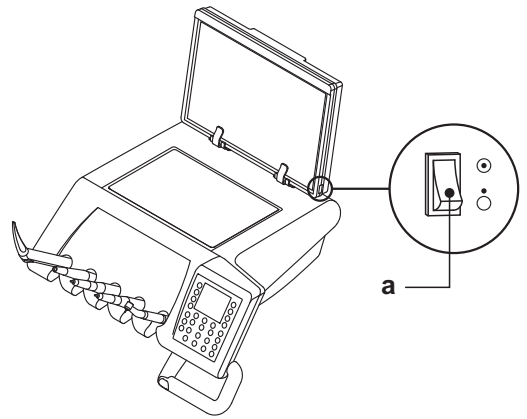
Use the switch [a] to turn on the x-ray film viewer:

= X-ray film viewer on.

= X-ray film viewer off.

8.5. H.P.A. filter (Hight Protected Air)

The H.P.A. filter [h] holds back any bacterial load present in the air line directed to the instrument sprays.

IMPORTANT INFORMATION!

- The filter element is not sterile when supplied;
- Use disposable gloves when replacing parts;
- The transparent ampule does not stand up to alcohol. Clean the ampule with neutral detergent and hot water.

Care

- Sterilize the element every month.

IMPORTANT INFORMATION!

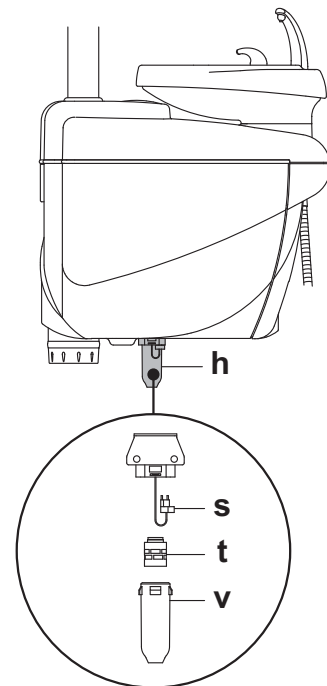
Sterilization must be done in steam autoclave at 135°C (2 bar) for at least 15 min.

Do not use dry sterilizers.

- Replace the element once a year or after 200 sterilization cycles.

Removing the filter element

- Move the dental chair completely up;
- Shut off the dental unit with the main switch (see paragraph 3.);
- Operate the syringe until the air pressure has been discharged;
- Remove the safety lock [s] found to the side of the ampule by pulling it down;
- Turn the transparent ampule [v] counter-clockwise to disengage it;
- Remove the filter element [t] by pulling it down.
- Sterilize/replace the filter element (code FH4200025).
- Put the filter element [t] back into place;
- Put the transparent ampule [v] back into place turning clockwise to engage it.
- Put the safety lock [s] back in place.

9. Maintenance

Preventive maintenance

CASTELLINI S.p.A., the manufacturer of the dental units, in accordance with applicable standards I.E.C. 60601-1 second ed., I.E.C. 62353 and directive MDD 93/42 and subsequent changes for medical devices underlines that the preventive maintenance checks for the dental unit specified in the Technical care manual and Maintenance and warranty handbook are to be carried out by authorised personnel at least once every 12 months.

ATTENTION!

The warranty is void if the equipment is serviced, repaired, altered or modified in any way by personnel who have not been duly authorized by CASTELLINI S.p.A..

Safety checks

In accordance with standard I.E.C. 62353, the safety checks specified in the Technical care manual and Maintenance and warranty handbook supplied with the dental unit are to be carried out at the intervals required by current local regulations. If no precise indications are given, CASTELLINI S.p.A., the manufacturer of the dental units, recommends checking them at least every 24 months at the time of installation and whenever electrical parts that are live are repaired/updated.

ATTENTION

The manufacturer shall not be held liable for any personal injury or equipment damage if the precautions given above are not observed.

9.1. Instrument maintenance

Maintenance instructions for the instruments are enclosed with each instrument.

ATTENTION!

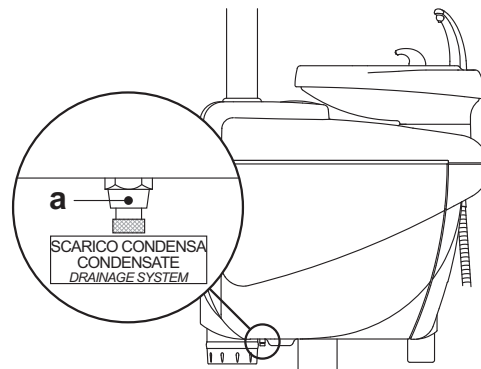
Maintenance of the instruments should be carried out with the equipment shut off.

9.2. Draining condensate

Perform this operation every day before starting work.

Proceed as follows:

- put a container under the cock [a] found below the hydrogroup,
- loosen the cock's knob
- once emptied, completely close the cock.

9.3. Cleaning the surgical suction filters

ATTENTION!

ALWAYS wear gloves when performing this operation to avoid contact with infected material.

At the beginning of each work day:

Put a VF CONTROL PLUS [v] tablet inside each filter [d].

At the end of each work day

Clean the filters as directed below:

- Open the door [c] that protects the filters.
- Take out the two filters [d].
- Clean/replace the filter (code 97290068).

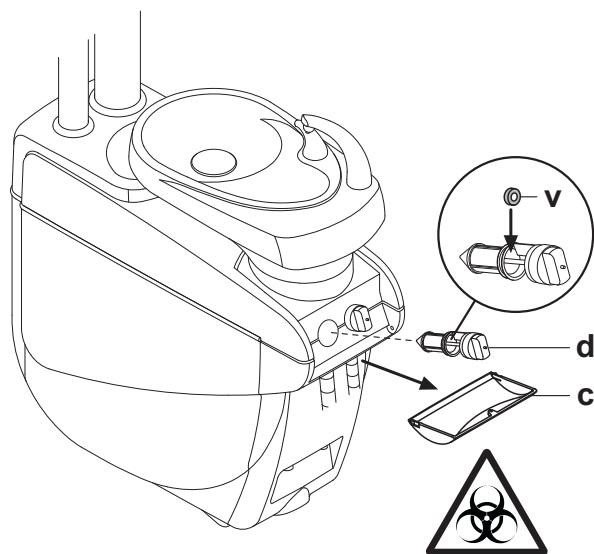
NOTE: The door [c] can be taken off and used as a tray for the filter removed.

- Put the two filters back in place being careful to remove any amalgam left at the seat of each filter.
- Close the door [c].

NOTE: in order to prevent liquids from dripping or secretions from the filter, these operations should be performed with the suction tube running.

9.4. Surgical suction

Surgical suction must be cleaned and disinfected with a disinfectant formulated for this purpose.

ATTENTION!

Use only a diluted solution of **STER 3 PLUS (Castellini)** as a liquid disinfectant.

As concerns the amount of disinfectant and percentage of solution follow the instructions provided with the product.

At the end of each surgical procedure.

- Perform an automatic flushing cycle (if present) or suck in about half a liter of the disinfecting solution with each of the suction tubes used.
- Sterilize the ends of the tube holder in autoclave with water steam at up to 135 °C (2 bar), minimum time 15 min.

At the end of each work day.

- Draw in 1 liter of water with each suction tube, alternating water and air (keep the suction tube alternately in and out of the water).
- Once rinsed with water, draw in approximately half a liter of the disinfecting solution with each of the suction tubes used.

NOTE: After disinfecting, it is good practise to draw in air only to dry the entire suction system (5 minutes).

Once a week

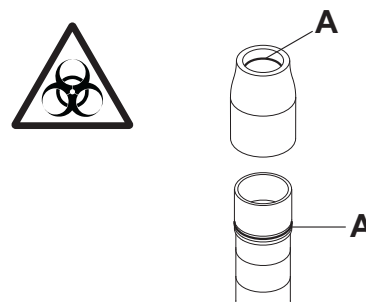
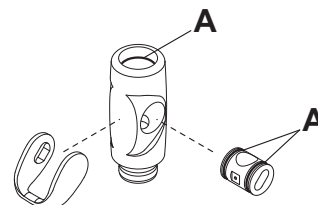
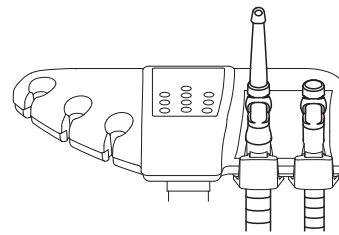
ATTENTION!

Perform this operation with the suction tube raised to keep suction on.

- Remove the suction tube body from its cord fitting and lubricate the o-rings [a] with **S1-Protection for o-rings**.

Once a year

Replace the suction tubes and ends of the tube holder.

9.5. CATTANI surgical separator

At the end of each surgical procedure.

- Suck in about half a liter of solution prepared with the selected disinfectant with each of the suction tubes used.
- Sterilize the ends of the tube holder in autoclave with water steam at up to 135°C (2 bar), minimum time 15 min.

At the end of each work day.

- Draw in 1 liter of water with each suction tube, alternating water and air (keep the suction tube alternately in and out of the water).
- Once rinsed with water, draw in approximately half a liter of the solution prepared with the selected disinfectant with each of the suction tubes used.

 **NOTE:** After disinfecting, it is good practise to draw in air only to dry the entire suction system (5 minutes).

Every 15 days.

- Clean the separator container and probes with a soft sponge and neutral detergent.
- Clean the drain valve for the separator's container with the device provided for this purpose.

Once a year.

- By technician: check the siphons and drains, check all the internal tubes and plastic and rubber parts subject to wear.

Before leaving the surgery empty for a few days.

- Start the aspirator and run it 20 - 30 minutes without sucking in liquids. The aspirator will dry itself completely. As a result, salt caused by moisture and basic substances will not form, salt that may cause fan seizure and motor blockage.

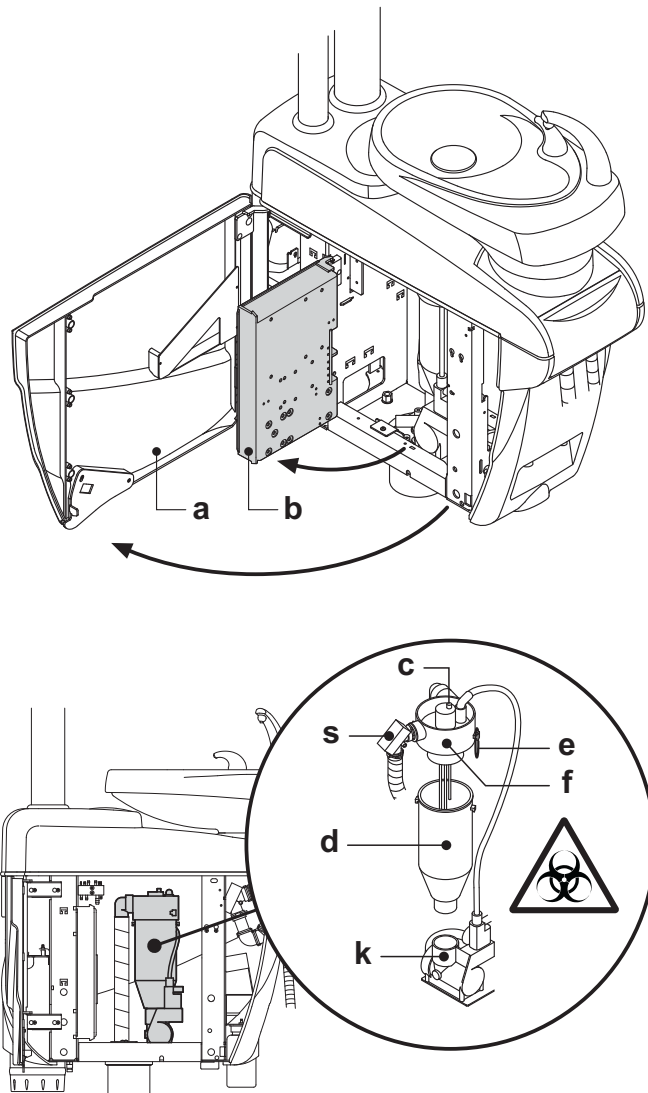
How to remove the separator's container

ATTENTION!

Gloves must be worn when carrying out the following operation to prevent contact with infected material.

- Move the dental chair fully up.
- Open the cover at the side of the hydrogroup [a] (see paragraph 7.7) and turn the electrical box [b].
- Turn the electrical box [b] and hydrogen peroxide tank [a] (if present).
- Completely empty the separator's container by pressing the timed button [c] found on the cover.
- If present, remove the valve (s) for central systems.
- Turn and raise the container until it is detached from the drain pump (k).
- Detach the container [d] from the cover [f] by pulling up the two side elastic bands [e].
- Once cleaned, put the container [d] back in place greasing the drain valve and o-ring with **S1-Protection for o-rings**.
- Put the electrical box [b] back in place and put the cover [a] back on the side of the hydrogroup (see paragraph 7.7).

9.6. Cleaning the turbine return air filter

Monthly check the oil container filter [g] present in the turbine's return air line.

If necessary, replace the filter element (code 97290014).

9.7. METASYS amalgam separator

The maintenance instructions for the METASYS amalgam separator are enclosed with the equipment if the equipment comes with this type of separator.

The separator's control device is located in the hydrogroup

ATTENTION!

Always wear gloves to prevent contact with infected material when cleaning the separator.

ATTENTION!

When disposing one-time use containers full of amalgam, observe current local and national laws.

9.8. DURR amalgam separator

The maintenance instructions for the DURR amalgam separator are enclosed with the equipment if the equipment comes with this type of separator.

The separator's control device is located in the hydrogroup

ATTENTION!

Always wear gloves to prevent contact with infected material when cleaning the separator.

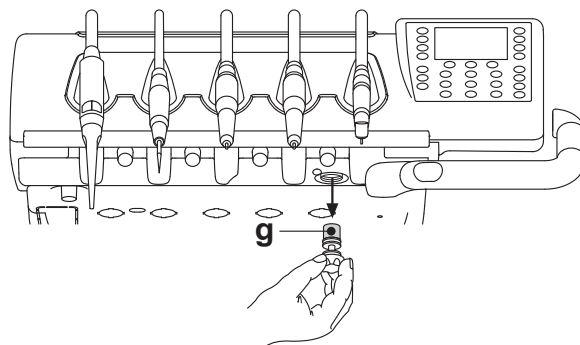
ATTENTION!

When disposing one-time use containers full of amalgam, observe current local and national laws.

9.9. CASTELLINI dental chair model NEW SKEMA

The dental chair NEW SKEMA does not need any special maintenance. In any case, operation should be checked by CASTELLINI authorized personnel once a year.

10. Packing

Packing dimensions.

Dental unit: 1570 x 780 x 1230 h. Maximum weight: 140 kg.

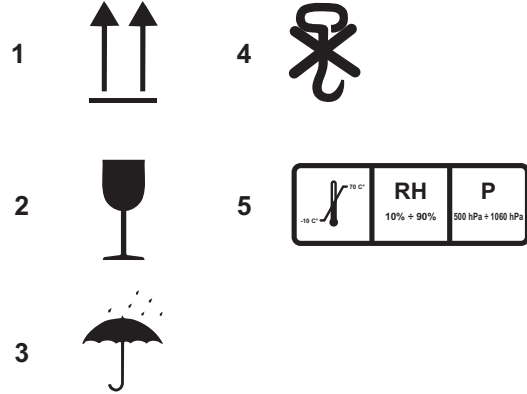
Dental chair: 1510 x 730 x 1000 h. Maximum weight: 150 kg.

10.1. Transport and storage

Storage, transport and unpacking instructions are given on the outside of the carton.

Carefully follow the instructions provided.

- 1) Transport and store only in the direction indicated by the arrows.
- 2) Avoid banging.
- 3) Protect against humidity.
- 4) Do not use hooks to lift.
- 5) Nameplate containing environmental conditions for storage.
 - temperature from -10° to 70 °C.
 - relative humidity from 10 to 90%.
 - atmospheric pressure from 500 to 1060 hPa.

11. Fault messages

M = Message shown on console display
C = Cause
R = Remedy

M: "H2O level low, fill tank"

C: The water in the independent water system's tank has dropped below the minimum acceptable level.

R: Fill the tank (see paragraph 7.2.).

M: "Put instruments back in place"

C: The system detected an instrument was already withdrawn while the disinfecting cycle was being set.

R: Make sure all the instruments are in place and then set the cycle again. If the fault message appears again, call technical support.

M: "Check instruments, repeat cycle"

C: The system detected the withdrawn instruments were altered during the flushing or disinfecting cycle.

R: Check the selected instruments and repeat the disinfecting (see paragraph 7.4.) or flushing (see paragraph 7.7.) cycle.

M: "H2O2 level low, fill tank"

C: The hydrogen peroxide in the relative tank has dropped below the minimum acceptable level.

R: Fill the hydrogen peroxide tank (see paragraph 7.4.).

M: "Open H2O spray cocks"

C: The system is not able to fill the lines with hydrogen peroxide during the disinfecting cycle.

R: Open the water spray cocks and repeat the disinfecting cycle (see paragraph 7.4.). If the message appears again, call technical support.

M: "Remove all instruments"

C: The system detected an internal malfunction during the disinfecting cycle.

R: Repeat the disinfecting cycle, selecting all the instruments. If the message appears again, call technical support.

M: "Empty M.W.B. system"

C: **M.W.B.** system malfunction.

R: Empty the tank inside the **M.W.B.** system and restart the system (see paragraph 5.1.1.7.). If the message appears again, call technical support.

M: "Turn on M.W.B."

C: The system is attempting to perform a task which requires that the **M.W.B.** system is turned on.

R: Turn on the **M.W.B.** system (see paragraph 7.2.).

M: "Check suction tubes, repeat cycle"

C: The system has detected that the suction tubes are not connected to the relative fittings during the flushing or disinfecting cycle.

R: Make sure the suction tubes are properly connected and repeat the disinfecting (see paragraph 7.4.) or flushing (see paragraph 7.5.) cycle. If the message appears again, call technical support.

M: "Withdraw at least one instrument"

C: An attempt has been made to start a disinfecting cycle without selecting any instruments or the cup.

R: Repeat the disinfecting cycle selecting at least one instrument or the cup. If the message appears again, call technical support.

M: "Instrument configured"

C: The instrument in the indicated position on the board has been automatically configured with the factory settings.

R: If the message appears again, call technical support.

M: "Put suction tubes back in place"

C: Suction tubes extracted when dental unit is turned on.

R: Make sure the suction tubes are correctly placed in their housings. If the message appears again, call technical support.

M: "Put instrument back in place"

C: Instrument extracted when dental unit is turned on.

R: Make sure all the instruments are correctly placed in their housings. If the message appears again, call technical support.

M: "Check suction tube filter"

C: Suction tube flushing cycle malfunction.

R: Make sure the filters are clean, the suction tubes are not closed and that the suction unit works correctly and then repeat the flushing cycle. If the message appears again, call technical support.

M: "Hydrogroup emergency device activated"

C: While performing an automatic movement, the assistant's board encountered an obstacle.

R: Clear the obstacle and press the button for the desired program again.

M: "Lower dental chair"

C: The bowl does not move because the dental chair is in its way.

R: Lower the dental chair so that it is no longer in the way.

M: "Check operating light fuses"

C: The operating light does not turn on because electric power is not supplied.

R: Call technical support.

M: "Maintenance required"

C: Scheduled maintenance required.

R: Call technical support to schedule the maintenance work.

M: "Footboard emergency device activated"

C: The dental chair encountered an obstacle.

R: Press the "seat up" button and clear the obstacle.

M: "Backrest emergency device activated"

C: The dental chair backrest encountered an obstacle.

R: Press the "seat up" button and clear the obstacle.

M: "Move bowl"

C: The dental chair does not move because the bowl is in the way.

R: Move the bowl so that it is no longer in the way of the dental chair.

M: "Seat emergency device activated"

C: The dental chair encountered an obstacle.

R: Press the "seat up" button and clear the obstacle.

M: "Dental chair blocked, put instrument back in place"

C: An attempt was made to move the dental chair with an instrument withdrawn.

R: Put instrument back in place and repeat dental chair movement.

M: "Dental chair blocked"

C: An attempt was made to move the dental chair while it was blocked

R: Remove dental chair blockage (see paragraph 4.).

IMPORTANT INFORMATION!

M: "XXXX, call technical support" (where XXXX represents a numerical code)

C: This type of message indicated a serious internal error.

R: Call technical support quoting the number of the error.

SKEMA 6 - OPERATING INSTRUCTION


12. Specifications.

Installation plan:	97042072	Water supply pressure:	3-5 bar
Technical manual:	97071123	Water delivery rate:	10 l/min.
Dental unit spare parts catalogue:	97023067	Water usage:	2 l/min.
Dental chair spare parts catalogue:	97023067	Water hardness:	60 mg/l
Maximum dental unit weight:	90 kg	Drain connection:	ø40 mm
Maximum dental chair weight:	115 kg	Drainage rate:	10 l/min.
Maximum load-bearing capacity of the NEW SKEMA dental chair	190 Kg.	Drain duct inclination:	10 mm/m
Voltage:	230 V~ /115 V~	Aspirator connection:	ø40 mm
Frequency:	50/60 Hz.	Vacuum (minimum):	65 mbar
Power absorbed:	1,5 kW	Vacuum delivery rate:	450 l/min.
Air connection:	1/2 Gas.	Mark of approval:	CE 0051
Air supply pressure:	6-8 bar	Electrical work in compliance with:	I.E.C. 60364-7-710
Air delivery rate:	82 l/min.		
Water connection:	1/2 Gas.		

FUSES			
Identification	Value	Protection	Position
<i>Dental unit.</i> Fuse F2	T 8 A	230 V~ : Dental unit power supply line.	Electrical box.
Fuse F4	T 8 A	115 V~ : Dental unit power supply line.	Electrical box.
Fuse F5	T 6,3 A	Secondary protection: Hydrogroup.	Electrical box.
Fuse F6	T 6,3 A	Secondary protection: Dental unit.	Electrical box.
	T 6,3 A	Secondary protection: Operatory light.	Electrical box.
<i>Dental chair.</i> Fuse FU1	T 4 A	230 V~ : Dental chair power supply line.	Electrical box..
<i>Quick-connect couplers.</i> Fuse	T 2 A	230 V~ : Electrical outlet power supply line.	Electrical box.
MONITOR power supply. Fuse	T 6,3 A	21 V~: MONITOR / WORKSTATION power line.	Dental chair card area.

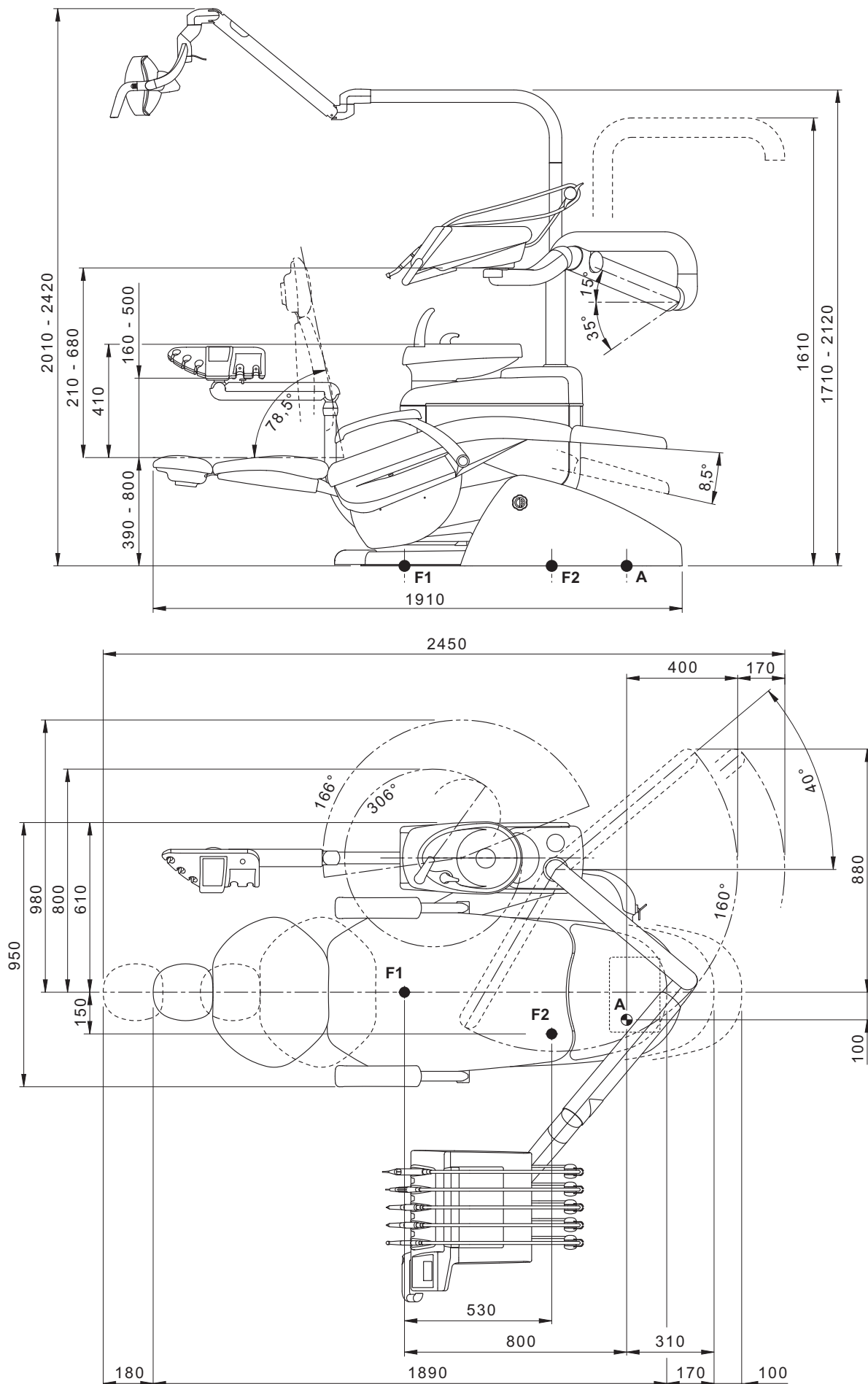12.1. Overall dimensions: SKEMA 6 RS


SKEMA 6 - OPERATING INSTRUCTION

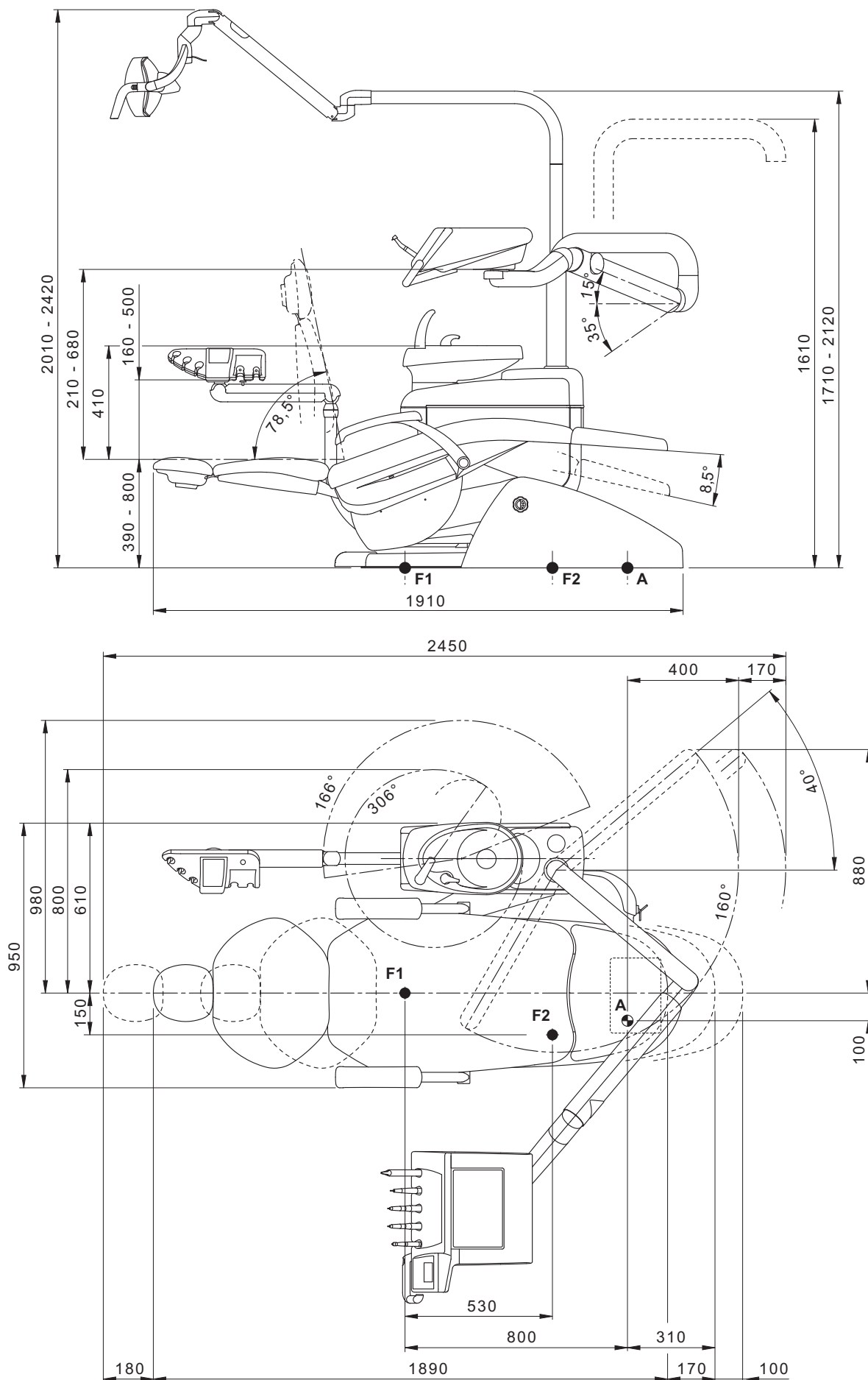12.2. Overall dimensions: SKEMA 6 CP

13. SKEMA 6 dental operatory disinfecting and maintenance plan

WHEN	PART	WHAT TO DO	SEE PARAGRAPH
Before starting work	Condensate drain cock	Drain condensate from air ducts	See paragraph 9.2
	Spray ducts (if the dental unit is turned off for more than 48 hours)	Disinfect water ducts of sprays with manual cycle (if the S.S.S. system is present)	See paragraph 7.2.1
		Disinfect water ducts of sprays with semi-automatic cycle (if the AUTOSTERIL system is present)	See paragraph 7.2
	Surgical suction filters	Place a VF CONTROL PLUS tablet inside each filter	See paragraph 9.3
	Single-use protective barriers	Put single-use barriers and flexible protection on dental chair and dental unit	/
	Strumenti.	Lubricate micromotor and turbine handpieces	See documentation enclosed with the instrument
After each treatment	Spray ducts	Flush water ducts of sprays (if the TIME FLUSHING system is present)	See paragraph 7.6
	Instruments	Sterilize or disinfect the outside	See documentation enclosed with handpiece
	Camera	Disinfect the outside Do not use harsh products or acids	See paragraph 5.8
	Surgical suction tubes	Perform an automatic flushing cycle (if present) or draw approximately half a liter of STER 3 PLUS solution with each suction tube	See paragraph 9.4
	Flexible protection	Sterilize the contaminated flexible protection	/
	Single-use barriers	Replace the single-use barriers	/
	Contaminated surfaces	Disinfect with STER 1 PLUS	/
At the end of the work day	Spray ducts	Disinfect water ducts of sprays with manual cycle (if the S.S.S. system is present)	See paragraph 7.2.1
		Disinfect water ducts of sprays with semi-automatic cycle (if the AUTOSTERIL system is present)	See paragraph 7.2
	Bowl filter	Clean the filter in running water The content must be disposed separately	See paragraph 7.1
	Bowl.	Clean with commercial detergents designed for ceramic materials. Do not use harsh products or acids.	See paragraph 7.1
	Surgical suction filter	Check the filter and replace it if the suction capacity is reduced (code 97290068)	See paragraph 9.3
	Surgical suction tubes.	Perform an automatic flushing cycle (if present) or draw approximately half a liter of STER 3 PLUS solution with each suction tube	See paragraph 9.4
	Hydraulic saliva ejector	Clean the saliva ejector filter	See paragraph 6.4
	Flexible protection	Sterilize the flexible protection present	/
	Single-use barriers	Remove the single-use barriers used	/
	Dental chair and dental unit surfaces	Clean surfaces with STER 1 PLUS	/

SKEMA 6 - OPERATING INSTRUCTION

When needed.	Removable instrument cords.	Clean with a suitable disinfectant carefully following the directions for use provided by the manufacturer. Spray the product on soft disposable paper towels. Do not use acids or harsh products.	See paragraph 5.
	Spray ducts	Disinfect water ducts of sprays with manual cycle (if the S.S.S. system is present)	See paragraph 7.2.1
		Disinfect water ducts of sprays with semi-automatic cycle (if the AUTOSTERIL system is present)	See paragraph 7.2
	CATTANI surgical separator	Clean the separator's container, drain valve and probes	See paragraph 9.5
	METASYS amalgam separator	Empty the separator bowl	See documentation enclosed with the separator
	DURR amalgam separator	Empty the separator bowl	See documentation enclosed with the separator
	Operatory light.	Clean the front glass and reflector mirror	See documentation enclosed with the equipment
	Monitor / Workstation with light pole.	Clean the surfaces	See documentation enclosed with the equipment
	Coated surfaces and dental chair upholstery.	Clean with a suitable disinfectant carefully following the directions for use provided by the manufacturer. Spray the product on soft disposable paper towels. Do not use acids or harsh products.	See paragraph 1.4
Weekly.	Suction tube holders	Lubricate the o-rings	See paragraph 9.4
Monthly.	Turbine return air filter.	Check the filter and replace it if necessary (code 97290014).	See paragraph 9.7.
Yearly.	Dental chair and dental unit	Contact the technical service department for general inspection.	/

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Castellini S.p.A. Castel Maggiore (BO) Italia